



*Graduation Project Submitted to:
Faculty of Computers and Artificial Intelligence,
Beni-Suef University, Egypt*



Melanosizer

Melanoma Detection device

Supervisor:

Dr. Mohammed Kayed

Eng: Mohamed Gamal

Presented by:

✚ Mohamed Essam Abdelmonem

✚ Shimaa Mostafa Mohamed

✚ Hader Mostafa Mohamed

✚ Youssef Abdelrazek Sayed

✚ Alaa Ahmed Ibrahim

Table Contents

| | |
|--|----|
| Purpose of documentation | 6 |
| Acknowledgements | 6 |
| Abstract | 6 |
| Contacts..... | 7 |
| 1 Chapter One (Introduction) | 7 |
| 1.1 Introduction..... | 7 |
| 1.2 Objectives..... | 9 |
| 1.3 Background..... | 9 |
| 1.4 Business impact | 10 |
| 1.5 Ways of marketing..... | 11 |
| 1.6 Summary | 13 |
| 2 Chapter two (Related work)..... | 13 |
| 2.1 Introduction..... | 13 |
| 2.2 Current and Existing systems | 16 |
| 2.2.1 Current systems..... | 16 |
| 2.2.2 Existing systems..... | 19 |
| 3 Chapter Three (System Analysis) | 20 |
| 3.1 Data Analysis | 21 |
| 3.2 Data Flow Diagrams..... | 21 |
| 3.3 System Requirements..... | 22 |
| 4 Chapter Four (Analysis and Design) | 23 |
| 4.1 Customer Requirements..... | 23 |
| 4.1.1 Main Features: | 23 |
| 4.1.2 Key Elements:..... | 23 |
| 4.1.3 Mapping table: | 23 |
| 4.2 Customer Requirements Specifications | 24 |

| | | |
|-------|--|----|
| 4.2.1 | System Context: - | 25 |
| 4.2.2 | CRS requirements:..... | 25 |
| 4.3 | Hardware Software Interface | 28 |
| 4.3.1 | Hardware specifications: | 28 |
| 4.3.2 | HSI System Context: | 31 |
| 4.3.3 | HW and SW Restrictions:..... | 31 |
| 4.3.4 | HSI Functional Requirements: | 32 |
| 4.4 | System Requirements Specifications | 33 |
| 4.4.1 | Software context: | 33 |
| 4.4.2 | SRS Requirements | 34 |
| 5 | Chapter Five (Implementation)..... | 38 |
| 5.1 | Deep learning model | 38 |
| 5.1.1 | Why CNN | 38 |
| 5.1.2 | CNN Features | 38 |
| 5.1.3 | Why TensorFlow..... | 39 |
| 5.1.4 | About used dataset | 40 |
| 5.1.5 | Balancing class weights | 41 |
| 5.2 | Image standardization | 45 |
| 5.3 | Benefits of using pretrained model | 46 |
| 5.4 | Loss and accuracy matrix and reasons behind them | 49 |
| 5.5 | Experiments..... | 51 |
| 5.5.1 | Some explanations before start experiments | 51 |
| 5.5.2 | Base Model..... | 52 |
| 5.5.3 | Experiment <i>one</i> | 56 |
| 5.5.4 | Experiment two | 60 |
| 5.5.5 | Experiment three..... | 64 |
| 5.5.6 | Experiment Four | 73 |

| | | |
|-------|--|-----|
| 5.5.7 | Experiment five | 77 |
| 5.5.8 | Experiments Six | 80 |
| 5.6 | Deployment..... | 89 |
| 5.6.1 | API deployment..... | 89 |
| 5.6.2 | Raspberry PI deployment | 90 |
| 5.7 | Embedded | 91 |
| 5.7.1 | Design Constraints:..... | 91 |
| 5.7.2 | Hardware Environment | 91 |
| 5.7.3 | Raspberry Pi TOUCH SCREEN LCD HDMI..... | 91 |
| 5.8 | Raspberry Pi computer with a Camera Module port..... | 93 |
| 5.9 | Software Environment..... | 97 |
| 5.9.1 | Operating system | 97 |
| 5.9.2 | Updating Raspberry Pi OS..... | 97 |
| 5.9.3 | Using APT | 98 |
| 5.9.4 | Using rpi-update..... | 98 |
| 6 | Chapter Six (Testing and Validation) | 98 |
| 6.1 | Module Testing..... | 98 |
| 6.1.1 | Model Testing..... | 99 |
| 6.1.2 | Web Server Testing | 99 |
| 6.2 | Overall Testing..... | 101 |
| 7 | Chapter Seven (Tools and Technologies) | 102 |
| 7.1 | Libraries and environment used in ML model..... | 102 |
| 8 | Chapter eight (Future work) | 104 |
| 8.1 | Conclusion | 104 |
| 8.2 | Future Work | 104 |
| 9 | References..... | 106 |

List of Figures

| | |
|---|-----|
| Figure 3. 1 Data Flow Diagram | 21 |
| Figure 4. 1 Context System Diagram | 25 |
| Figure 4. 2 Hardware Components | 31 |
| Figure 4. 3 System flow diagram | 33 |
| Figure 5. 1 benign and malignant counter..... | 42 |
| Figure 5. 2 computing class weights..... | 45 |
| Figure 5. 3 TP vs. FP rate at different classification thresholds..... | 50 |
| Figure 5. 4 Predictions ranked in ascending order of logistic regression score..... | 50 |
| Figure 5. 5 base model summary | 54 |
| Figure 5. 6 Model accuracy | 55 |
| Figure 5. 7 Model Loss | 55 |
| Figure 5. 8 Model summary | 64 |
| Figure 5. 9 Raspberry Pi TOUCH SCREEN LCD HDMI | 91 |
| Figure 5. 10 Raspberry Pi to the back of the Touch Display..... | 92 |
| Figure 5. 11 camera module port..... | 93 |
| Figure 5. 12 Raspberry Pi Camera Module | 94 |
| Figure 5. 13 shows step 1,2..... | 94 |
| Figure 5. 14 show steps 3,4..... | 95 |
| Figure 6. 1 False Negative values matrix | 99 |
| Figure 6. 2 AUC matrix | 99 |
| Figure 6. 3 Swagger UI documentation of API | 100 |
| Figure 6. 4 Example of skin image prediction..... | 100 |
| Figure 6. 5 Screen on input image..... | 101 |
| Figure 6. 6 Screen on output prediction | 102 |

Purpose of documentation

Imagine that many of the tools that today we use its development team did not attach documentation with it, this leads to the inability to reuse some of these tools.

This is inconsistent with our goal of this project, which is to be useful to both the consumers and the developers working on similar projects.

The purpose of this documentation is to record all aspects of the project and to clarify the faces of software development life cycle such as analysis, design of software or hardware using manuals, listings, diagrams, and other hard- or soft-copy written and graphic materials, integration, testing and so on, this will achieve our target to make this project maintainable and reusable.

Acknowledgements

Our deepest gratitude goes to my supervisor **Dr. Mohamed Kayed**, for his guidance, advice and constructive comments on our project, requirement specification and for his valuable support. Besides, he also guided us by telling us how to build a good system. He also provides a lot of feedback which increases the overall quality of this project. We would not have reached this phase, if it were not for his permanent support, advice, and guidance

not forgotten, our appreciation to **Dr. Mohamed Gamal**, for providing us with technical help to get better implementation which without it the experiments may be impossible to realize.

We would also like to thank our true friends, family and colleagues for giving us encouraging words to keep up trying to reach results when we get stuck in problems

Abstract

Our project is an embedded device that work with data science model to help doctors in early and accurate detection of melanoma. The early detection can make treatment more effective and save more lives.

classifiers based on convolutional neural networks (CNNs) were shown to classify images of skin cancer on par with dermatologists, could enable lifesaving and fast diagnoses, and we want this process to be inside the hospital to let the doctors take model output with other medical analysis and diagnose the type of skin cancer that the patient infected with, so we go to embedded systems to produce a medical device.

Approaches that use a CNN and optimize its parameters to the classification of skin lesions are the most common ones used and they display the best performance with the currently available limited dataset, CNN's display a high performance as skin lesion classifiers and use raspberry PI as CPU with camera and touch screen to make the device more usable.

Our project aims to make the medical field more powerful to save more lives. This will achieve by producing a model with good results and always work to improve these results to make the device as the main actor in its field.

Contacts

Team member

| Name | Email |
|--------------------------|------------------------------|
| Mohamed Essam Abdelmonem | m.esamelbanna@gmail.com |
| Youssef Abdelrazek Sayed | yussefabdelrazik@gmail.com |
| Shimaa Mostafa Mohamed | shimaa.mostafaa.07@gmail.com |
| Alaa Ahmed Ibrahim | roo7.elfahham@gmail.com |
| Hader Mostafa Mohamed | haderms2000@gmail.com |

Supervisor

| | | |
|--------------------|-----------------------|----------------|
| Dr. Mohammed Kayed | kayed@fcis.bsu.edu.eg | +2 01091666401 |
|--------------------|-----------------------|----------------|

1 Chapter One (Introduction)

1.1 Introduction

The world's most common cancer is a relentless disease that strikes one in five people by age 70. The good news is that 99 percent of all cases are curable if they

are diagnosed and treated early enough. But in order to stop skin cancer, we have to spot it on time.

Skin cancer — the abnormal growth of skin cells — most often develops on skin exposed to the sun. But this common form of cancer can also occur on areas of your skin not ordinarily exposed to sunlight. There are four main types of skin cancer — basal cell carcinoma, squamous cell carcinoma, merkel cell cancer and melanoma.

- Basal cell carcinoma. Basal cells are the round cells found in the lower epidermis. About 80% of skin cancers develop from this type of cell. These cancers are described as basal cell carcinomas. Basal cell carcinoma most often develops on the head and neck, although it can be found anywhere on the skin. It is mainly caused by sun exposure or develops in people who received radiation therapy as children. This type of skin cancer usually grows slowly and rarely spreads to other parts of the body.
- Squamous cell carcinoma. Most of the epidermis is made up of flat, scale-like cells called squamous cells. Around 20% of skin cancers develop from these cells, and these cancers are called squamous cell carcinomas. Squamous cell carcinoma is mainly caused by sun exposure, so it may be diagnosed on many regions of the skin. It can also develop on skin that has been burned, damaged by chemicals, or exposed to x-rays. Squamous cell carcinoma is commonly found on the lips; at sites of a long-standing scar; and on the skin outside the mouth, anus, and a woman's vagina. About 2% to 5% of squamous cell carcinomas spread to other parts of the body.
- Merkel cell cancer. Merkel cell cancer is a highly aggressive, or fast-growing, rare cancer. It starts in hormone-producing cells just beneath the skin and in the hair follicles. It is usually found in the head and neck region. Merkel cell cancer may also be called neuroendocrine carcinoma of the skin.
- Melanoma. There are scattered cells called melanocytes where the epidermis meets the dermis. These cells produce the pigment melanin, which gives skin its color. Melanoma starts in melanocytes, and it is the most serious type of skin cancer. It accounts for about 1% of all skin cancers.

As we mentioned above that melanoma is the least common skin cancer but it, specifically, is responsible for 75% of skin cancer deaths. The American Cancer Society estimates over 100,000 new melanoma cases will be diagnosed in 2020. It's also expected that almost 7,000 people will die from the disease. As with other

cancers, early and accurate detection—potentially aided by our project—can make treatment more effective.

Currently, dermatologists evaluate every one of a patient's moles to identify outlier lesions or “ugly ducklings” that are most likely to be melanoma. Existing AI approaches have not adequately considered this clinical frame of reference. Dermatologists could enhance their diagnostic accuracy if detection algorithms take into account “contextual” images within the same patient to determine which images represent a melanoma. If successful, classifiers would be more accurate and could better support dermatological clinic work.

As we mentioned above, melanoma is the most dangerous type of skin cancer and early and accurate detection of melanoma can make treatment more effective, so we choose to work on detecting this type in an attempt to save more lives. Our idea is providing a better support to dermatological clinic work by developing embedded device that takes a skin picture and classify it using data science model. This device helps doctors in early and accurate detection, which is the most important thing to achieve early recovery.

1.2 Objectives

The goal of the project is to help doctors detect whether skin cancer is benign or malignant and detect it early so that it can be treated in a simpler way and in a shorter time than usual. This early diagnosis and treatment save many more lives because many deaths occur as a result of late diagnosis and treatment.

1.3 Background

Skin cancer is the most prevalent type of cancer. Melanoma, specifically, is responsible for 75% of skin cancer deaths, despite being the least common skin cancer. The American Cancer Society estimates over 100,000 new melanoma cases will be diagnosed in 2020. It's also expected that almost 7,000 people will die from the disease. As with other cancers, early and accurate detection—potentially aided by data science—can make treatment more effective.

Currently, dermatologists evaluate every one of a patient's moles to identify outlier lesions or “ugly ducklings” that are most likely to be melanoma. Existing AI approaches have not adequately considered this clinical frame of reference. Dermatologists could enhance their diagnostic accuracy if detection algorithms take

into account “contextual” images within the same patient to determine which images represent a melanoma. If successful, classifiers would be more accurate and could better support dermatological clinic work.

1.4 Business impact

In this section we will discuss the potential impacts of an interruption to critical business operations, due to disasters, accidents or emergencies.

In our project we have two types of these business operations:

1. Impact in business operations related to hardware components.

This impact will affect our ability to produce more devices from our own devices. To solve this problem, we may resort to two ways of thinking.

- First, we may limit our reliance on certain hardware components. This will make the user of the device dependent on some of their smart devices, and we will implement an integration method between these devices and our device to try to reduce their dependence on some hardware components.
- Second, what if these factors (disasters, accidents, emergencies) affect many companies and a global crisis occurs, this will definitely affect the number of smart devices in the world and we may not be able to implement this integration so in this case we may resort to making our virtual model that enable you to take a picture with any camera and we will work hardly on our software to give you the best result that will help doctor to give you the best diagnosis.

2. Impact in business operations related to Cloud storage.

This effect will affect the way we deploy the model to the cloud to give old and new devices the advantage of an immediate effect to improve the model, so we will be forced to abandon this feature and deploy the model to the processor of the device. The device or the development event in the model.

The other solution to this impact depends on any level of business we are in. If it is a high level, we will work on developing our own cloud servers to serve the old and new devices, this will be a huge advantage to us that will enable us to avoid future impacts even if it related to hardware component impacts.

1.5 Ways of marketing

Steps to product marketing were Searching for opportunities, defining a target market, developing a unique offer, and selecting channels are all important steps in product marketing, we will explain in this section how we will work on marketing using the above-mentioned steps.

Here are the seven steps we need to follow to market the product:

1. Conduct market research

Information should be collected about customers who use the device and what their buying patterns and location are.

In addition to this research, it is possible to know the first sales of the device, the trend of the market to the product, and to know the actions of competitors

2. Profiling our target markets

It is necessary to know the target group, which is the category of skin cancer patients, as it is effective but expensive.

And since the main reason for using the device is the ability to know whether the patient has skin cancer or not

The device will be used in hospitals, centers and institutes for the treatment of cancer

3. Identifying our unique selling proposition (USP)

USP is the reason that distinguishes us from other competitors and why our customers buy from us and not from the competitors.

As it examines the skin and determines the part affected by the tumor and knows whether the tumor is benign or malignant, it gives the doctor the opportunity to diagnose the patient more accurately and give him the appropriate treatment.

4. Choosing our marketing avenues

Through the many targeted ways to market the device, marketing is done through social media, print ads, networking events and a business website.

5. Set your goals and budget

Marketing objectives make us get what we want to achieve through marketing activities. Our objectives are smart: specific, measurable, achievable, relevant, and time-based.

As we set a budget for our marketing activities such as:

- Website development and maintenance
- Search engine optimization strategy
- Brand Design
- Printing promotional materials (business cards, brochures, banners, etc.)
- Advertisement
- Donations and sponsorship
- Hiring employees to carry out marketing activities

6. Nurture your loyal customers

We take care of distinguished customers and offer them offers as exceptional customer service is provided that can keep people coming back and differentiate me from my competitors.

Strategies to build loyalty in customers include:

- communicating regularly with customers through social media, blogs or e-news
- providing after-sale follow up
- delivering on your promises
- going the ‘extra mile’ and providing benefits that exceed initial expectations
- using feedback and complaints as an opportunity to improve services
- listening to customers
- training staff in customer service and basic sales processes.

7. Monitor and review

We regularly monitor and review marketing activities to determine if they are achieving the desired result, such as increasing sales.

Initially, we review the marketing plan every three months to ensure that our activity supports our strategy. Once the business is more established, we review our plan when introducing a new product or service, if a new competitor enters the market.

1.6 Summary

Skin cancer is the most prevalent type of cancer. Melanoma, specifically, is responsible for 75% of skin cancer deaths, despite being the least common skin cancer. The American Cancer Society estimates over 100,000 new melanoma cases will be diagnosed in 2020. It's also expected that almost 7,000 people will die from the disease. As with other cancers, early and accurate detection—potentially aided by data science—can make treatment more effective.

Currently, dermatologists evaluate every one of a patient's moles to identify outlier lesions or “ugly ducklings” that are most likely to be melanoma. Existing AI approaches have not adequately considered this clinical frame of reference. Dermatologists could enhance their diagnostic accuracy if detection algorithms take into account “contextual” images within the same patient to determine which images represent a melanoma. If successful, classifiers would be more accurate and could better support dermatological clinic work.

2 Chapter two (Related work)

2.1 Introduction

We will talk about the literature work in two fields traditional medical analysis and using computer vision technologies:

- Traditional medical ways

Melanoma is often called "the most serious skin cancer" because it tends to spread.

- Melanoma can develop within a mole that you already have on your skin or appear suddenly as a dark spot on the skin that looks different from the rest.
- Early diagnosis and treatment are crucial.
- Knowing the [ABCDE warning signs of melanoma](#) can help you find an early melanoma.



The melanomas are detected through DERMOSCOPY

Dermoscopy of superficial melanoma

By the time in situ and invasive superficial spreading melanoma (SSM) is recognized as a changing or distinctive lesion by the patient or their doctor, it is often large (>6mm). Characteristically, superficial melanoma is asymmetrical and irregular in shape and structure.

Superficial melanomas usually have one or more of the following dermoscopic features:

- Blue-white veil
- Multiple brown dots
- Pseudopods
- Radial streaming
- Scar-like depigmentation
- Peripheral black dots/globules
- Multiple (5-6) colours, especially red and blue
- Broad network
- Focal sharply cut-off border
- Negative network
- Irregular vascularity

Uniform blue-white structures may be observed over some blue naevi and haemangiomas but in melanoma they are focal, asymmetrical and irregular.

Scar-like depigmentation due to regression of melanoma results in irregular white areas that must be distinguished from the uniform peripheral loss of pigment seen in benign halo naevi. It arises in about 50% of melanomas.

Negative network, although a feature of melanoma, may also be found in some benign melanocytic lesions (especially Spitz naevus) and seborrhoeic keratoses.

Melanoma may be recognized when there are only 2-3 colours in the lesion on dermoscopy (or 1-2 clinically). Deeper melanomas reveal more colours.



Melanoma



Irregular shape/structure

Computer vision technology

The research of skin cancer detection based on image analysis has advanced significantly over the years. Many different techniques have been tried. The International Skin Imaging Collaboration (ISIC) event of 2018 has become a de facto benchmark in skin cancer detection by hosting a challenge contest. It is also reported that a mobile app can be used to detect skin cancer. In all these efforts researchers have tried to improve the accuracy of diagnosis by employing different classification algorithms and techniques. Image classification took to new bounds when convolutional neural network (CNN) structure was introduced by Fukushima (1988) and later Le-Cunn (1990).

Because of Some of the structural features may be subtle in early melanoma We present an embedded software system that implements based on data science model to classify images of skin (naevus) on par with dermatologists, could enable lifesaving and fast diagnoses, even outside the hospital, Real-time multimodal registration to enable dual-camera spatio-temporal feature extraction in a skin cancer screening application.

We test the system on a combination of image sequences on data science model and which takes via a camera, then Image registration is performed by matching common features between each frame of organ image to each frame of another image sequence of camera to detect naevus is benign or malignant.

2.2 Current and Existing systems

In this section we will discuss firstly Describe the system currently available mostly within our chosen market segment (Current systems), and systems that exists (not necessarily within your chosen market segment) which is similar to your intended project (Existing systems).

2.2.1 Current systems

Our project is a melanoma detection device, that is, an embedded system device with a data science model that will be developed to detect whether skin cancer is benign or malignant.

In this section, we will discuss relevant works in our field (in our market). Hence, we will focus on data science models that have been published for the same purpose. Many or all these models are deployed in mobile apps which we will talk about next.

A handful of smartphone apps and devices claim to aid early detection and keep you on track with regular self-exams. You can capture photos of suspicious moles or marks and track them yourself or send them off to a dermatologist for assessment. Either way, these apps can be helpful, but they do have limitations.

A handful of skin cancer detection apps popped up allowing you to analyze your skin with your smartphone and artificial intelligence algorithms.

Some send photos to a dermatologist; some provide instant feedback and others offer helpful reminders about self-checking your skin and scheduling a doctor's appointment. Here's what you need to know about using your smartphone to detect skin cancer.

1. SkinVision

The Amsterdam-based company developed a smartphone app to easily evaluate risk factors for skin cancer and keep track of potentially problematic moles. The app uses

a deep learning algorithm to analyze your mole photo and assess whether it is high-risk within a minute.

2. UMSkinCheck

This University of Michigan Medical School app is free. UMSkin Check allows users to do a complete skin cancer exam and track changes over time. This app provides helpful advice on how to perform a skin exam. The app stores your baseline photos for comparison. It also furnishes prompts to remind you to check your skin regularly.

3. Miiskin

Miiskin is a paid "aka premium" app. It uses high-res photography to take photos of large parts of your body. The app allows the user to compare individual moles over time to detect changes.

You can take a photo of any new spots or moles and mark where they are on the body. You can add more notes and view them side by side to see if there have been any changes. Miiskin uses mole mapping to analyze your skin. Dermatologists perform mole maps as part of a clinical full-body skin exam, using digital dermoscopy (magnified digital photography) to catch suspicious lesions they may not catch with their own eyes. The Miiskin app also provides information about melanoma, risk factors and what to look for when doing your at-home check.

The app isn't free, it costs \$3.49 per month or \$34.99 for a full year.

4. MoleScope

Developed by Maryam Sadeghi during her PhD research in computer science several years ago. This is a high-resolution camera compatible with many different smartphones. This camera uses high magnification and special lighting to take more detailed and better-quality photos than other skin cancer apps. It also contains many features that other apps do, such as skin mapping, image management, and regular reminders. Though Mole-Scope itself won't analyze or diagnose your moles, you can use the ABCD guide in the app to keep tabs on any suspicious moles. The app helps

you document your moles with photos and sends them to a dermatologist, who can assess them using the ABCD method.

5. Mole Monitor

The Mole Monitor app enables you to track moles. It has the same features in which you can photograph moles, record their location, and add notes. You can view your mole images side by side (with the other apps, you have to swipe between the photos). There is also an overlay feature to help you compare mole size, but it is not as accurate as a doctor's image with a ruler imposed upon it.

You can export photos of your mole to email, but with an added benefit. You can answer questions about your risk factors (family history, personal history of sunburns) and the app generates a pdf report with your risk factors, photos and notes. You can email it to yourself or to your doctor's office before your appointment.

But the app itself is not as user-friendly as the other mole-tracking apps. The camera requires you to hold the phone at a certain distance from your skin, wait for a circle to appear on the screen, and then pull the phone back before it will take a photo. The idea is to improve image capture and the images were not any better than the other apps.

It's a free app, but you can purchase upgrades with extra imaging features.

6. Dermatology A to Z

The American Academy of Dermatology developed the Dermatology A to Z, specifically designed to serve consumers looking for skin health information. The app gives users evidence-based, dermatologist-approved health information, insights on diseases affecting skin, hair, and nails, and the latest medical and cosmetic treatments. Utilizing the smartphone's GPS tracking system, the app can show the UV Index in real time to fight against dangers of ultraviolet radiation as well as find the nearest dermatologist in the area.

7. DermaChecker

Derma-checker is a free skin cancer app that analyzes skin spots and provides an assessment using artificial intelligence within seconds. The algorithm helps to identify atypical moles which may show signs of skin cancer. Derma-Checker checks the skin spots for any signs of skin cancer. It analyzes the photo and gives a full assessment on the mole. This makes it possible to detect melanoma at an early stage when it is the most treatable.

8. NgoziYangu

The mobile app, NgoziYangu (a Swahili word that means my skin), has functions to acquire images either directly through the mobile camera or by uploading an image that had already been taken and stored. It also allowed our research assistant to add patient information (eg, sex, contacts, date of birth, and so on) and skin lesion characteristics.

9. MoleMapper

Mole Mapper is the result of a cancer biologist's efforts to help his wife. Oregon Health & Science University collaborated with Apple and Sage Bionetworks to develop this app. It's available at no cost. The Oregon Health & Science University guides physicians to help monitor suspicious lesions without monthly in-person visits. It allows users to take photos and gather measurements of any moles on their body.

Finally, Smartphone apps cannot diagnose skin cancer. It can help you with at home checks, and it can help you track moles and lesions so that you can tell your doctor about any changes or catch changes sooner for an earlier diagnosis. But current algorithm-based smartphone apps cannot be relied on to detect all cases of melanoma. Test performance is likely to be poorer than reported here when used in clinically relevant populations and by the intended users of the apps.

2.2.2 Existing systems

Skin cancer diagnosis includes:

- A physical exam: A physician, often a dermatologist, examines your skin to look for suspicious growths. Basal and squamous cell carcinomas tend to look red and flaky. Melanomas are often larger and multi-colored.
- Mole mapping: This procedure uses full-body photographs of your skin to monitor changes in existing moles and to detect new ones. Physicians use this cutting-edge preventive measure at your annual skin exam. It's especially helpful for tracking moles for high-risk individuals, such as those with a family history of melanoma.
- Biopsy: We remove a small amount of tissue from your skin, for examination under a microscope. A biopsy can determine if cancer is present. During a biopsy, we give you an anesthetic so you will be comfortable. Biopsy procedures include:
 - Shave biopsy, using a tool similar to a razor
 - Punch biopsy, using a circular tool that removes a small section of skin
 - Excisional biopsy, using a scalpel to remove an entire lump
- Sentinel lymph node mapping: Our specialists use this minimally invasive melanoma diagnostic tool to remove a sentinel node near cancerous tissues. A pathologist who specializes in skin cancer determines whether the node contains cancer cells and if the cancer has spread. Right now, we are using this for melanoma and other aggressive skin cancers including Merkel cell carcinoma. Our researchers are pioneering this technique for squamous cell cancer as well.
- Skin cancer signs and symptoms: Skin cancer can occur anywhere on your body, either in normal skin or in a mole that becomes cancerous. In men, skin cancer most often appears on the face or trunk. For women, this type of cancer most often appears in the lower legs. For both men and women, skin cancer can occur on skin that has not been exposed to sunlight.
Skin cancer can affect people of any color. In people with darker skin, skin cancer tends to occur on the palms of the hands, the soles of the feet, or under the fingernails or toenails.

3 Chapter Three (System Analysis)

In this chapter we explain the system in detail that be the set of functionalities that the system it's meant to encompass, that will be provided to the users. It defines what the system will do and what the system won't do.

3.1 Data Analysis

A diagram of the data flow is called tools, a style of modeling Important analysis and building information can be noted that a major user who logs on to the system and start registration the main picture then The algorithm is implemented using a combination of embedded software and dedicated hardware units on a heterogeneous reconfigurable system-on-a-chip ,Control between the API and machine learning and the raspberry pi connects to the raspberry pi came to record main picture thus to performs feature detection and extraction, while the software estimates the transformation parameters and maps each visible video , also connects raspberry pi to power supply.


3.2 Data Flow Diagrams

The system consists of Raspberry pi 4 Model B 4GB RAM UK running a flavor of OPERATING SYSTEM (RASPBRRY).

Have several modules for efficient functioning over a small embeddedcomputer.it is hosting a CD CARD to store all its files and the Raspberry Pi OS operating system and allow the user various controls of the raspberry pi via access it. A Raspberry pi 4 Model B based processing controller unit board and UCB _c power supply provides connectivity to various cd card running the RS Raspberry pi NOIR Camera, Raspberry pi 7capacitive Touch Screen LCD HDMI Buttons and Buzzer. CD CARD is directly connected to the processing controller unit board and relays

Figure 3.1:(Data Flow Diagram

information to the Application programming interface that connect to machine learning.

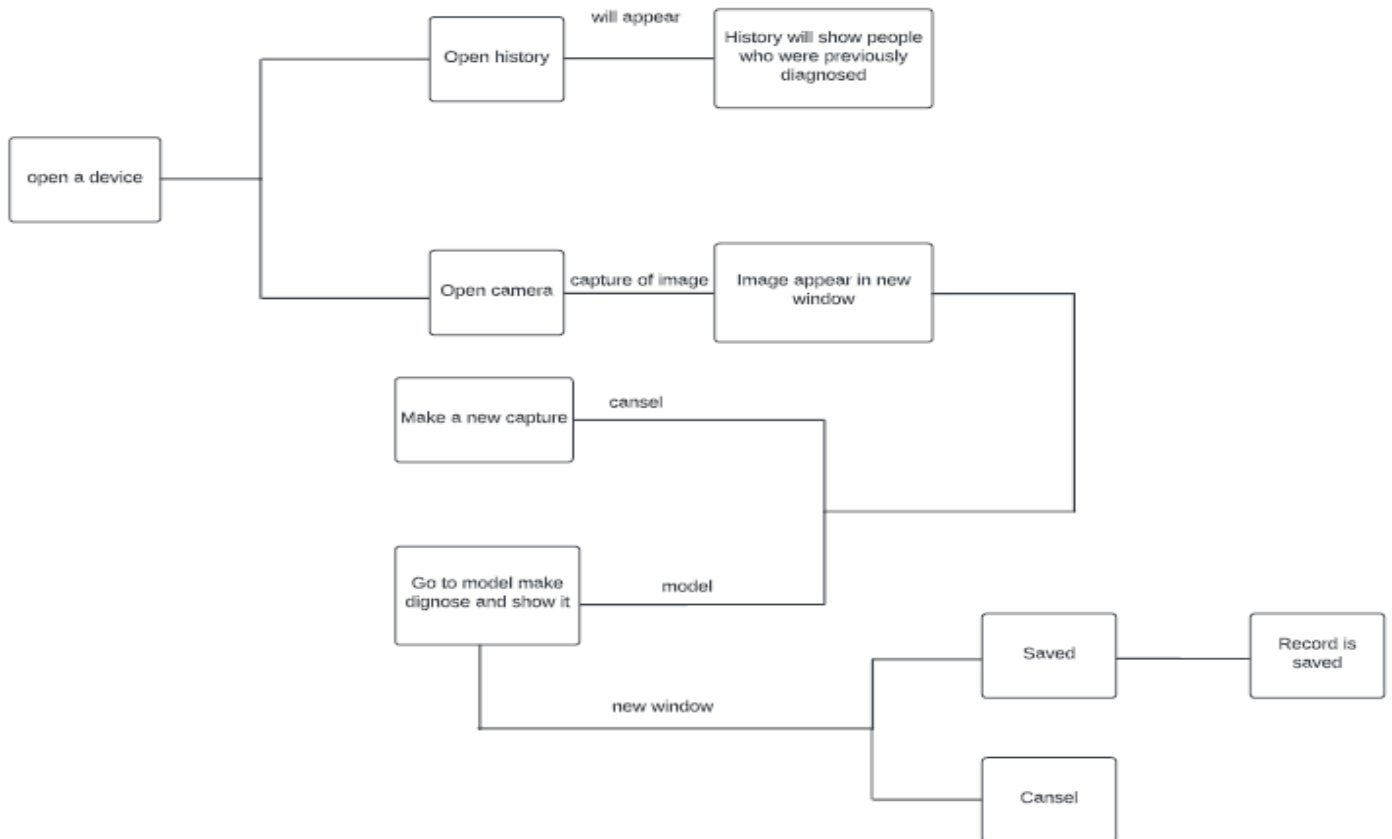

The data is received from a user side raspberry pi 4 via raspberry pi camera that takes a picture that is recorded & displayed at lcd that convert the data to the API and start to detect and extract the picture feature and it back into raspberry pi to display the result at touch lcd.

3.3 System Requirements

- (Clients, Customers and Users)
- **Developers:** Responsible for development and porting software firmware and future modifications to the system

- **dermatologists:** Uses the data from the real time imagery to classify images of skin cancer on par with dermatologists to know if it is disease malignant or not
- **Students/Scientists/Researchers:** Analyses the data as well as the Video stream for scientific evaluation and be responsible to extracting the features
- **Patients:** could enable lifesaving and fast diagnoses, even outside the hospital.

4 Chapter Four (Analysis and Design)

In this chapter we will present a detailed description of our system features, core components, system and customer requirements. It will give in-depth information about our system that we will develop.

4.1 Customer Requirements

Customer requirements are related to what is needed from a customer or client perspective.

4.1.1 Main Features:

- Predict skin cancer is benign or malignant.
- The Skin image and prediction printed on the device screen.

4.1.2 Key Elements:

- The device consists of camera, camera arm, screen and the processing unit.
- The device has four buttons (power button, camera arm button, button for tacking a picture, button for specifying that the image is good to run the model on it.
- The device has a power management system.
- The doctor uses the device camera to take the skin image.
- The image is processed and classified.
- The image and its prediction presented on the screen.

4.1.3 Mapping table:

| Requirement ID | Covers |
|----------------------|---|
| Covers_M_CR_001_V1.0 | The device opened using the power button. |
| Covers_M_CR_002_V1.0 | Use touch screen to increase usability. |
| Covers_M_CR_003_V1.1 | The doctor chooses whether choose open camera for take an image or redirect to previous records. |
| Covers_M_CR_004_V1.0 | A real time camera will show on the screen. |
| Covers_M_CR_005_V1.1 | The doctor can take a skin picture using capture button and specify if the image is good to send it to the model or he will take another one. |
| Covers_M_CR_006_V1.0 | The image sends to the ML model to detect whether it is benign or malignant. |
| Covers_M_CR_007_V1.0 | The image and its prediction presented on the screen. |
| Covers_M_CR_008_V1.1 | The prediction will be saved into the records. |
| Covers_M_CR_009_V1.1 | The doctor can list the previous records for all patients with their classifications. |

4.2 Customer Requirements Specifications

The customer requirements specifications describe the customer requirements in detailed steps to explain expected cases or scenarios.

4.2.1 System Context: -

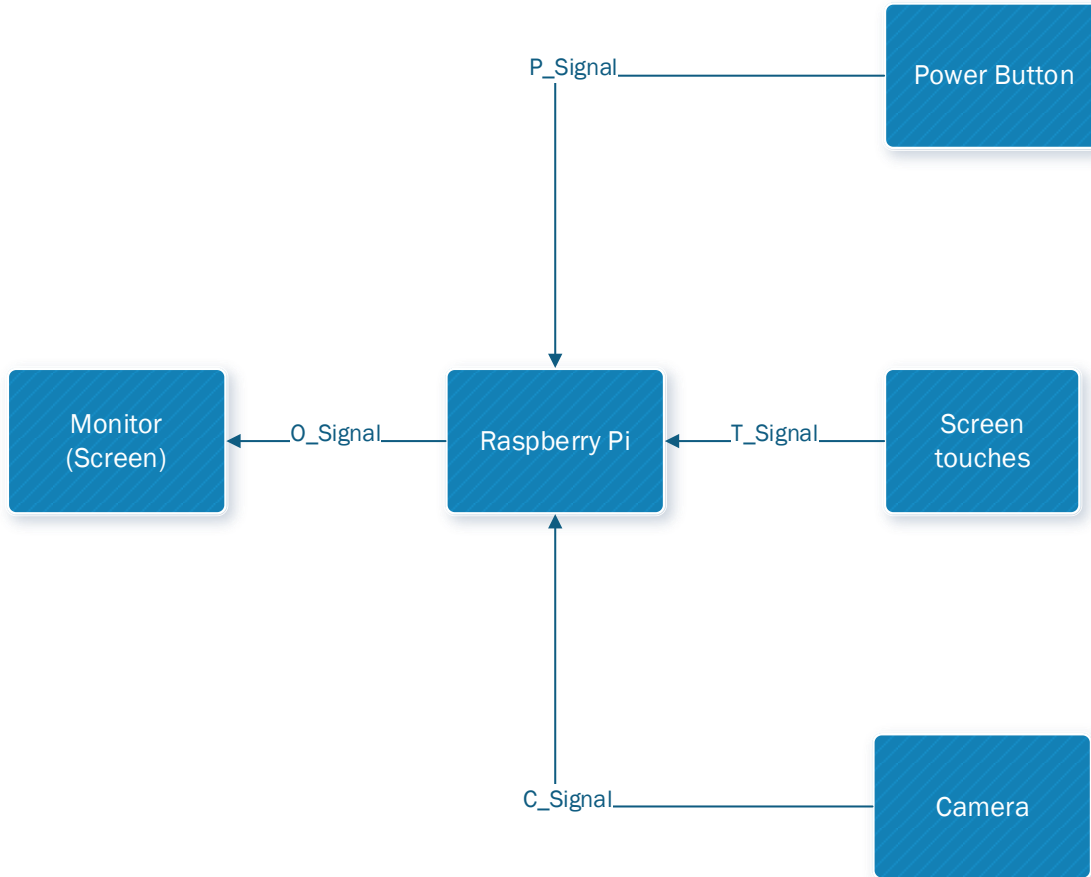


Figure 4.1: Context System Diagram

4.2.2 CRS requirements:

| | | | |
|-------------|--|--------|----------------------|
| Req_ID | Req_M_CRS_001- v1.1 | Covers | Covers_M_CR_001_V1.0 |
| Description | Connect device with power supplier outputs 5v and 2A | | |

| | | | |
|-------------|--|--------|----------------------|
| Req_id | Req_M_CRS_002- v1.1 | Covers | Covers_M_CR_001_V1.0 |
| Description | Press the power button to open the device. | | |

| | | | |
|-------------|---|--------|----------------------|
| Req_id | Req_M_CRS_003- v1.1 | Covers | Covers_M_CR_002_V1.1 |
| Description | The doctor talk action by pressing on a touch screen outputting a T_Signal. | | |

| | | | |
|-------------|--|--------|----------------------|
| Req_id | Req_M_CRS_004- v1.1 | Covers | Covers_M_CR_003_V1.1 |
| Description | the device opens on a screen with two options: <ul style="list-style-type: none"> • Open camera to take skin cancer picture from a patient. • Open previews records. | | |

| | | | |
|-------------|---|--------|----------------------|
| Req_id | Req_M_CRS_005- v1.1 | Covers | Covers_M_CR_004_V1.0 |
| Description | If the doctor chooses to “open camera” button Raspberry PI will take T_Signal and detect the action (open the camera) then takes C_Signals and show a real time view on the screen. | | |

| | | | |
|-------------|---|--------|----------------------|
| Req_id | Req_M_CRS_006- v1.1 | Covers | Covers_M_CR_005_V1.1 |
| Description | when the doctor press on capture button the Raspberry PI will takes T_Signal to detect the action (taking a picture) and capture the image. | | |

| | | | |
|-------------|---|--------|----------------------|
| Req_id | Req_M_CRS_007- v1.1 | Covers | Covers_M_CR_005_V1.1 |
| Description | The captured image will show on the screen with two options: <ul style="list-style-type: none"> • Retake another picture. • Send picture the model to detect whether it is benign or malignant. | | |

| | | | |
|-------------|--|--------|----------------------|
| Req_id | Req_M_CRS_008- v1.1 | Covers | Covers_M_CR_005_V1.1 |
| Description | If the doctor press on retake button the Raspberry PI will takes “ T_Signal ” to detect the action (back to the Previous window). | | |

| | | | |
|-------------|--|--------|----------------------|
| Req_id | Req_M_CRS_009- v1.1 | Covers | Covers_M_CR_005_V1.1 |
| Description | The Raspberry PI takes “ C_Signal ” and convert it into Frames. | | |

| | | | |
|-------------|---|--------|----------------------|
| Req_id | Req_M_CRS_010- v1.1 | Covers | Covers_M_CR_005_V1.1 |
| Description | The Raspberry PI display these Frames into camera window on the Screen. | | |

| | | | |
|-------------|---|--------|----------------------|
| Req_id | Req_M_CRS_011- v1.1 | Covers | Covers_M_CR_005_V1.1 |
| Description | If the doctor press on send button the Raspberry PI takes T_Signal to detect action (the image is good to send to the model). | | |

| | | | |
|-------------|--|--------|----------------------|
| Req_id | Req_M_CRS_0012- v1.1 | Covers | Covers_M_CR_006_V1.0 |
| Description | The Raspberry PI will send image to the model using API Post method. | | |

| | | | |
|-------------|---|--------|----------------------|
| Req_id | Req_M_CRS_013- v1.1 | Covers | Covers_M_CR_006_V1.0 |
| Description | the image will be resized and standardized. | | |

| | | | |
|-------------|---|--------|----------------------|
| Req_id | Req_M_CRS_014- v1.1 | Covers | Covers_M_CR_006_V1.0 |
| Description | After standardizing and resizing the image, it will be sent to the model and the prediction produced. | | |

| | | | |
|-------------|--|--------|----------------------|
| Req_id | Req_M_CRS_015- v1.1 | Covers | Covers_M_CR_006_V1.0 |
| Description | The API return prediction to the Raspberry PI. | | |

| | | | |
|-------------|---|--------|----------------------|
| Req_id | Req_M_CRS_016- v1.1 | Covers | Covers_M_CR_007_V1.0 |
| Description | The Raspberry PI present the image with its prediction on the screen. | | |

| | | | |
|-------------|--|--------|-----------------------|
| Req_id | Req_M_CRS_017- v1.1 | Covers | Covers_SS_CR_008_V1.1 |
| Description | After displaying the image, the doctor will have two options: <ul style="list-style-type: none"> • Save the image. • Cancel. | | |

| | | | |
|-------------|--|--------|----------------------|
| Req_id | Req_M_CRS_018- v1.1 | Covers | Covers_M_CR_008_V1.1 |
| Description | If the doctor press on Save button the Raspberry PI takes the T_Signal to detect the action (save the image), and save diagnosis to the records. | | |

| | | | |
|-------------|---|--------|----------------------|
| Req_id | Req_M_CRS_019- v1.1 | Covers | Covers_M_CR_008_V1.1 |
| Description | After saving window is displayed, doctor will has two options: <ul style="list-style-type: none"> • Add new record. • Save to existing records. | | |

| | | | |
|--------|---------------------|--------|----------------------|
| Req_id | Req_M_CRS_020- v1.1 | Covers | Covers_M_CR_008_V1.1 |
|--------|---------------------|--------|----------------------|

| | | | |
|-------------|---|--|--|
| Description | When the doctor press on “Add new record” the Raspberry PI takes the T_Signal to detect the action (displaying text box and create record button) | | |
|-------------|---|--|--|

| | | | |
|-------------|--|--------|----------------------|
| Req_id | Req_M_CRS_021- v1.1 | Covers | Covers_M_CR_008_V1.1 |
| Description | The doctor enters patient name on the text box and then press “Save” button to save new record with this patient name. | | |

| | | | |
|-------------|--|--------|----------------------|
| Req_id | Req_M_CRS_022- v1.1 | Covers | Covers_M_CR_008_V1.1 |
| Description | When the doctor press on “save to existing record” button a list of existing record names will be displayed with searching option. | | |

| | | | |
|-------------|---|--------|----------------------|
| Req_id | Req_M_CRS_023- v1.1 | Covers | Covers_M_CR_008_V1.1 |
| Description | After the doctor chooses a specific record the image with its prediction will be saved to the chosen record by pressing “Save” button | | |

| | | | |
|-------------|---|--------|----------------------|
| Req_id | Req_M_CRS_024- v1.1 | Covers | Covers_M_CR_009_V1.1 |
| Description | When the doctor press on “open previous records” button, Raspberry PI will takes T_Signal and detect the action (open the records). | | |

| | | | |
|-------------|---|--------|----------------------|
| Req_id | Req_M_CRS_025- v1.1 | Covers | Covers_M_CR_009_V1.1 |
| Description | The doctor can search for specific patient record using his name. | | |

4.3 Hardware Software Interface

4.3.1 Hardware specifications:

This project contains the following hardware components:

4.3.1.1 Raspberry pi 4 Model B

- The Raspberry Pi 4 Model B (Pi4B) is the first of a new generation of Raspberry Pi computers supporting more RAM and with significantly enhanced CPU, GPU and I/O performance; all within a similar form factor, power envelope and cost as

the previous generation Raspberry Pi 3B+. The Pi4B is available with either 1, 2 and 4 Gigabytes of LPDDR4 SDRAM.

- Quad core 64-bit ARM-Cortex A72 running at 1.5GHz
- 1, 2 and 4 Gigabyte LPDDR4 RAM options
- H.265 (HEVC) hardware decode (up to 4Kp60)
- H.264 hardware decode (up to 1080p60)
- VideoCore VI 3D Graphics
- Supports dual HDMI display output up to 4Kp60
- 802.11 b/g/n/ac Wireless LAN
- Bluetooth 5.0 with BLE
- 1x SD Card
- 2x micro-HDMI ports supporting dual displays up to 4Kp60 resolution
- 2x USB2 ports
- 2x USB3 ports
- 1x Gigabit Ethernet port (supports PoE with add-on PoE HAT)
- 1x Raspberry Pi camera port (2-lane MIPI CSI)
- 1x Raspberry Pi display port (2-lane MIPI DSI)
- 28x user GPIO supporting various interface options:
 - Up to 6x UART
 - Up to 6x I2C – Up to 5x SPI
 - 1x SDIO interface
 - 1x DPI (Parallel RGB Display)
 - 1x PCM – Up to 2x PWM channels
 - Up to 3x GPCLK outputs
- **Electrical Specification:** Stresses above those listed in Table 1 may cause permanent damage to the device. This is a stress rating only; functional operation of the device under these or any other conditions above those listed in the operational sections of this specification is not implied. Exposure to absolute maximum rating conditions for extended periods may affect device reliability.

| Symbol | Parameter | Minimum | Maximum | Unit |
|--------|------------------|---------|---------|------|
| VIN | 5V Input Voltage | -0.5 | 6.0 | V |

Table 1: Absolute maximum

4.3.1.2 Raspberry Pi Camera Module

- **Description:** High-Definition camera module compatible with all Raspberry Pi models. Provides high sensitivity, low crosstalk and low noise image capture in an ultra-small and lightweight design. The camera module connects to the Raspberry Pi board via the CSI connector designed specifically for interfacing to cameras. The

CSI bus is capable of extremely high data rates, and it exclusively carries pixel data to the processor.

- **Image Sensor:** Sony IMX 219 PQ CMOS image sensor in a fixed-focus module.
- **Resolution:** 8-megapixel
- Still picture resolution: 3280 x 2464
- **Max image transfer rate:** 1080p: 30fps (encode and decode) 720p: 60fps
- **Connection to Raspberry Pi:** 15-pin ribbon cable, to the dedicated 15-pin MIPI Camera Serial Interface (CSI-2).
- Image control functions:
 - Automatic exposure control
 - Automatic white balance
 - Automatic band filter
 - Automatic 50/60 Hz luminance detection
 - Automatic black level calibration
- Temp range:
 - Operating: -20° to 60°
 - Stable image: -20° to 60°

4.3.1.3 7-inch Touchscreen Display B

- 7 inch Capacitive Touch Screen LCD, HDMI interface, with a resolution of 800×480 and a capacitive touch panel, which supports Raspberry Pi and can also be used as a computer monitor.
- 800×480 hardware resolution.
- 5-points capacitive touch control
- When used with Raspberry Pi, supports Raspberry Pi OS / Ubuntu / Kali and Retropie
- Support backlight control, more power saving.

4.3.1.4 5V Cooling Fan

- Operating voltage: 5V
- Current: 0.2 A
- Brushless DC fan

4.3.1.5 Button

- Mode of Operation: Tactile feedback
- Power Rating: MAX 50mA 24V DC
- Insulation Resistance: 100Mohm at 100v
- Operating Force: 2.55±0.69 N
- Contact Resistance: MAX 100mOhm
- Operating Temperature Range: -20 to +70 °C
- Storage Temperature Range: -20 to +70 °C

4.3.2 HSI System Context:

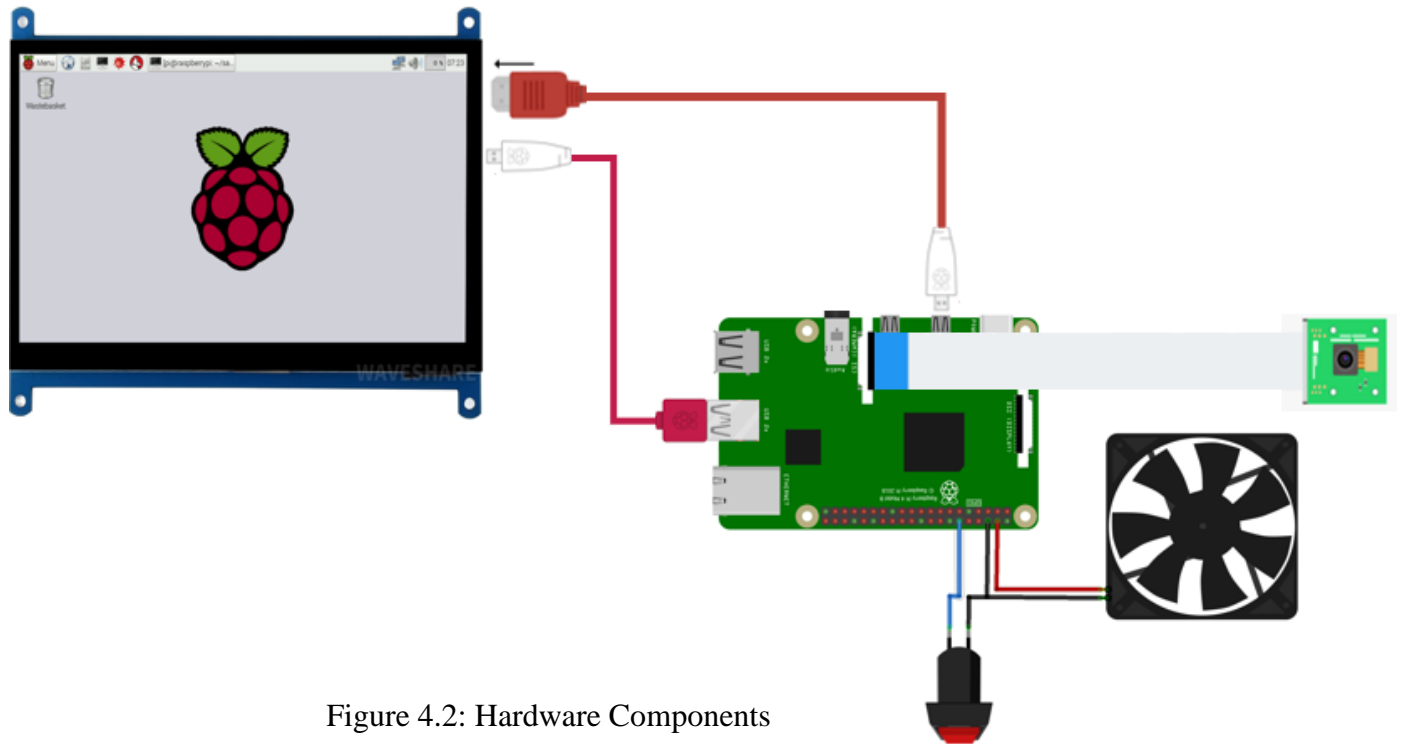


Figure 4.2: Hardware Components

4.3.3 HW and SW Restrictions:

4.3.3.1 Raspberry Pi 4 Model B:

- Power Requirements the Pi4B requires a good quality USB-C power supply capable of delivering 5V at 3A. If attached downstream USB devices consume less than 500mA, a 5V, 2.5A supply may be used.

4.3.3.2 Raspberry Pi Camera Module and 7 inch Touchscreen Display:

To avoid malfunction or damage to this product, please observe the following:

- Before connecting the device, shut down your Raspberry Pi computer and disconnect it from external power.
- If the cable becomes detached, pull the locking mechanism forward on the connector, insert the ribbon cable ensuring the metal contacts face towards the circuit board, then push the locking mechanism back into place.
- This device should be operated in a dry environment at 0–50°C.

- Do not expose it to water or moisture, or place on a conductive surface whilst in operation.
- Do not expose it to excessive heat from any source.
- Care should be taken not to fold or strain the ribbon cable.
- Care should be taken when screwing in parts or fitting a tripod. A cross-thread can cause irreparable damage and void the warranty.
- Take care whilst handling to avoid mechanical or electrical damage to the printed circuit board and connectors.
- Avoid handling the printed circuit board whilst it is powered and only handle by the edges to minimize the risk of electrostatic discharge damage.
- Store in a cool, dry location.
- Avoid rapid changes of temperature, which can cause moisture build up in the device, affecting image quality

4.3.4 HSI Functional Requirements:

| | | | |
|--------------------|--|---------------|----------------------|
| Req_ID | Req_ HSI_M_001-V1.0 | Covers | Covers_M_CR_001_V1.0 |
| Description | Connect Raspberry pi to source power and wait for the program to run | | |

| | | | |
|--------------------|---|---------------|--|
| Req_ID | Req_ HSI_M_002-V1.0 | Covers | Covers_M_CR_005_V1.1 Covers_M_CR_004_V1.0 |
| Description | using Home Button in app to open capture window for taking images | | |

| | | | |
|--------------------|---|---------------|----------------------|
| Req_ID | Req_ HSI_M_003-V1.0 | Covers | Covers_M_CR_009_V1.1 |
| Description | using History button in app to open Previews images | | |

| | | | |
|--------------------|--|---------------|----------------------|
| Req_ID | Req_ HSI_M_004-V1.0 | Covers | Covers_M_CR_006_V1.0 |
| Description | In Home window, use Send image to Machine learning model or recapture a good one | | |

| | | | |
|--------------------|--|---------------|----------------------|
| Req_ID | Req_ HSI_M_005-V1.0 | Covers | Covers_M_CR_008_V1.1 |
| Description | Using Save button to save result of diagnosis as image from Machine learning model | | |

| | | | |
|--------------------|--|---------------|--|
| Req_ID | Req_ HSI_M_006-V1.0 | Covers | Covers_M_CR_002_V1.0 Covers_M_CR_004_V1.0 |
| Description | Display the Result or captures image on the Screen with touch usage. | | |

| | | | |
|--------------------|---|---------------|--|
| Req_ID | Req_ HSI_M_007-V1.0 | Covers | |
| Description | Shut down button to turn off the raspberry pi | | |

4.4 System Requirements Specifications

4.4.1 Software context:

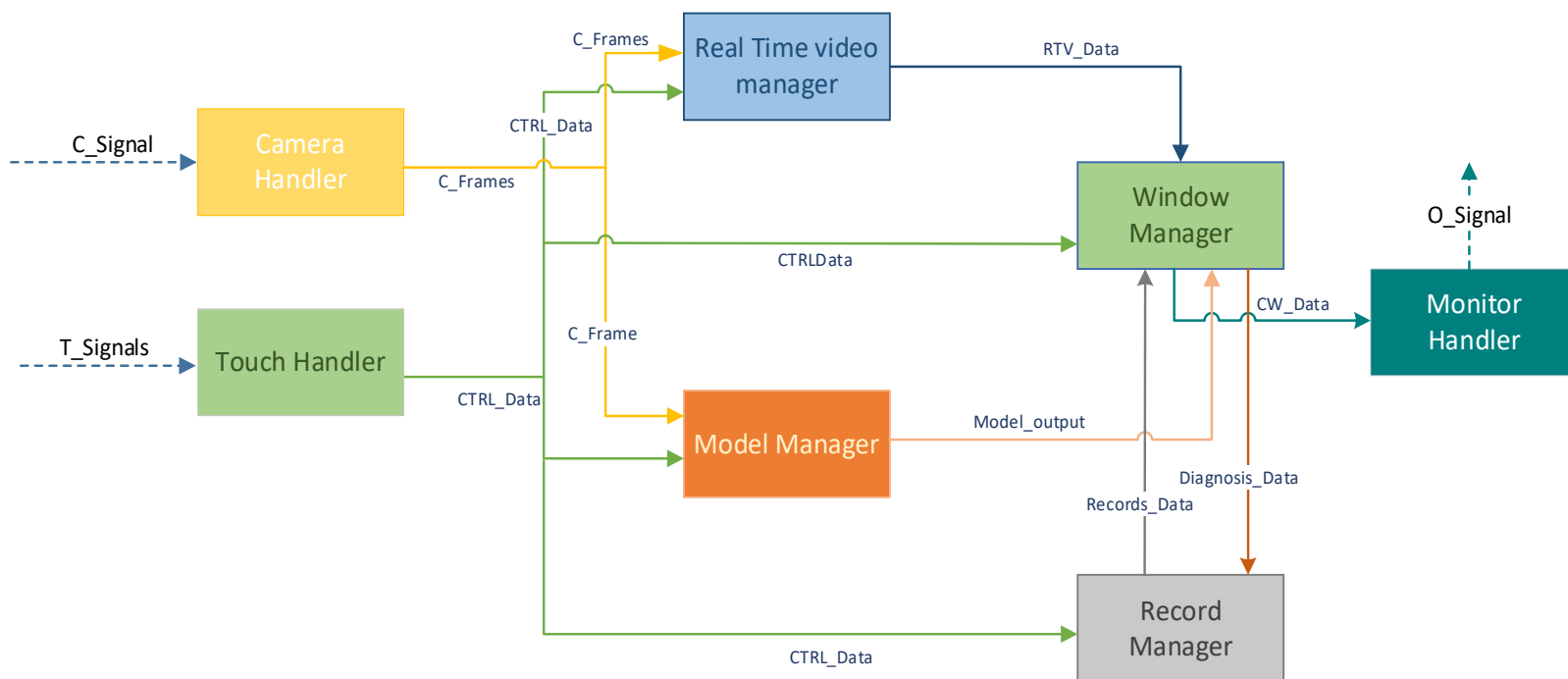


Figure 4.3: System flow diagram

4.4.2 SRS Requirements

Touch Handler

| | | | |
|--------------------|--|----------------|---------------------|
| Req_ID | Req _M_SRS_Touches_Handler_01-V1.0 | Covers | Req_M_CRS_003- v1.1 |
| Description | In Touch handler, SW shall read the touch signal “T_Signal “ | | |
| Inputs | T_Signal | Outputs | ----- |

| | | | |
|--------------------|---|----------------|---------------------|
| Req_ID | Req _M_SRS_Touches_Handler_02-V1.0 | Covers | Req_M_CRS_003- v1.1 |
| Description | The doctor press on the touch screen “T_Signal” and the touch handler specify action in “CTRL_Data” | | |
| Inputs | T_Signal | Outputs | CTRL_Data |

Camera Handler

| | | | |
|--------------------|--|----------------|---------------------|
| Req_ID | Req _M_SRS_Camera_Handler_02-V1.0 | Covers | Req_M_CRS_009- v1.1 |
| Description | In Camera handler, SW shall read the camera signal “C_Signal “ | | |
| Inputs | C_Signal | Outputs | ----- |

| | | | |
|--------------------|--|----------------|---------------------|
| Req_ID | Req _M_SRS_Camera_Handler_02-V1.0 | Covers | Req_M_CRS_009- v1.1 |
| Description | The “C_Signal” will convert into Frames. | | |
| Inputs | C_Signal | Outputs | C_Frames |

Real Time Video Manager

| | | | |
|--------------------|---|----------------|---------------------|
| Req_ID | Req _M_SRS_RTV_Handler_01-V1.0 | Covers | Req_M_CRS_011- v1.1 |
| Description | The model manager shall read “CTRL_Data” and “C_Frames” | | |
| Inputs | C_Frames, CTRL_Data | Outputs | ---- |

| | | | |
|--------------------|---|----------------|---------------------|
| Req_ID | Req _M_SRS_RTV_Handler_02-V1.0 | Covers | Req_M_CRS_007- v1.1 |
| Description | The real time video manager shall detect the action that “CTRL_Data” Specifies. | | |
| Inputs | C_Frames, CTRL_Data | Outputs | ---- |

| | | | |
|--------------------|--|----------------|---------------------|
| Req_ID | Req _M_SRS_RTV_Handler_03-V1.0 | Covers | Req_M_CRS_010- v1.1 |
| Description | IF “CTRL_Data” Specifies that open camera is the next action, the camera frames will be presented on the screen. | | |
| Inputs | C_Frames, CTRL_Data | Outputs | RTV_Data |

Model Manager

| | | | |
|--------------------|--|----------------|---------------------|
| Req_ID | Req _M_SRS_Camera_Handler_01-V1.0 | Covers | Req_M_CRS_011- v1.1 |
| Description | The model manager shall read “CTRL_Data” and “C_Frame” | | |
| Inputs | C_Frame, CTRL_Data | Outputs | ---- |

| | | | |
|--------------------|---|----------------|---------------------|
| Req_ID | Req _M_SRS_Camera_Handler_02-V1.0 | Covers | Req_M_CRS_011- v1.1 |
| Description | The model manager shall detect the action that “CTRL_Data” Specifies. | | |
| Inputs | C_Frame, CTRL_Data | Outputs | ---- |

| | | | |
|--------------------|--|----------------|---------------------|
| Req_ID | Req _M_SRS_Camera_Handler_03-V1.0 | Covers | Req_M_CRS_012- v1.1 |
| Description | IF “CTRL_Data” Specifies that send image to the model is the next action, the captured camera frame will be sent to the API and Present the output throw window manager. | | |
| Inputs | C_Frame, CTRL_Data | Outputs | ---- |

| | | | |
|--------------------|---|----------------|---------------------|
| Req_ID | Req_M_SRS_Camera_Handler_04-V1.0 | Covers | Req_M_CRS_012- v1.1 |
| Description | Send the “ model_output ” as output to window manager. | | |
| Inputs | C_Frame, CTRL_Data | Outputs | Model_output |

Record Manager

| | | | |
|--------------------|---|----------------|---------------------|
| Req_ID | Req_M_SRS_Record_Handler_01-V1.0 | Covers | Req_M_CRS_018- v1.1 |
| Description | The record manager shall read “ CTRL_Data ” and “ Diagnosis_Data ”. | | |
| Inputs | Diagnosis_Data, CTRL_Data | Outputs | ---- |

| | | | |
|--------------------|---|----------------|---------------------|
| Req_ID | Req_M_SRS_Record_Handler_02_V1.0 | Covers | Req_M_CRS_018- v1.1 |
| Description | The record manager shall detect the action that “ CTRL_Data ” Specifies. | | |
| Inputs | Diagnosis_Data, CTRL_Data | Outputs | ---- |

| | | | |
|--------------------|---|----------------|---------------------|
| Req_ID | Req_M_SRS_Record_Handler_03_V1.0 | Covers | Req_M_CRS_018- v1.1 |
| Description | IF “ CTRL_Data ” Specifies that save diagnosis is the next action, the captured image with its diagnosis will be saved into records. | | |
| Inputs | Diagnosis_Data, CTRL_Data | Outputs | ---- |

| | | | |
|--------------------|---|----------------|---------------------|
| Req_ID | Req_M_SRS_Record_Handler_04_V1.0 | Covers | Req_M_CRS_022- v1.1 |
| Description | IF “ CTRL_Data ” Specifies that record history is the next action, the saved records will be presented throw window manager. | | |
| Inputs | Diagnosis_Data, CTRL_Data | Outputs | Records_Data |

Window Manager

| | | | |
|---------------|----------------------------------|---------------|---------------------|
| Req_ID | Req_M_SRS_Window_Handler_01-V1.0 | Covers | Req_M_CRS_018- v1.1 |
|---------------|----------------------------------|---------------|---------------------|

| | | | |
|--------------------|---|----------------|------|
| Description | The window manager shall read “CTRL_Data”, “RTV_Data”, “Model_output” and “Record_Data” | | |
| Inputs | CTRL_Data, RTV_Data, Model_output, Record_Data | Outputs | ---- |

| | | | |
|--------------------|--|----------------|----------------------|
| Req_ID | Req_M_SRS_Window_Handler_02_V1.0 | Covers | Covers_M_CR_003_V1.1 |
| Description | The window manager shall detect the action that “CTRL_Data” Specifies. | | |
| Inputs | CTRL_Data, RTV_Data, Model_output, Record_Data | Outputs | ---- |

| | | | |
|--------------------|---|----------------|----------------------|
| Req_ID | Req_M_SRS_Camera_Handler_03_V1.0 | Covers | Covers_M_CR_010_V1.1 |
| Description | IF “CTRL_Data” Specifies that open camera is the next action, the window manager present “RTV_Data” into camera window. | | |
| Inputs | CTRL_Data, RTV_Data, Model_output, Record_Data | Outputs | CW_Data |

| | | | |
|--------------------|---|----------------|----------------------|
| Req_ID | Req_M_SRS_Camera_Handler_04_V1.0 | Covers | Covers_M_CR_024_V1.1 |
| Description | IF “CTRL_Data” Specifies that open records action is the next action, the window manager present “Records_Data” into record history window. | | |
| Inputs | CTRL_Data, RTV_Data, Model_output, Records_Data | Outputs | CW_Data |

| | | | |
|--------------------|--|----------------|----------------------|
| Req_ID | Req_M_SRS_Camera_Handler_05_V1.0 | Covers | Covers_M_CR_016_V1.1 |
| Description | IF “CTRL_Data” Specifies that send image to the model action is the next action, the window manager present “Model_output” diagnosis window. | | |
| Inputs | CTRL_Data, RTV_Data, Model_output, Records_Data | Outputs | CW_Data |

| | | | |
|---------------|----------------------------------|---------------|----------------------|
| Req_ID | Req_M_SRS_Camera_Handler_06_V1.0 | Covers | Covers_M_CR_018_V1.1 |
|---------------|----------------------------------|---------------|----------------------|

| | | | |
|--------------------|---|----------------|----------------|
| Description | IF “CTRL_Data” Specifies that save diagnosis action is the next action, the window manager present the image with its diagnosis to the records. | | |
| Inputs | CTRL_Data, RTV_Data, Model_output, Record_Data | Outputs | Diagnosis_Data |

5 Chapter Five (Implementation)

5.1 Deep learning model

This section describes in detail the model implementation. The model has been written in python and developed in Kaggle, google collaborator (colab) in some experiments and features testing and finally MS Azure trying to get more computation power.

We use the convolution neural network to build our model and train this model on the official datasets of the SIIM-ISIC Melanoma Classification over multiple Challenges and use TensorFlow as a framework to implement the model.

Some techniques are used like making class weights balanced instead of unbalanced, transfer learning, fine tuning, learning rate techniques (constant LR, schedule LR or CLRs), standardize or normalize images or some regularization techniques. Execute a lot of experiments by changing all the hyperparameter like the mentioned techniques and other hyperparameters.

5.1.1 Why CNN

A convolutional neural network is very good at classifying images. This is one of the widely known and well-publicized facts, but why is it so, and why is it so in dealing with skin cancer images?

5.1.2 CNN Features

The number of parameters in a neural network grows rapidly with the increase in the number of layers. This can make training for a model computationally heavy (and sometimes not feasible). Tuning so many of parameters can be a very huge task. The time taken for tuning these parameters is diminished by CNNs.

CNNs are fully connected feed forward neural networks. CNNs are very effective in reducing the number of parameters without losing on the quality of models. Images

have high dimensionality (as each pixel is considered as a feature) which suits the above-described abilities of CNNs. CNNs are trained to identify the edges of objects in any image.

All the layers of a CNN have multiple convolutional filters working and scanning the complete feature matrix and carry out the dimensionality reduction. This enables CNN to be a very apt and fit network for image classifications and processing.

Skin Cancer Images require:

When work in the problem of skin cancer you must:

- Make the huge number of parameters under control this is very useful when dealing with skin cancer images as you must have control on every feature inside the image.
- Treat each pixel as a feature so, when reducing the number of parameters, it must be ensured that this process occurs without losing the quality of the model to reveal that the image is benign or malignant.
- Take caution when dealing with edges as edges are very important factor in classifying skin cancer images (edges shape, texture, density, ...etc.). so, identifying and detecting images is very important.

From these points and the characteristics of CNN identified above, we have concluded that CNN is the best solution in our case.

5.1.3 Why TensorFlow

Whether you're an expert or a beginner, TensorFlow is an end-to-end platform that makes it easy for you to build and deploy ML models. This discussed in the next three points:

- **Easy model building**

TensorFlow offers multiple levels of abstraction so you can choose the right one for your needs. Build and train models by using the high-level Keras API, which makes getting started with TensorFlow and machine learning easy.

If you need more flexibility, eager execution allows for immediate iteration and intuitive debugging. For large ML training tasks, use the Distribution Strategy API for distributed training on different hardware configurations without changing the model definition.

- **Robust ML production anywhere**

TensorFlow has always provided a direct path to production. Whether it's on servers, edge devices, or the web, TensorFlow lets you train and deploy your model easily, no matter what language or platform you use.

Use TensorFlow Extended (TFX) if you need a full production ML pipeline. For running inference on mobile and edge devices, use TensorFlow Lite. Train and deploy models in JavaScript environments using TensorFlow.js.

- **Powerful experimentation for research**

Build and train state-of-the-art models without sacrificing speed or performance. TensorFlow gives you the flexibility and control with features like the Keras Functional API and Model Subclassing API for creation of complex topologies. For easy prototyping and fast debugging, use eager execution.

TensorFlow also supports an ecosystem of powerful add-on libraries and models to experiment with, including Ragged Tensors, TensorFlow Probability, Tensor2Tensor and BERT.

5.1.4 About used dataset

The dataset is combination of the following:

- **ISIC 2020 challenge dataset**

The dataset contains 33,126 dermoscopic training images of unique benign and malignant skin lesions from over 2,000 patients. Each image is associated with one of these individuals using a unique patient identifier. All malignant diagnoses have been confirmed via histopathology, and benign diagnoses have been confirmed using either expert agreement, longitudinal follow-up, or histopathology. A thorough publication describing all features of this

dataset is available in the form of a pre-print that has not yet undergone peer review.

The dataset was generated by the International Skin Imaging Collaboration (ISIC) and images are from the following sources: Hospital Clínic de Barcelona, Medical University of Vienna, Memorial Sloan Kettering Cancer Center, Melanoma Institute Australia, University of Queensland, and the University of Athens Medical School.

- **ISIC 2019 challenge dataset**

25,331 images are available for training across 8 different categories (Melanoma, Melanocytic nevus, Basal cell carcinoma, Actinic keratosis, Benign keratosis (solar lentigo / seborrheic keratosis / lichen planus-like keratosis), Dermatofibroma, Vascular lesion, Squamous cell carcinoma, None of the others). Additionally, the test dataset contains an additional outlier class not represented in the training data, which developed systems must be able to identify.

5.1.5 Balancing class weights

As we mentioned above our dataset consist of all categories (types) of skin cancer and all of them divided into malignant (melanoma) which represent 4922 image (8.60% of total) and benign (other types) which represent 52,302 image that causes imbalance class weights and we will discuss it in three parts what, why and how.

- **What is imbalance Class Weights?**

Class imbalance is a problem that occurs in machine learning classification problems. It merely tells that the target class's frequency is highly imbalanced, i.e., the occurrence of one of the classes is very high compared to the other classes present. In other words, there is a bias or skewness towards the majority class present in the target. Suppose we consider a binary classification where the majority target class has 10000 rows, and the minority target class has only 100 rows. In that case, the ratio is 100:1, i.e., for every 100-majority class, there is only one minority class present. This problem is what we refer to as class imbalance. Some of the general areas where we can

find such data are fraud detection, churn prediction, medical diagnosis, e-mail classification, etc.

We will be working on our dataset to understand class imbalance properly. In the below image, you can see the number of images in each class in our target variable.

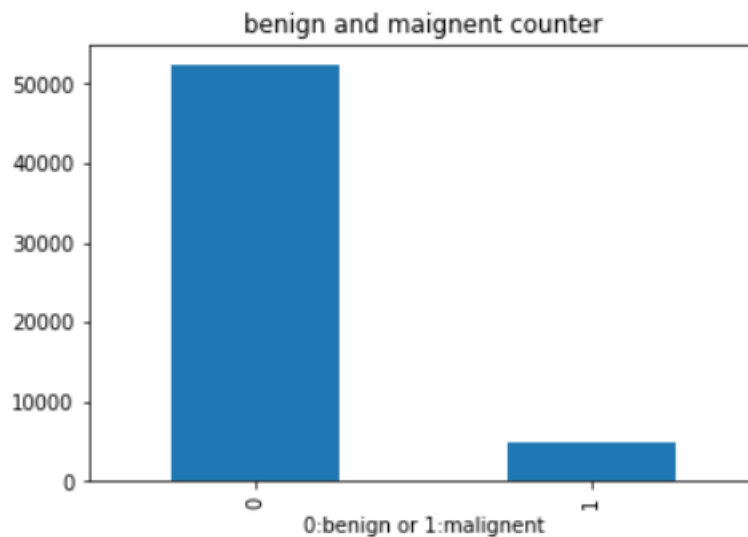


Figure 5.1: benign and malignant counter

Here,

0: signifies that the patient didn't have a melanoma.

1: signifies that the patient had a melanoma.

From the counter, we can see that there are only 4922 image (8.6 %) of patients who had a melanoma. So, this is a classic class imbalance problem.

- **Why is it essential to deal with class imbalance?**

So far, we got the intuition about class imbalance. But why is it necessary to overcome this, and what problems does it create while modeling with such data?

Most machine learning algorithms assume that the data is evenly distributed within classes. In the case of class imbalance problems, the extensive issue is that the algorithm will be more biased towards predicting the majority class (Malignant our case). The algorithm will not have enough data to learn the patterns present in the minority class (benign).

Let's take a real-life example to understand this better.

Consider you have shifted from your hometown to a new city, and you been living here for the past month. When it comes to your hometown, you will be very familiar with all the locations like your home, routes, essential shops, tourist spots, etc. because you had spent your whole childhood there. But when it comes to the new city, you would not have many ideas about where each location exactly is, and the chances of taking the wrong routes and getting lost will be very high. Here, your hometown is our majority class, and the new city is the minority class.

Similarly, this happens in class imbalance. The model has adequate information about the majority class but insufficient information about our minority class. That is why there will be high misclassification errors for the minority class.

When we train our model without solving this problem, we got 8% wrong prediction images which is so good but when we use AUC matrix we got 42% error rate this discussing the impact of imbalance class weights.

- **How to use class weights to solve this imbalance classes**

Most machine learning algorithms are not very useful with biased class data. But we can modify the current training algorithm to consider the skewed distribution of the classes. This can be achieved by giving different weights to both the majority and minority classes. The difference in weights will influence the classification of the classes during the training phase. The whole purpose is to penalize the misclassification made by the minority class by setting a higher-class weight and at the same time reducing weight for the majority class.

To make this a bit clear, we will be reviving the city example we considered earlier.

Please think of it this way that the last month you have spent in the new city, instead of going out when it is needed, you spent the whole month exploring the city. You spent more time understanding the city routes and places the entire month. Giving more time to research will help you to understand the new city better, and the chances of getting lost will reduce. And this is precisely how class weights work. During the training, we give more weightage to the minority class in the cost function of the algorithm so that it could provide a higher penalty to the minority class and the algorithm could focus on reducing the errors for the minority class.

There is a threshold to which you should increase and decrease the class weights for the minority and majority class respectively. If you give very high-class weights to the minority class, chances are the algorithm will get biased towards the minority class, and it will increase the errors in the majority class.

In our case we can increase the class weights of minority class (malignant) as if we tend to zero false negative this is very useful in our case even if false positive matrix increases a little. This comes from the logic that if a patient has a malignant tumor, we must diagnose it correctly, this is a higher priority than the diagnosis of a benign condition.

When using CNN in TensorFlow the fit method has an in-built parameter “class_weight” which helps us optimize the scoring for the minority class just the way we have learned so far.

By default, the value of class_weight=None, i.e., both the classes have been given equal weights. Other than that, we can either give it as ‘balanced’ or we can pass a dictionary that contains manual weights for both the classes.

To be more precise, the formula to calculate a balance class weight in our case is:

```

total_img = train["target"].size

malignant = np.count_nonzero(train["target"])
benign = total_img - malignant
print('Examples:\n    Total: {}\n    Positive: {} ({:.2f}% of total)\n'.format(
    total_img, malignant, 100 * malignant / total_img))

weight_for_0 = (1 / benign)*(total_img)/2.0
weight_for_1 = (1 / malignant)*(total_img)/2.0

class_weight = {0: weight_for_0, 1: weight_for_1}

print('Weight for class 0: {:.2f}'.format(weight_for_0))
print('Weight for class 1: {:.2f}'.format(weight_for_1))

```

Examples:

```

    Total: 57224
    Positive: 4922 (8.60% of total)

Weight for class 0: 0.55
Weight for class 1: 5.81

```

Figure 5.2: computing class weights

5.2 Image standardization

Standardization is one kind of scaling. We need to scale when the features are of different scales, units, ranges etc. But in image all the feature columns are nothing but the intensities. Hence, all of them are already scaled in the same range [0-255]. So, what's the necessity of re-scaling them? But we can still standardize the data and which is helpful in our case.

Consider how a neural network learns its weights. C(NN)s learn by continually adding gradient error vectors (multiplied by a learning rate) computed from backpropagation to various weight matrices throughout the network as training examples are passed through.

The thing to notice here is the "multiplied by a learning rate".

If we didn't scale our input training vectors, the ranges of our distributions of feature values would likely be different for each feature, and thus the learning rate would cause corrections in each dimension that would differ (proportionally speaking) from

one another. We might be overcompensating a correction in one weight dimension while undercompensating in another.

This is non-ideal as we might find ourselves in a oscillating (unable to center onto a better maxima in cost(weights) space) state or in a slow moving (traveling too slow to get to a better maxima) state.

It is of course possible to have a per-weight learning rate, but it's yet more hyperparameters to introduce into an already complicated network that we'd also have to optimize to find. Generally learning rates are scalars.

Thus, we try to normalize images before using them as input into NN (or any gradient based) algorithm.

5.3 Benefits of using pretrained model

Pre-Trained models most of them have a good architecture and this architecture are trained on a very powerful machines talking a lot of time with huge datasets, so this improves its weights to a very good stage.

So, in our case we have not large dataset with 57224 image and not having the ability to run our model in a powerful machine we use the online free machine, but this does not give you all what you need.

There are two techniques to use pretrained models with them this is discussed below:

- **Transfer learning**

Transfer learning means taking the relevant parts of a pre-trained machine learning model and applying it to a new but similar problem. This will usually be the core information for the model to function, with new aspects added to the model to solve a specific task. we have to identify which areas of the model are relevant to the new task, and which parts will need to be retrained. For example, a new model will keep the processes that allow the machine to identify objects but retrain the model to identify a different specific object.

A machine learning model which identifies (recognize) a certain subject within a set of images is a prime candidate for transfer learning. The bulk of the model which deals with how to select different subjects can be kept. The part of the algorithm

which highlights a specific subject to categorize is the element that will be retrained. In this case, there's no need to rebuild and retrain a machine learning algorithm from scratch.

In supervised machine learning, models are trained to complete specific tasks from labelled data during the development process. Input and desired output are clearly mapped and fed to the algorithm. The model can then apply the learned trends and pattern recognition to new data. Models developed in this way will be highly accurate when solving tasks in the same environment as its training data. It will become much less accurate if the conditions or environment changes in real-world application beyond the training data. The need for a new model based on new training data may be required, even if the tasks are similar.

Transfer learning is a technique to help solve this problem. As a concept, it works by transferring as much knowledge as possible from an existing model to a new model designed for a similar task. For example, transferring the more general aspects of a model which make up the main processes for completing a task. This could be the process behind how objects or images are being identified or categorized. Extra layers of more specific knowledge can then be added to the new model, allowing it to perform its task in new environments.

Benefits of transfer learning in our case:

The main benefits of transfer learning include the saving of resources and improved efficiency when training new models.

- ✓ **Saving training data**

A huge range of data is usually required to accurately train a machine learning algorithm. Labelled training data takes time, effort and expertise to create. Transfer learning cuts down on the training data required for new machine learning models, as most of the model is already previously trained.

- ✓ **Efficiently train multiple models**

Machine learning models designed to complete complex tasks can take a long time to properly train. Transfer learning means that you don't have to start from scratch each time a similar model is required. The resources and time put into training a machine learning algorithm can be shared across different

models. The whole training process is made more efficient by reusing elements of an algorithm and transferring the knowledge already held by a model.

- Fine-Tune

Specifically, fine-tuning is a process that takes a model that has already been trained for one given task and then tunes the model to make it perform a second similar task.

Building and validating our model can be a huge task in its own right, depending on what data we're training it on.

This is what makes the fine-tuning approach so attractive. If we can find a trained model that already does one task well, and that task is similar to ours in at least some remote way, then we can take advantage of everything the model has already learned and apply it to our specific task.

If we have a model that has already been trained to recognize (categorize) images and we want to fine-tune this model to our project, we can first import our original model.

Generally, in fine tune we may remove last layer that classify the output may in thousand class like image_net dataset so we add our own classifying layer

In some experiments, we may want to remove more than just the last single layer, and we may want to add more than just one layer. This will depend on how similar the task is for each of the models.

Layers at the end of our model may have learned features that are very specific to the original task, whereas layers at the start of the model usually learn more general features like edges, shapes, and textures.

After we've modified the structure of the existing model, we then want to freeze the layers in our new model that came from the original model.

By freezing, we mean that we don't want the weights for these layers to update whenever we train the model on our new data for our new task. We want to keep all

of these weights the same as they were after being trained on the original task. We only want the weights in our new or modified layers to be updating.

After we do this, all that's left is just to train the model on our new data. Again, during this training process, the weights from all the layers we kept from our original model will stay the same, and only the weights in our new layers will be updating.

5.4 Loss and accuracy matrix and reasons behind them

❖ Loss function

For loss we use binary cross entropy this according to a lot of literature work and papers we read after this we find the following benefits of it:

- cross entropy is equivalent to fitting the model using maximum likelihood estimation. This on the other hand can be interpreted as minimizing the dissimilarity between the empirical distribution of training data and the distribution induced by the model.
- why maximizing the likelihood. Maximum Likelihood estimators have nice asymptotical properties (they are the best estimators in terms of convergence rates), they are consistent and statistically efficient.
- Given the prevalence of gradient-based learning algorithms (especially in deep learning), the logarithm in the cross-entropy will undo any exponential behavior that is often given through popular activation/output units like the sigmoid/softmax. Thus, the logarithm will avoid that the gradient saturates for extreme values. Large gradients are important for gradient-descent algorithms to make sufficient progress in each iteration.

This make us use it in the first time, but we tried a lot of loss functions like Dice loss, Focal Loss...etc. and we found that binary is actually best for us.

❖ Accuracy Matrix

We use AUC according to evaluation in the Kaggle competition to test our results on it. We will discuss why we use it regardless the Kaggle competition.

We will start by explaining what is ROC-curve?

An ROC curve (receiver operating characteristic curve) is a graph showing the performance of a classification model at all classification thresholds. This curve plots two parameters:

- True Positive Rate
- False Positive Rate

True Positive Rate (TPR) is a synonym for recall and is therefore defined as follows:

$$TPR = \frac{TP}{TP + FN}$$

False Positive Rate (FPR) is defined as follows: $FPR = \frac{FP}{FP + TN}$

An ROC curve plots TPR vs. FPR at different classification thresholds. Lowering the classification threshold classifies more items as positive, thus increasing both False Positives and True Positives. The following figure shows a typical ROC curve.

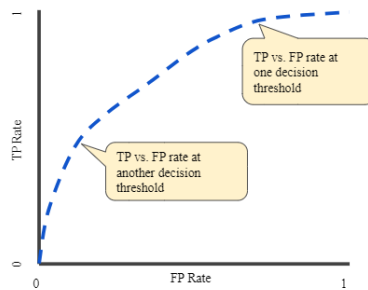


Figure 5.3: TP vs. FP rate at different classification thresholds

To compute the points in an ROC curve there's an efficient, sorting-based algorithm that can provide this information for us, called AUC.

AUC: Area Under the ROC Curve, That is, AUC measures the entire two-dimensional area underneath the entire ROC curve (think integral calculus) from (0,0) to (1,1).

AUC provides an aggregate measure of performance across all possible classification thresholds. One way of interpreting AUC is as the probability that the model ranks a random positive example more highly than a random negative example. For example, given the following examples, which are arranged from left to right in ascending order of logistic regression predictions:

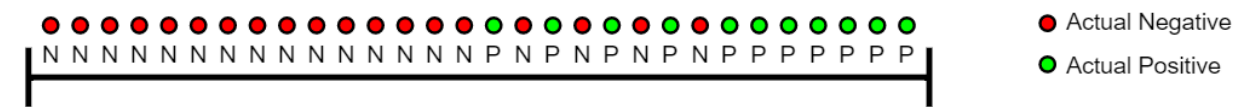


Figure 5.4: Predictions ranked in ascending order of logistic regression score.

5.5 Experiments

In the previous sections we discussed techniques that wouldn't be used at the same time because of the combination of them and more of changing hyperparameters in the experiments as we will discuss next.

5.5.1 Some explanations before start experiments

How the experiment will be evaluated

We will evaluate our experiments using validation accuracy matrices and by submitting prediction to “SIIM-ISIC Melanoma Classification” competition on Kaggle.

Ways of deciding the next step and what we will do

When we start building our model, we can either go about it randomly or intentionally. A random approach to model building is tempting to take. It requires less effort and still brings about results. A random approach would be when after each iteration of a model, you decide to try something new to make it better without much thought as to why. There is no systematic way to improve your model. It is simply this: an idea popped into your head that might improve the AUC-score so why not try it? It's easy and it works.

The intentional approach to building a model is using error analysis. Error analysis requires you to dig into the results of your model after each iteration. You look at the data and predictions on an observational level and form hypotheses as to why your model failed on certain predictions. Then you test your hypothesis by changing the model in a way that might fix that error and begin the next iteration. Each iteration of modeling becomes more time consuming with error analysis, but the final result is better and will likely arrive faster.

In error analysis the important factors are human level and train, dev and test accuracy. There is no human level in our problem so we will consider it to reach to 94.9% accuracy as the higher accuracy in the SIIM ISIC competition.

So, we choose to do error analysis after each experiment and depending on the results of error analysis we will try what to do next, but you may see that we sometimes do some experiments to run at the same time its changes was taken from error analysis trying to make the process faster.

Explanation about the distribution of train/div/test set

In all experiments you will see that train may be from a distribution that different from div and test set that they are always from the same distribution instead of any preprocessing like normalizing images this must be the same for the three sets.

5.5.2 Base Model

In this section, the model will be explained in detailed to avoid explaining the model add each experiment.

Let's describe the base model in the next steps:

- Import all libraries that will be used in the experiments from Tensorflow, Keras, SKLearn and Matplotlib, and import numpy, pandas and os.
- Use “marking.csv” from the dataset to create pandas' data frame that contains training images names and their labels.
- Train_test_split is used to split train data frame into train_df and validate_df to use validation data to evaluate the model at each step.
Use test_size equal 0.2 making the number of validation images equal 11445 that close to number of test set images and use a certain random state to see changes in experiments accuracy without outside effect.
- Print a figure using plot in matplotlib that explain the difference between the number of images in each of the two classes (benign, malignant).
- We use class weights to try to avoid the effect of this huge difference.
There are a lot of techniques to avoid this difference rather than class weights but using it backing two seeing a lot of previous work and experiments and the difficult to collect malignant images sometimes it takes years to collect the number of images that we want to balance our dataset.
- In this step we used the images data generator to collect data set images into a tensor data structure that is optimized to work with GPU in parallel by run batches of the data set at the same time, we assign image size to (224,224) and batch size to 32 batch.

In image data generator we rescale our images to make pixel values in range from zero to one and add augmentation to present it in the dataset at this step we use only zoom augmentation with ratio equal 0.3.

- In validation dataset we use image data generator with the same image size and batch size but only rescale the images without augmentation as it must be from the same distribution of the test set.
- The model explained in the model summary in figure 5.5 and the model compiled using binary_crossentropy loss function, Adam optimizer with LR equal 10^{-8} as small number in the first and as accuracy matrices we used Accuracy Under the Curve (AUC) matrix and False Negative matrix as the False classified malignant as benign as more hurtful for use.

| Layer (type) | Output Shape | Param # |
|--------------------------------|-----------------------|-----------|
| input_2 (InputLayer) | [(None, 224, 224, 3)] | 0 |
| conv2d_1 (Conv2D) | (None, 222, 222, 32) | 896 |
| conv2d_2 (Conv2D) | (None, 222, 222, 32) | 9248 |
| max_pooling2d (MaxPooling2D) | (None, 111, 111, 32) | 0 |
| conv2d_3 (Conv2D) | (None, 109, 109, 64) | 18496 |
| conv2d_4 (Conv2D) | (None, 109, 109, 64) | 36928 |
| max_pooling2d_1 (MaxPooling2D) | (None, 54, 54, 64) | 0 |
| conv2d_5 (Conv2D) | (None, 52, 52, 128) | 73856 |
| conv2d_6 (Conv2D) | (None, 52, 52, 128) | 147584 |
| max_pooling2d_2 (MaxPooling2D) | (None, 26, 26, 128) | 0 |
| conv2d_7 (Conv2D) | (None, 24, 24, 256) | 295168 |
| conv2d_8 (Conv2D) | (None, 24, 24, 256) | 590080 |
| max_pooling2d_3 (MaxPooling2D) | (None, 12, 12, 256) | 0 |
| flatten (Flatten) | (None, 36864) | 0 |
| dense (Dense) | (None, 4096) | 150999040 |
| dense_1 (Dense) | (None, 1) | 4097 |
| Total params: 152,175,393 | | |
| Trainable params: 152,175,393 | | |
| Non-trainable params: 0 | | |

Figure 5.5: base model summary

- Define callback to save model and model weights of the best validation accuracy after training the model.
- Fit the model with the train and validation datasets to start training the number of epochs depend on that the maximum number that can be assign must not exceeds the nine hours of training on GPU as constraints in Kaggle notebooks

the pass the number of training steps (`total_train//batch_size`), validation steps (`total_validate//batch_size`), callback variable and classes weights.

- Display loss, AUC accuracy, and False negative accuracy in the following figures.

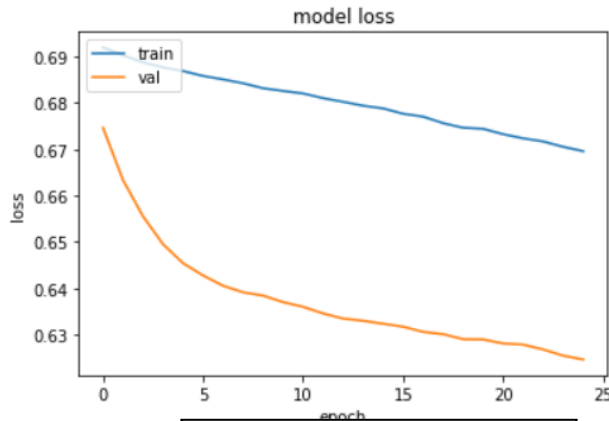


Figure 5.6: Model Loss

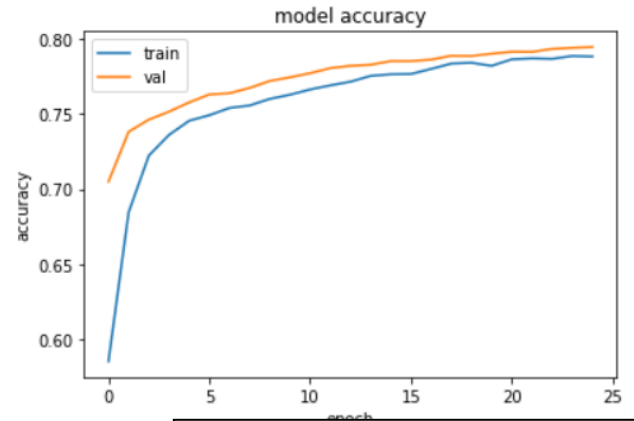


Figure 5.7: Model accuracy

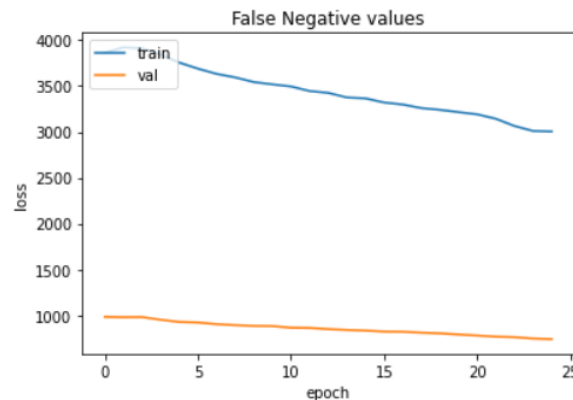


Figure 5.8: False Negative Values

- List images names in test image directory in `test_filenames` list and put this list in test data frame.
- For each image in `test_filenames` import the image using `Image.open` from `PIL` and resize it to `(224,224,3)`, convert it to numpy array and rescale it by dividing each pixel on `255.0` then using this array predict the probability of the image that identify it's benign or malignant append this prediction to "target" list.

- Add each prediction in the target list to test data frame in a column called target then save the data frame as a submission.csv to submit it into the competition to get test accuracy, Test accuracy of the base model is 0.7415 as private score and 0.7104 as public score.

This is the base model that used in the next experiments.

5.5.3 Experiment one

After analyzing the previous figure, we see that train is close to validation accuracy, so we face high bias low variance problem, depending on this we have to increase training one of the techniques to do so is reach to the optimal learning rate.

So, we try to reach to a learning rate close to the optimal by try to change learning rate.

First, we change learning rate to 10^{-7} , the result of this change is in the following figures:

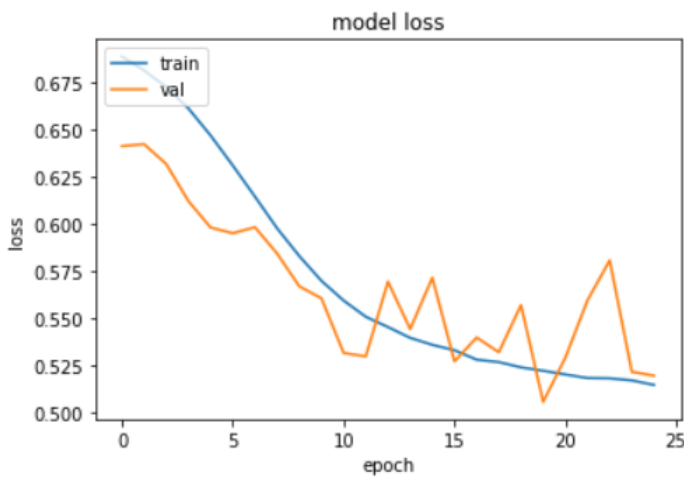


Figure 5.9: Model Loss

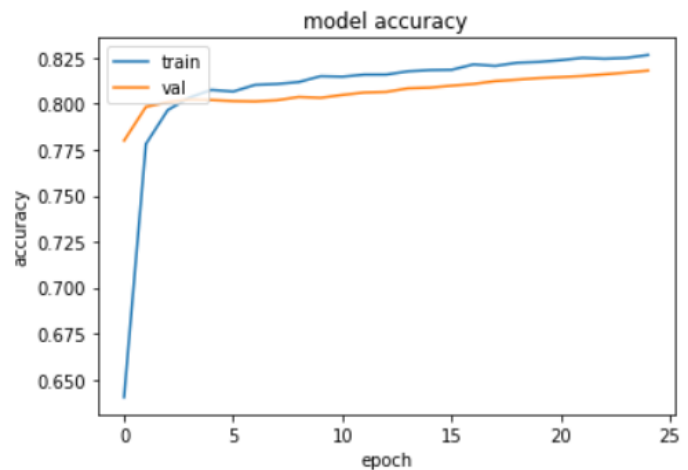


Figure 5.10: Model accuracy

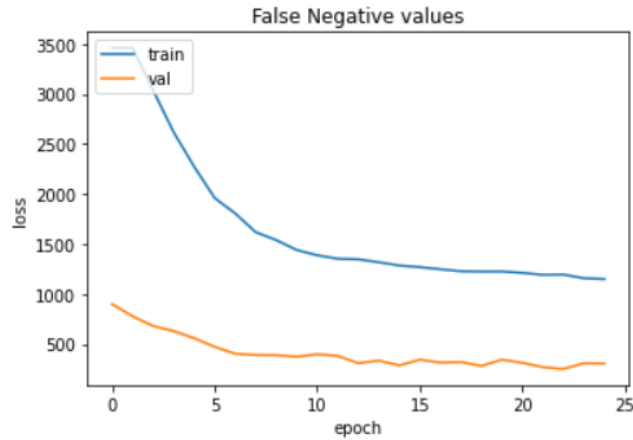


Figure 5.11: Model False negative values

After changing learning rate to 10^{-7} the accuracy increased by 3% in both train and validation, loss function decreases in training but become not stable in validation comparing to the previous model but it's a bigger decrease that what happen in the last change and false negative values decreased by 2000 values in the training and 500 in the validation.

From the previous the model has a better improve than the previous model but still in the same state of high bias low variance, so we will increase the learning rate again to 10^{-6} and the following figures present what happen.

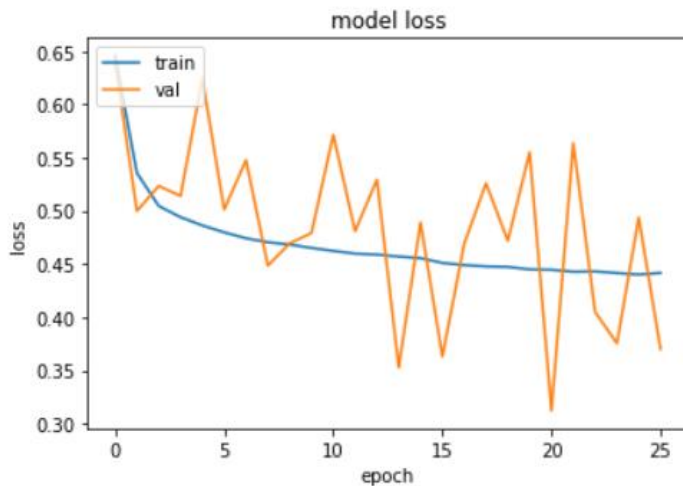


Figure 5.12: Model Loss

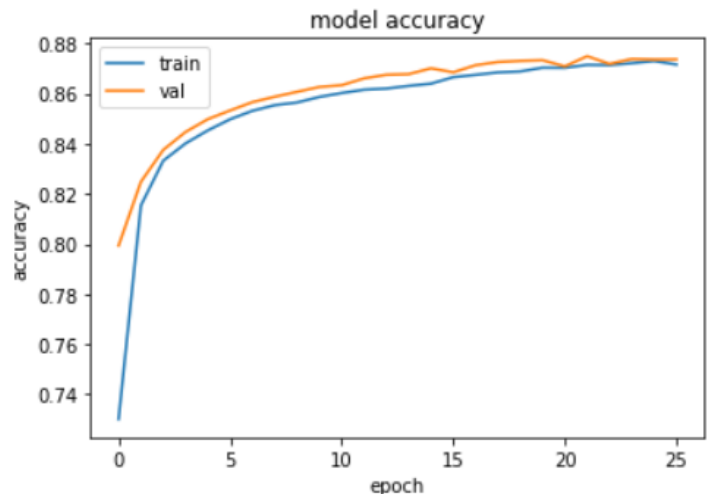


Figure 5.13: Model accuracy

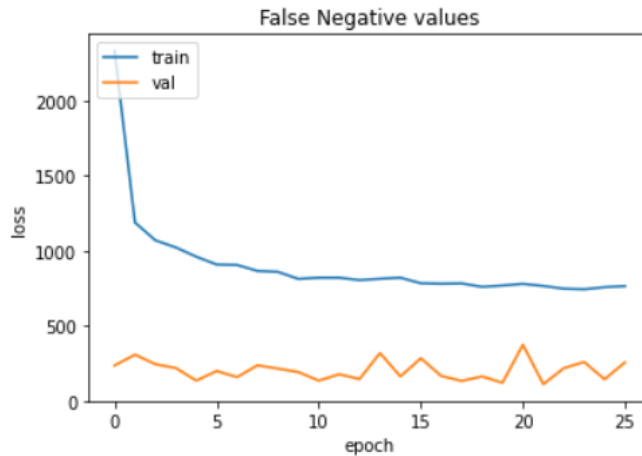


Figure 5.14: Model False negative

After changing learning rate to 10^{-6} the accuracy increased by 5% in both train and validation, loss function decreases in training and validation but become not stable in validation comparing to the previous model and false negative values decreased by a little in training values and almost decreased by 110 in the validation values.

From the previous the model has a better improve than the previous change but still in the same state of high bias low variance, so we will increase the learning rate again to 10^{-5} and the following figures present what happen.

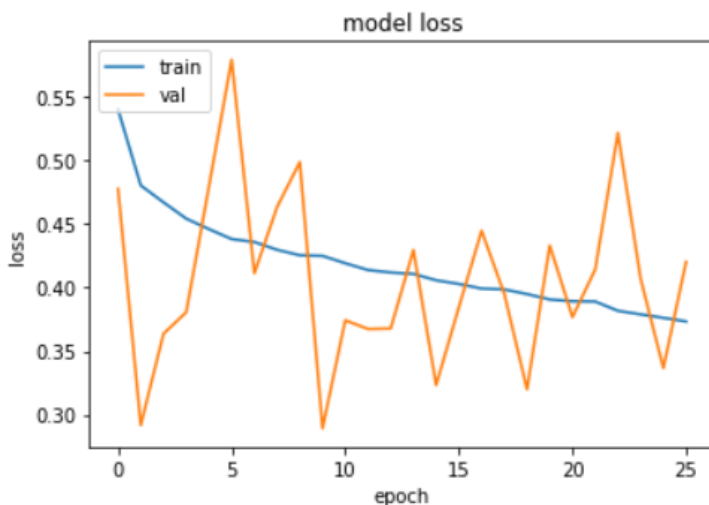


Figure 5.15: Model Loss

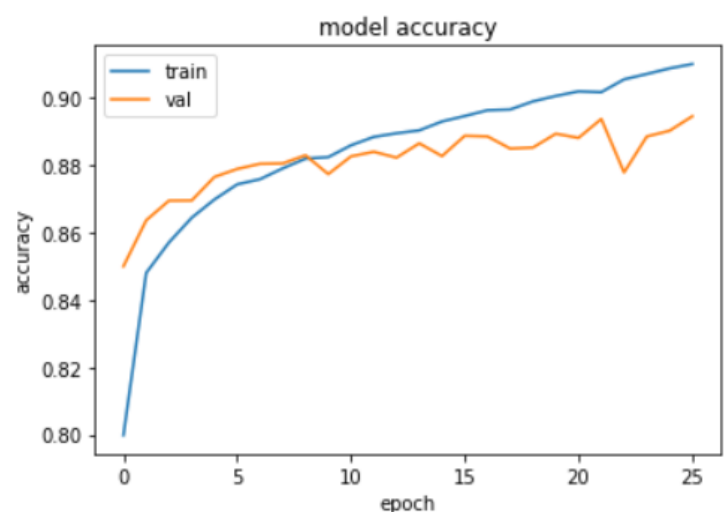


Figure 5.16: Model Accuracy

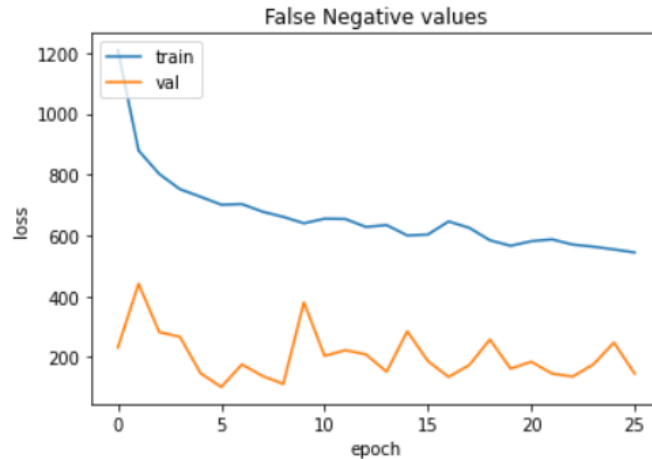


Figure 5.17: Model False negative values

After changing learning rate to 10^{-5} the accuracy increased by 4% in training and 2.5% in validation, loss function decreases in training and validation almost with the same difference as the previous change and false negative values almost doesn't change from the previous model.

From the previous the model has a better improve than the previous change but the difference between train and validation accuracy increases a little but still in the same state high bias and low variance, so we will increase the learning rate again to 10^{-4} and see whether the model will improve or not let's see in the next figures.

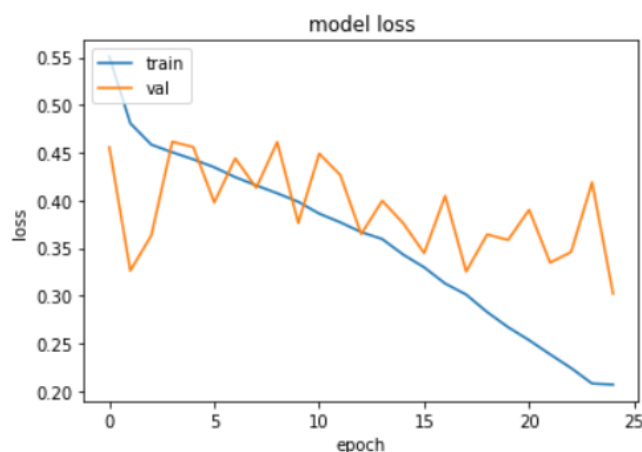


Figure 5.18: Model Loss

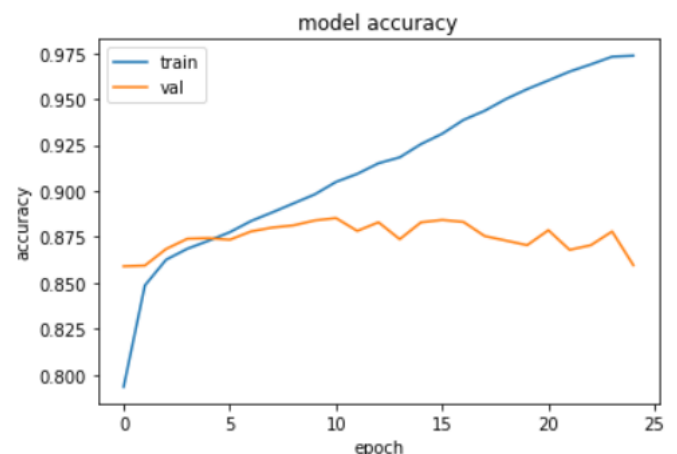


Figure 5.19: Model Accuracy

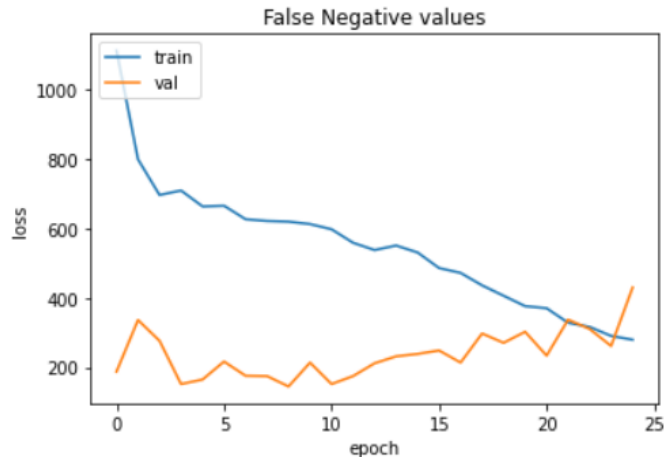


Figure 5.20: Model False negative values

After changing learning rate to 10^{-4} the accuracy increased by 6% in training and decreased by 1% in validation, loss function decreases a lot in training and almost is the same in validation and false negative decreases a lot in training and increases with small values in validation.

From the previous the model has a better improve in training but getting bad in validation, so we are facing overfitting at this step (low bias high variance).

The decrease of learning rate in this change makes the model overfit the training data so we will stop at this step and go to the next experiment.

After submitting the test prediction of each of the previous changes the test accuracy after changing learning rate to 10^{-7} is 0.7489 as private score and 0.7633 as public score, and with 10^{-6} learning rate is 0.8239 as private score and 0.8288 as public score, and with 10^{-5} learning rate is 0.8104 as private score and 0.8139 as public score, and with 10^{-4} learning rate is 0.7605 as private score and 0.7556 as public score.

5.5.4 Experiment two

In this experiment we will use additional augmentation trying to reduce overfitting that happen with 10^{-4} learning rate and this augmentation will increase the amount of data this increase in data may improve the accuracy of the model.

So, next we will discuss what happen when adding horizontal and vertical flip to the augmentation with learning rates 10^{-6} , 10^{-5} , 10^{-4} .

We will start with learning rate 10^{-6} , the following figures present what happen.

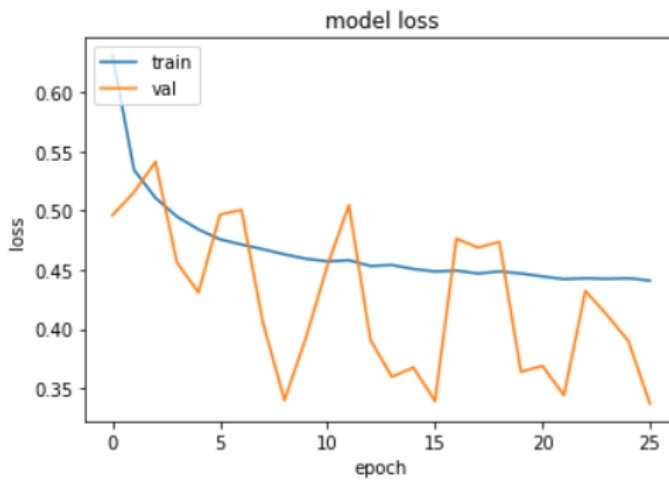


Figure 5.21: Model Loss

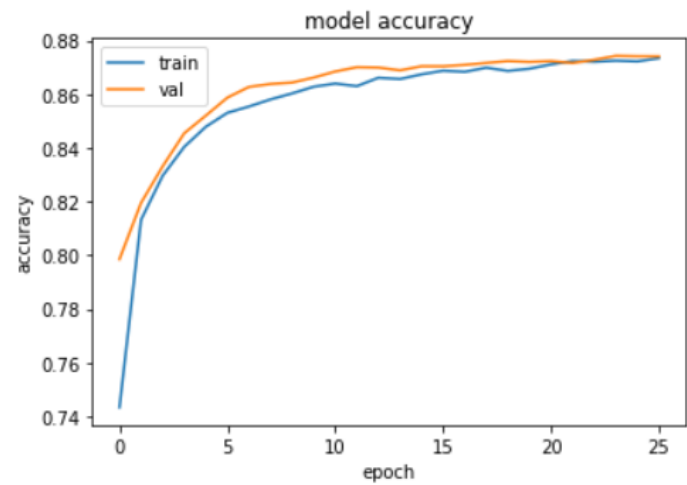


Figure 5.22: Model Accuracy

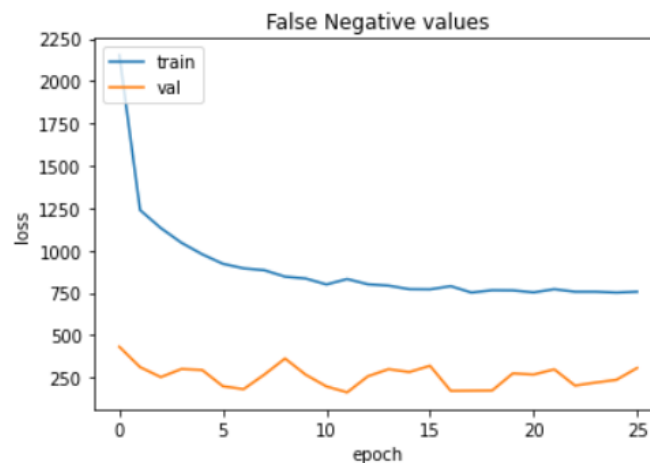


Figure 5.23: Model False negative values

After adding augmentation with learning rate 10^{-6} the training and validation accuracy almost the same without augmentation and also the loss function but the False Negative values increased in validation and almost the same in training.

From the previous after adding augmentation with learning rate 10^{-6} there is no improve in the model and it still in the same state high bias low variance.

Next, we will see what happen when adding augmentation with learning rate 10^{-5}

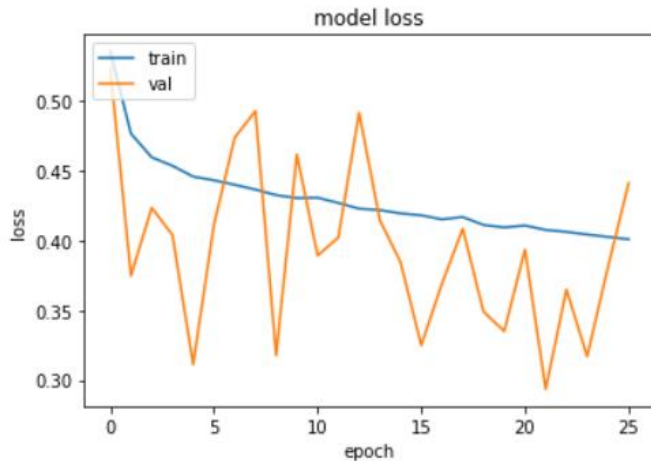


Figure 5.24: Model Loss

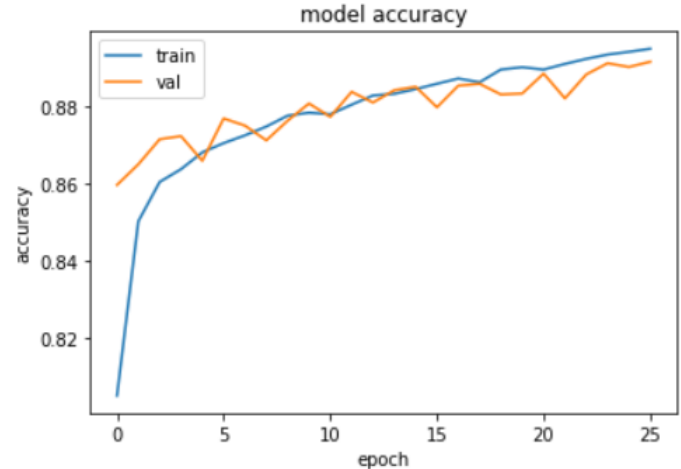


Figure 5.25: Model Accuracy

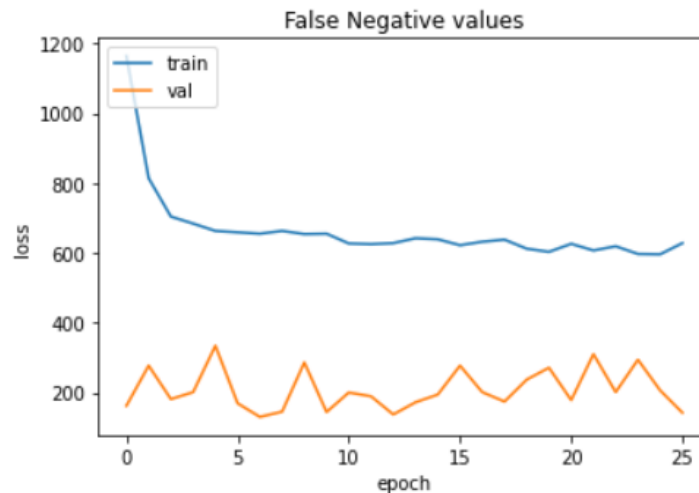


Figure 5.26: Model False negative values

After adding augmentation with learning rate 10^{-5} the training accuracy decreases by 1% and validation accuracy almost the same as without augmentation, the loss function increased a little bit we can say that it is not changed, and the False Negative values increased in training and validation.

From the previous after adding augmentation with learning rate 10^{-5} there is improve in the model by decreasing overfitting and it is in the state high bias low variance.

Next, we will see what happen when adding augmentation with learning rate 10^{-4}

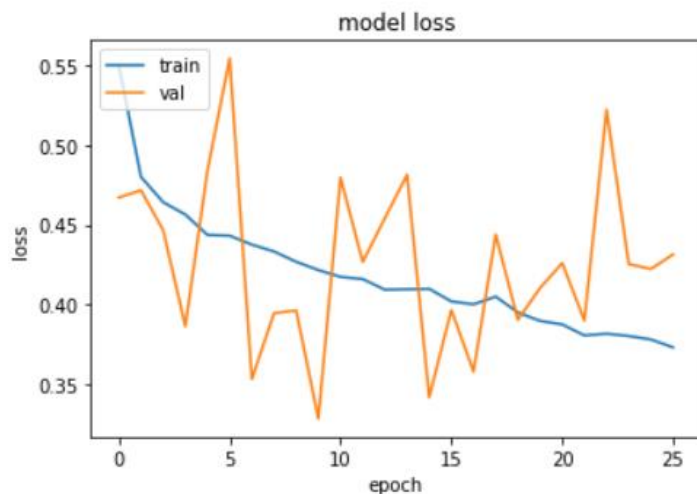


Figure 5.27: Model Loss

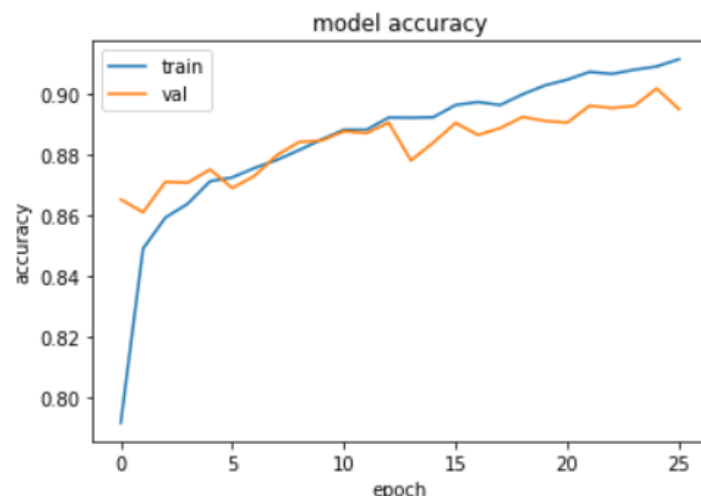


Figure 5.26: Model Accuracy

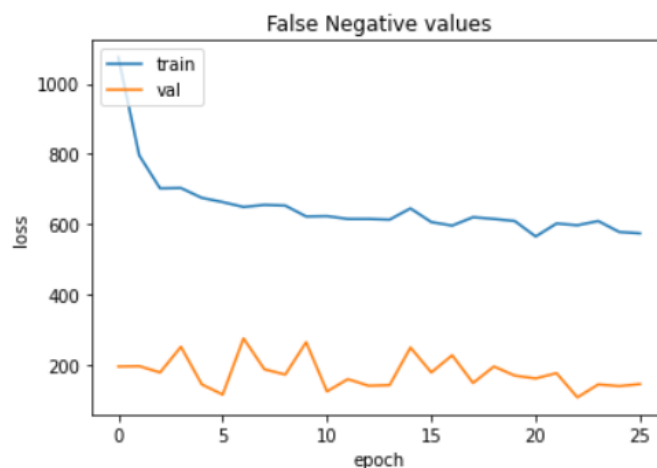


Figure 5.29: Model False negative values

After adding augmentation with learning rate 10^{-4} the training accuracy decreases by 6% and validation accuracy increases by 4%, the loss function increased in both training and validation, and the False Negative values increased in training and decreased in validation.

From the previous after adding augmentation with learning rate 10^{-4} there is improve in the model by decreasing overfitting and it is become in the state high bias low variance.

Adding augmentation improves models by decreasing overfitting in the next experiment we will try to increase the mode accuracy.

After submitting the test prediction of each of the previous changes the test accuracy after adding augmentation with learning rate 10^{-6} 0.8191 as private score and 0.8325 as public score, and with 10^{-5} learning rate is 0.8287 as private score and 0.8507 as public score, and with 10^{-4} learning rate is 0.8086 as private score and 0.8356 as public score.

5.5.5 Experiment three

In this experiment we will try two techniques to increase accuracy, we will make the neural network bigger and use the weights of the model after each change as initialization to the next to train the model more hours to see the effect of this.

Explanation: we save the model with the higher validation accuracy to get the model weights in a point with least overfitting and run the next model using these weights as initialization with some changes if the model starts to overfit the training data.

We start with rescaling, no augmentation and use a bigger model, the model summary presented in figure 5.30.

| Layer (type) | Output Shape | Param # | | | |
|---------------------------------|-----------------------|---------|---------------------------------|---------------------|-----------|
| ===== | | | | | |
| input_2 (InputLayer) | [(None, 224, 224, 3)] | 0 | conv2d_20 (Conv2D) | (None, 28, 28, 512) | 1180160 |
| conv2d_13 (Conv2D) | (None, 224, 224, 64) | 1792 | conv2d_21 (Conv2D) | (None, 28, 28, 512) | 2359808 |
| conv2d_14 (Conv2D) | (None, 224, 224, 64) | 36928 | conv2d_22 (Conv2D) | (None, 28, 28, 512) | 2359808 |
| max_pooling2d_5 (MaxPooling 2D) | (None, 112, 112, 64) | 0 | max_pooling2d_8 (MaxPooling 2D) | (None, 14, 14, 512) | 0 |
| conv2d_15 (Conv2D) | (None, 112, 112, 128) | 73856 | conv2d_23 (Conv2D) | (None, 14, 14, 512) | 2359808 |
| conv2d_16 (Conv2D) | (None, 112, 112, 128) | 147584 | conv2d_24 (Conv2D) | (None, 14, 14, 512) | 2359808 |
| max_pooling2d_6 (MaxPooling 2D) | (None, 56, 56, 128) | 0 | conv2d_25 (Conv2D) | (None, 14, 14, 512) | 2359808 |
| conv2d_17 (Conv2D) | (None, 56, 56, 256) | 295168 | max_pooling2d_9 (MaxPooling 2D) | (None, 7, 7, 512) | 0 |
| conv2d_18 (Conv2D) | (None, 56, 56, 256) | 590080 | flatten (Flatten) | (None, 25088) | 0 |
| conv2d_19 (Conv2D) | (None, 56, 56, 256) | 590080 | dense (Dense) | (None, 4096) | 102764544 |
| max_pooling2d_7 (MaxPooling 2D) | (None, 28, 28, 256) | 0 | dense_1 (Dense) | (None, 4096) | 16781312 |
| | | | dense_2 (Dense) | (None, 1) | 4097 |
| | | | ===== | | |
| | | | Total params: | 134,264,641 | |
| | | | Trainable params: | 134,264,641 | |
| | | | Non-trainable params: | 0 | |

Figure 5.30: Model summary

The effect of the above changes after 55 epochs, shown in the next figures.

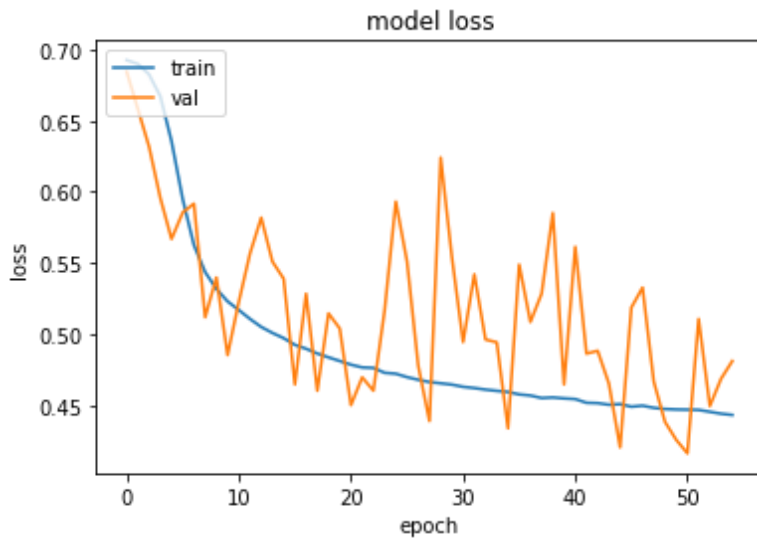


Figure 5.31: Model Loss

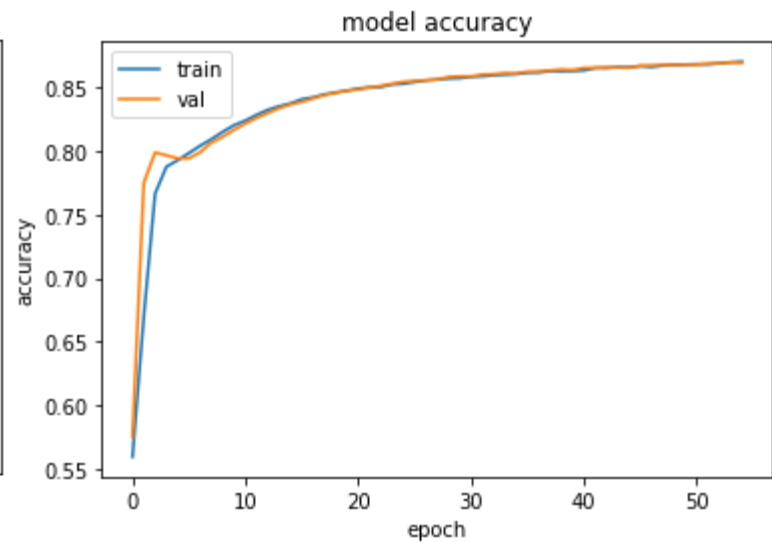


Figure 5.32: Model Accuracy

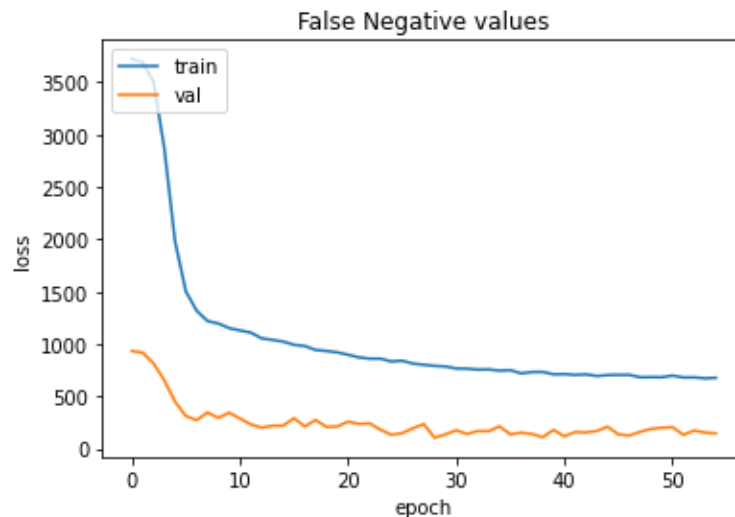


Figure 5.33: Model False negative values

From the figures the model performance very good as start because there is no overfitting it do very well in both train and validation, the model loss slightly bigger and the False Negative values is almost in the same range.

Every two changes of the next changes trained in parallel to see the effect of them on the accuracy matrices and the loss.

- Firstly, use the previous model weights and run additional 42 epoch, the effect shown in the next figures.

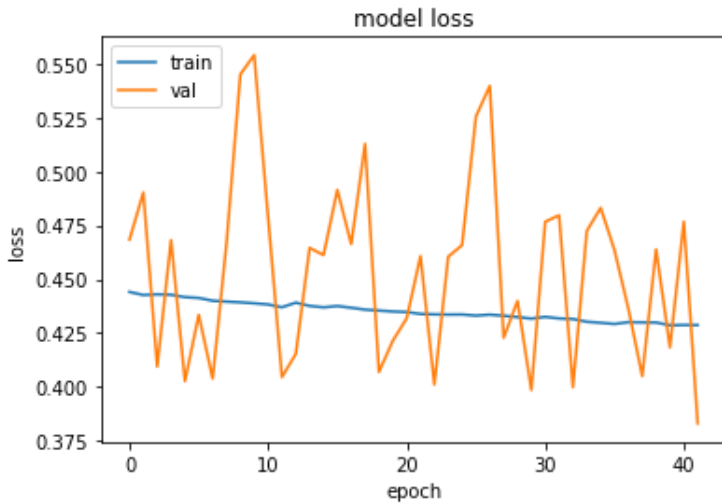


Figure 5.34: Model Loss

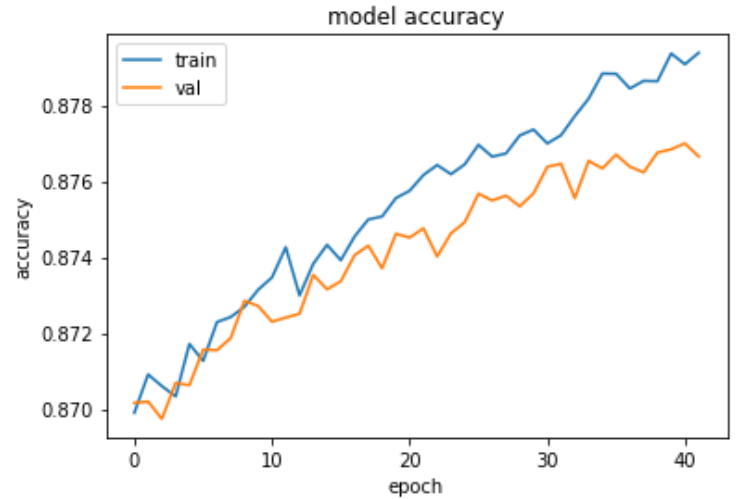


Figure 5.35: Model Accuracy

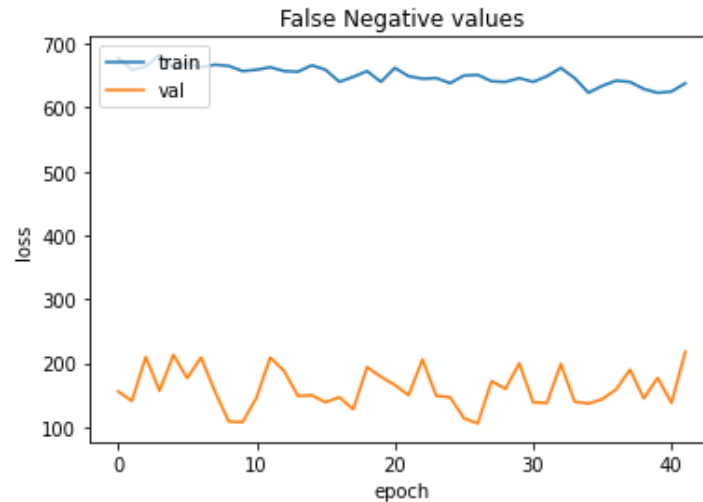


Figure 5.36: Model False negative values

The performance of the model slightly improves, the validation loss not stable but finally it is decreased, and False Negative values decreases with good difference comparing to the previous.

From the above the model doesn't improve comparing to the number of epochs and GPU hours.

- Secondly, add three dense layers with dropout layers to the model and run with its weights as initialization.

The effect of this change shown in the next figures.

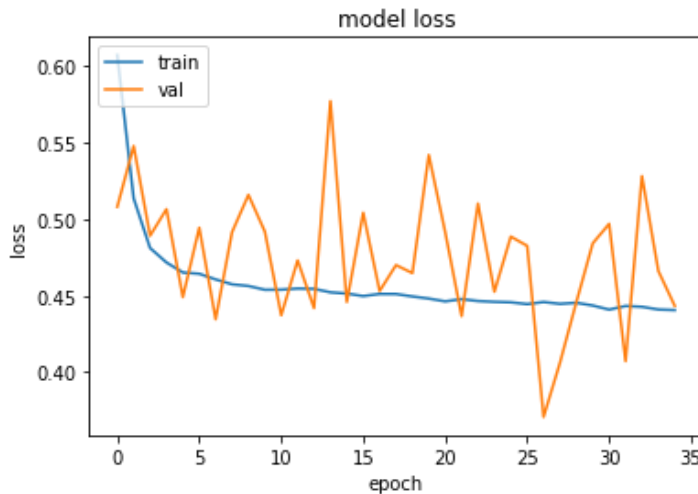


Figure 5.36: Model Loss

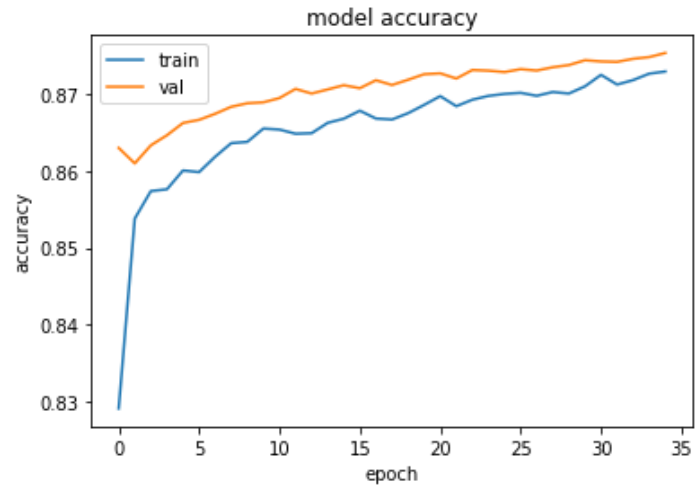


Figure 5.38: Model Accuracy

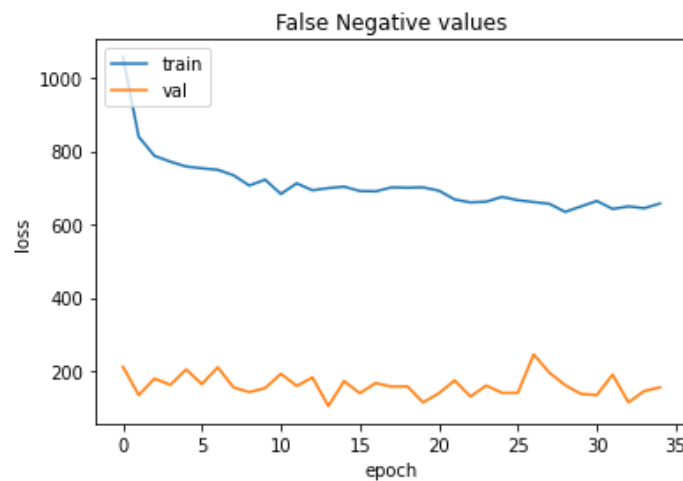


Figure 5.39: Model False negative values

The model does well in the validation accuracy than the training accuracy but comparing to its parallel change the validation accuracy almost the same, The Loss and False Negative values decreased.

From the above the model improves in overfitting but totally it doesn't improve comparing to the number of epochs and GPU hours.

- Firstly, add more augmentation to the first model of the previous two changes and use 10^{-5} learning rate and like all in this experiment weights of the previous model will be used as initialization to the model.

The next figures present what happen after train 37 epoch with this model.

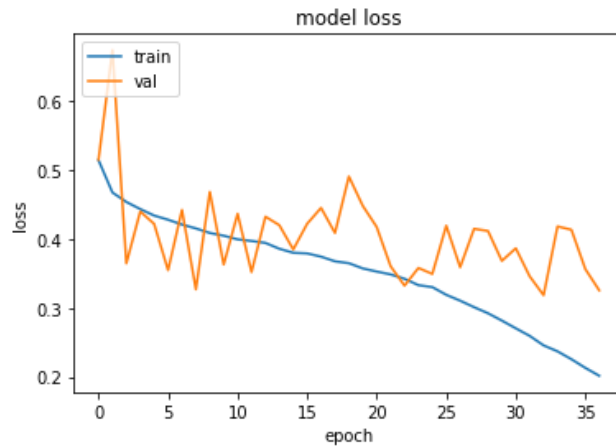


Figure 5.40: Model Loss

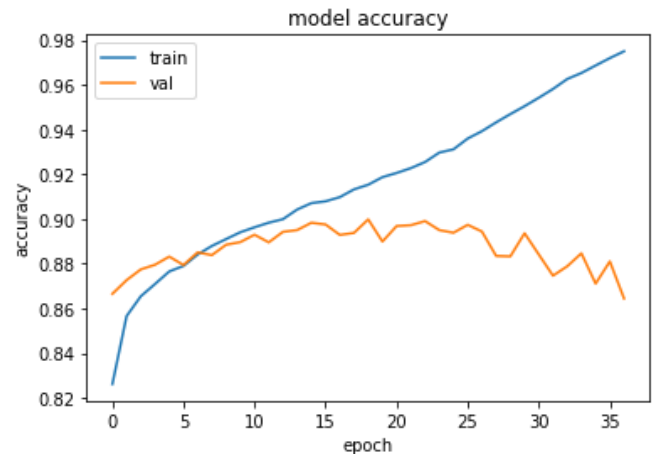


Figure 5.41: Model Accuracy

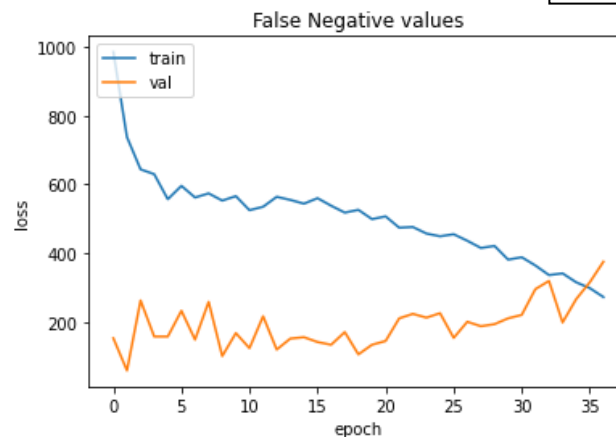


Figure 5.42: Model False negative values

From the figures model Accuracy start increasing in both train and validation but after ten epochs increases in training but still in the same range in validation, loss firstly decreases in both but after some epochs training accuracy decreases and validation still in the same values, and the False Negative values start with decreasing in validation but after almost 20 epochs start increasing but training continued decreasing.

From the above overfit the data, but since we save the model with the best validation accuracy, it will not affect the next model much.

- Secondly, use the first change but with 10^{-4} learning rate.

The next figures present what happen after train 47 epoch with this model.

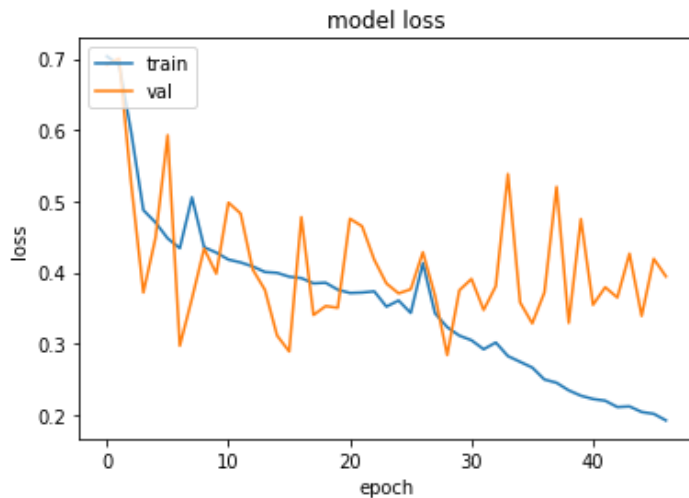


Figure 5.43: Model Loss

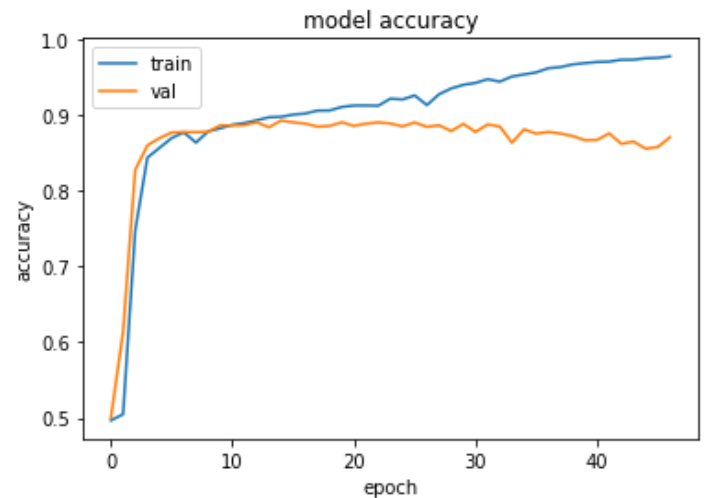


Figure 5.44: Model Accuracy

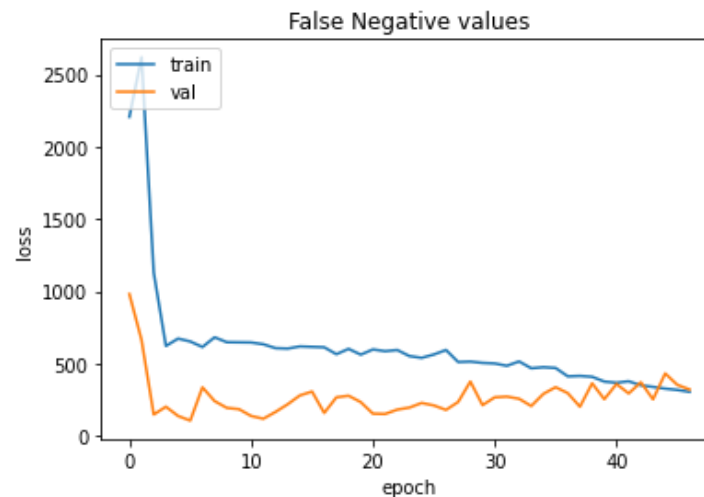


Figure 5.45: Model False negative values

From the figures model Accuracy start increasing in both train and validation but after almost fifteen epoch increases in training but slightly decreased in validation, loss firstly decreases in both but after some epochs training accuracy decreases and validation still in the same values but the loss here is better than the previous change, and the False Negative values start with decreasing in validation but after almost 20 epochs start slightly increasing but training continued decreasing.

From the above overfit the data, comparing to the previous model this is better in loss but the previous has better accuracy so we will continue with both models and see.

In the next two changes add two dropout layers in the last two blocks in the model and add zoom augmentation with range [0.5,1.0].

- Firstly, with learning rate 10^{-5} .

The next figures present what happen after train 30 epoch with this model.

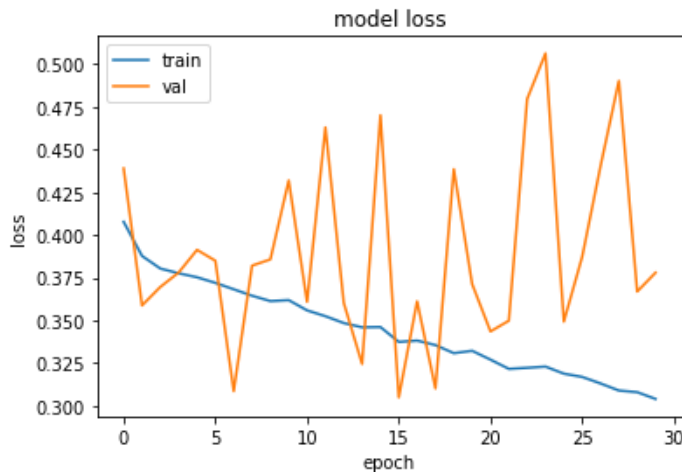


Figure 5.46: Loss

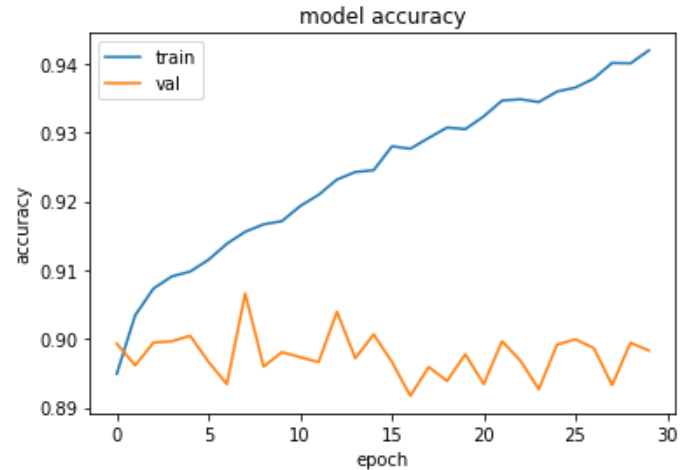


Figure 5.47: Model Accuracy

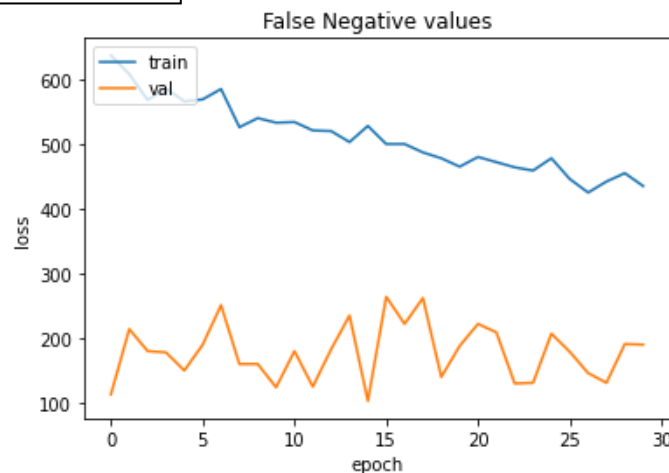


Figure 5.48: Model False negative values

From the figures training accuracy increased in almost straight line but validation in most time it decreased, it increased slightly in some experiments, the Loss function in training decreases in almost a straight line but in validation it frequently changes sometimes takes small values and sometime a very high values, and False Negative values it decreases in both but in validation it isn't stable.

- Secondly, with 10^{-4} learning rate

The next figures present what happen after train 32 epoch with this model.

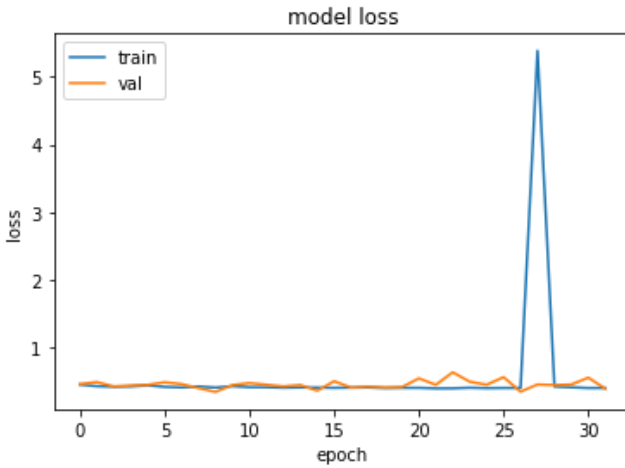


Figure 5.49: Model Loss

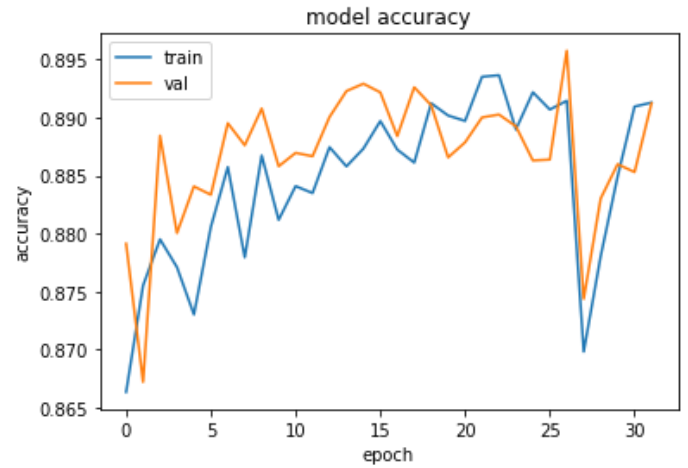


Figure 5.50: Model Accuracy

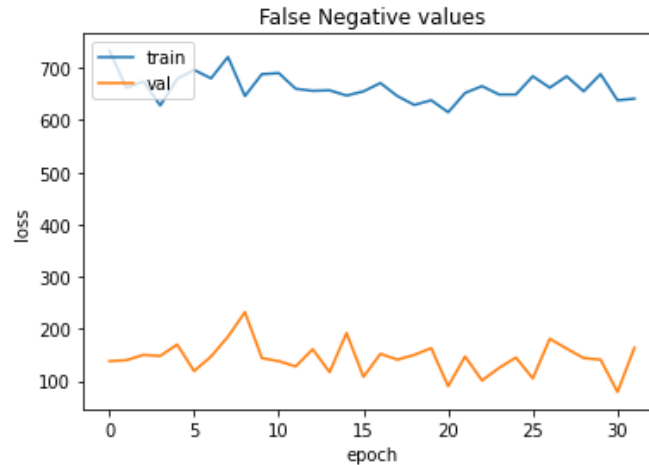


Figure 5.51: Model False negative values

From the figures train and validation accuracy increases in the same range but both are not stable but at most of the epochs the validation accuracy is higher, the loss is very good in both, training increases suddenly but quickly back to the small values, and False Negative values in validation is very good but in training if compared to the last model, it's not very good.

Next, we used the weights of the first change (the highest validation accuracy) and train it after standardizing the dataset images using training images STD and mean.

The next figures present what happen after train 28 epoch with this model.

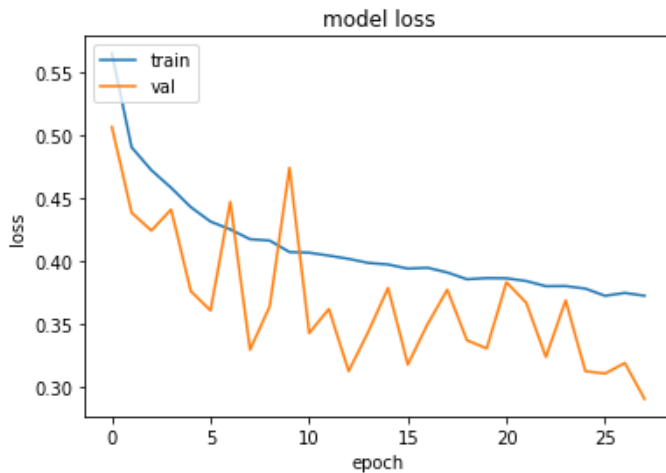


Figure 5.52: Model Loss

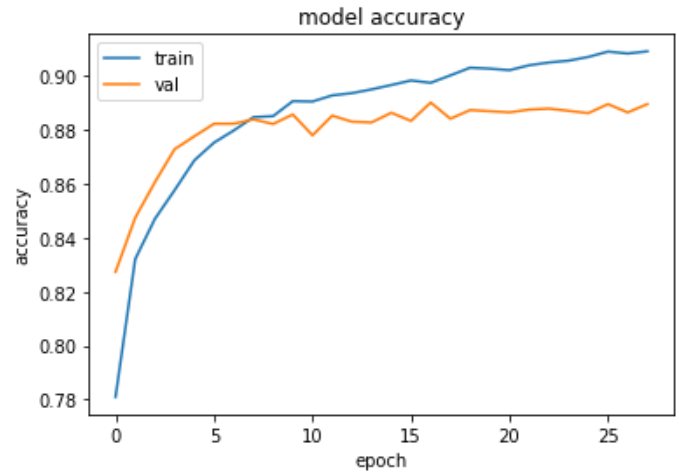


Figure 5.53: Model Accuracy

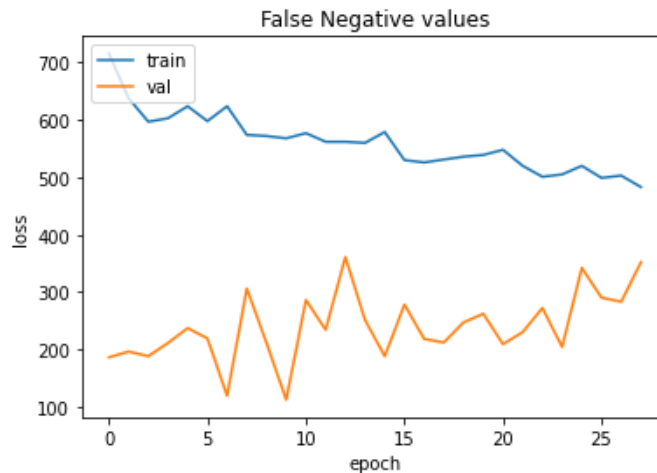


Figure 5.54: Model False negative values

From the figures the model accuracy increases but after epoch nine the validation almost fixed and training continues increasing, Loss is not very well in both train and validation, but it isn't bad, and False Negative values did like Accuracy.

After this change we did a lot of experiments but the performance in them is like what shown in the next figures.

As shown in the figures training and validation are very different.

So, to solve this problem we must connect the output of network layers not to only the next layer but to of the layers in the same block of layers.

There are some networks that have this advantage like ResNet and DenseNet.

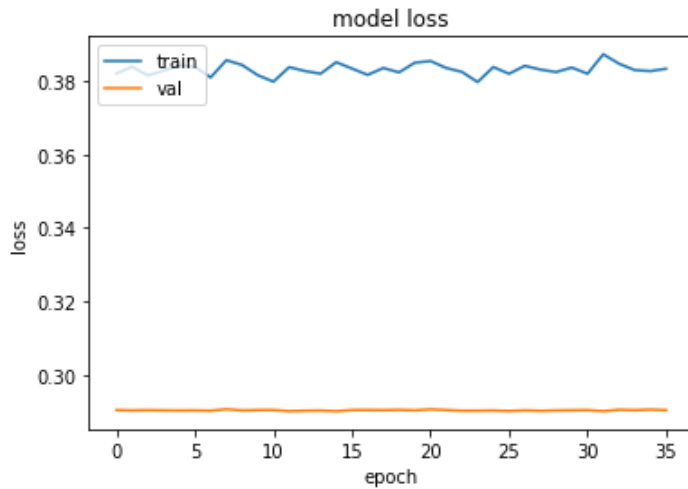


Figure 5.55: Model Loss

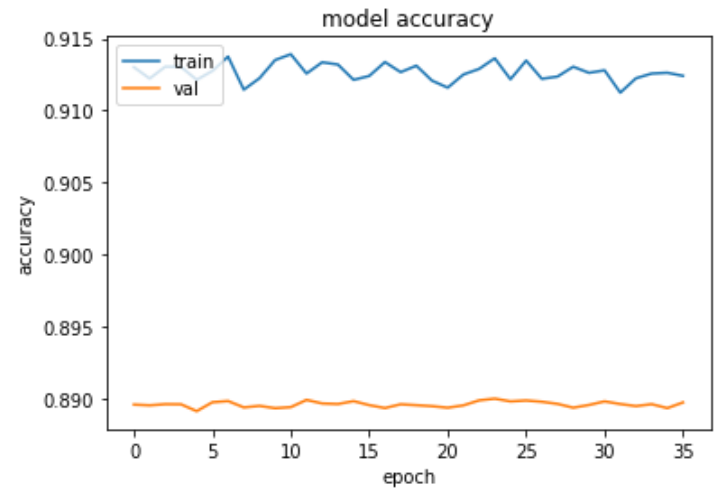


Figure 5.56: Model Accuracy

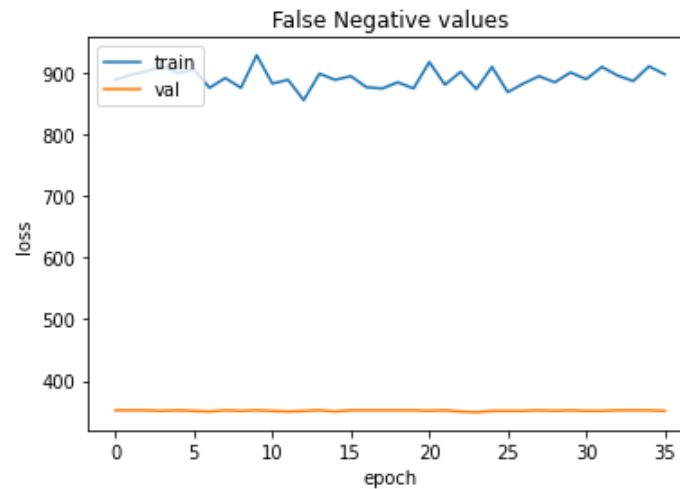


Figure 5.57: Model False negative values

5.5.6 Experiment Four

In this experiment we will use ResNet50 architecture with constant and cycling learning rate and other hyperparameter tuning.

1. Use the model with augmentation used in the last experiment, 10^{-5} constant learning rate, and standardization.

The next figures present what happen after train 24 epoch with this model.

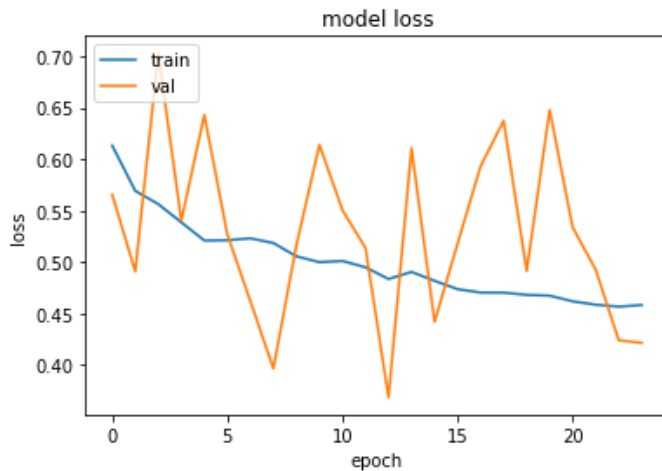


Figure 5.58: Model Loss

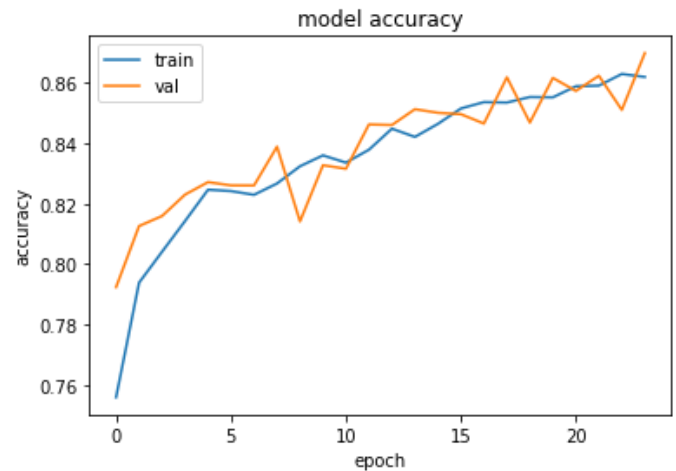


Figure 5.59: Model Accuracy

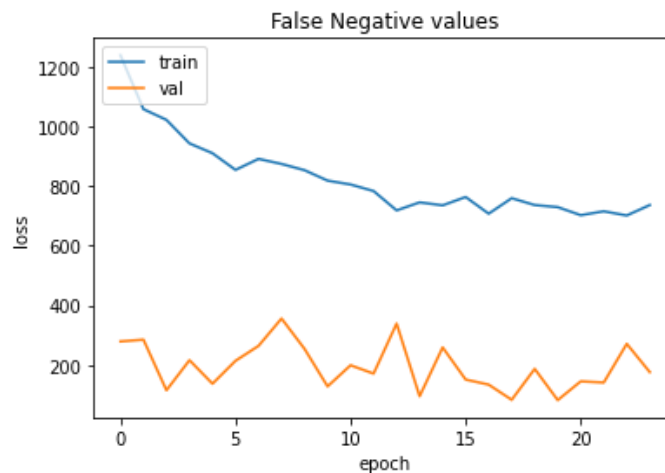


Figure 5.60: Model False negative values

The model performance is good and not overfit training data so we will only change zoom augmentation value to range and run again.

As shown in the figures (61,62,63) after 35 epoch the model did very well in Accuracy until epoch 15 and then it overfit the training data, in loss the results as high in both training and validation, and in False Negative values it did well in validation but not very good in training, but it is acceptable.

But until now it is not satisfied, we change a lot of hyper parameter and run the model many times on high number of GPU hours and still in the same range of validation accuracy next we will try using cycling learning rate and see the difference.

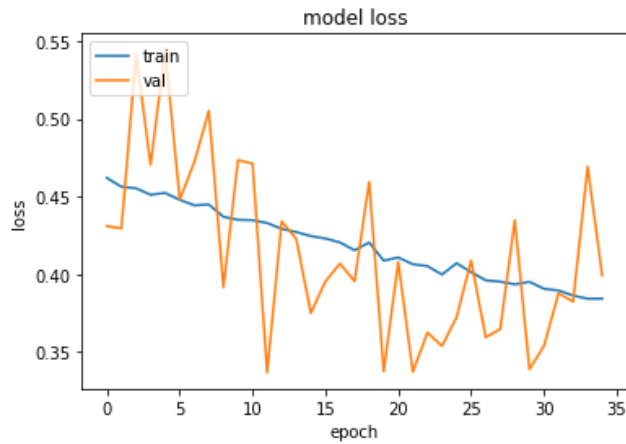


Figure 5.60: Model Loss

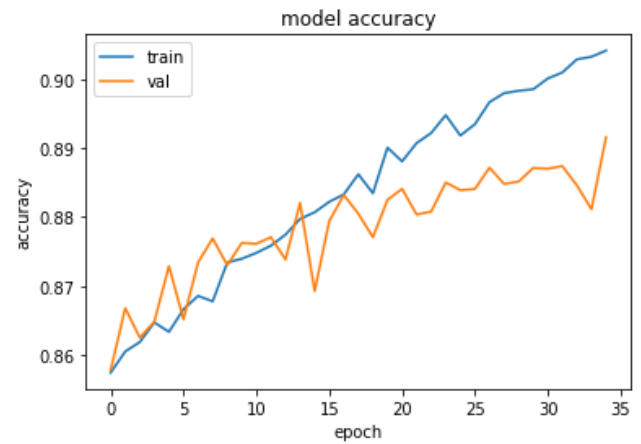


Figure 5.62: Model Accuracy

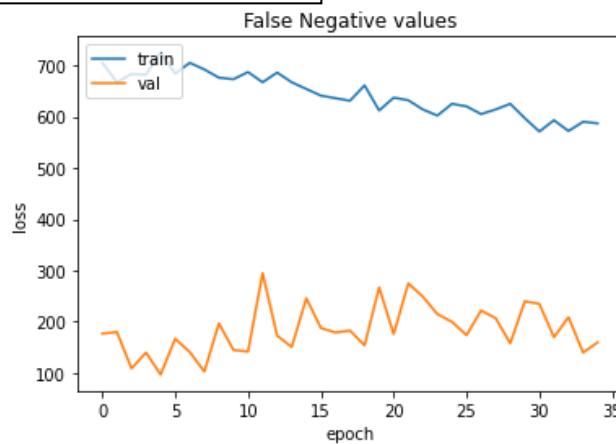


Figure 5.63: Model False negative values

2. Using cycling learning rate with 10^{-2} as max and 10^{-5} as min this numbers defined after a lot of experiments to define the best values.

The next figures present what happen after train 28 epoch with this model.

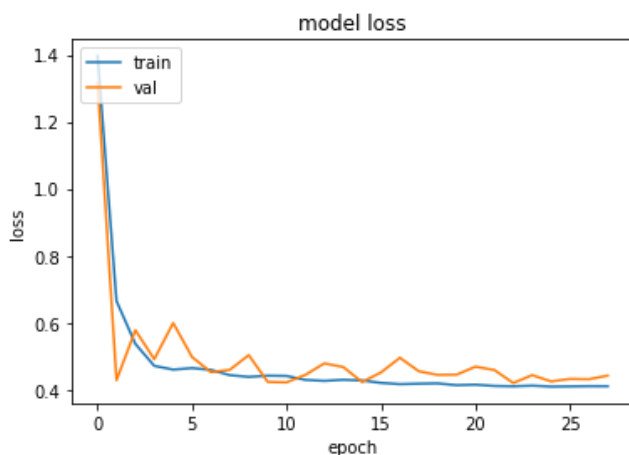


Figure 5.64: Model Loss

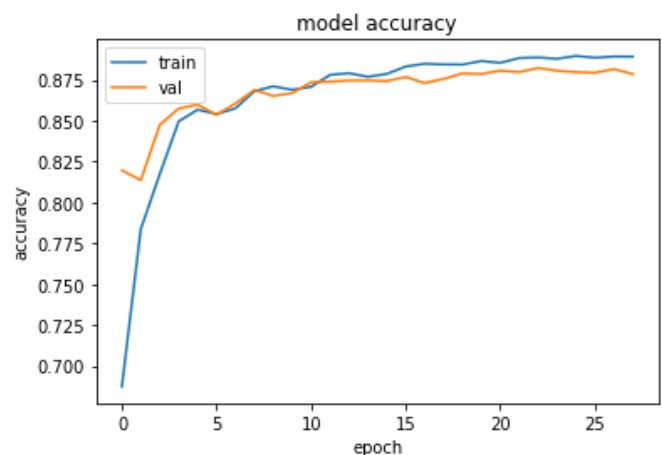


Figure 5.65: Model Accuracy

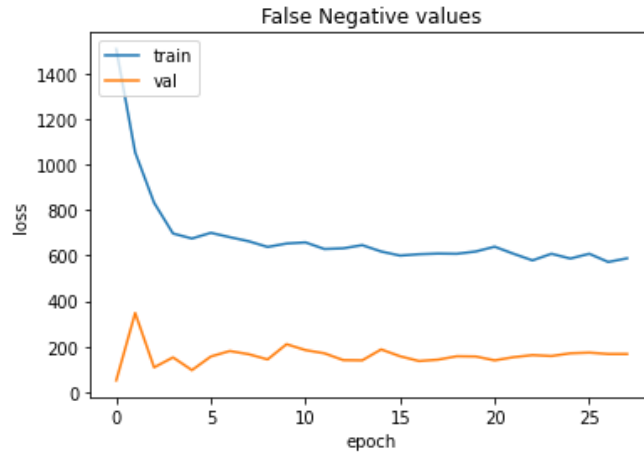


Figure 5.66: Model False negative values

The model performance is very good in all accuracy, loss and false negative values comparing to the model with constant learning rate.

Next, we will add some linear layers with dropout to run the model a lot 71 epoch on three times to see the effect of running the model on a lot of GPU hours.

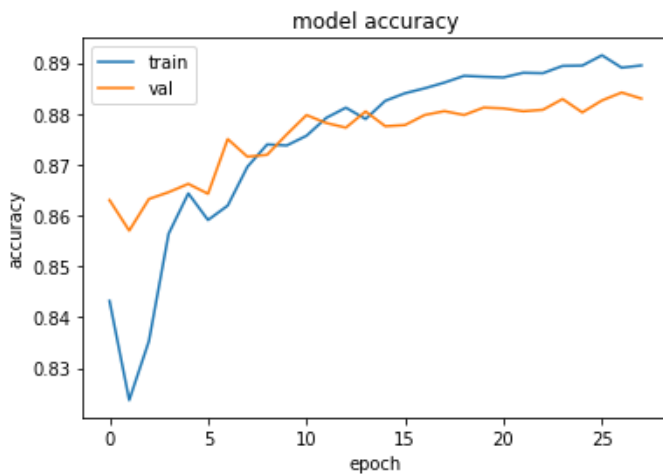


Figure 5.67: Model Accuracy of first 28

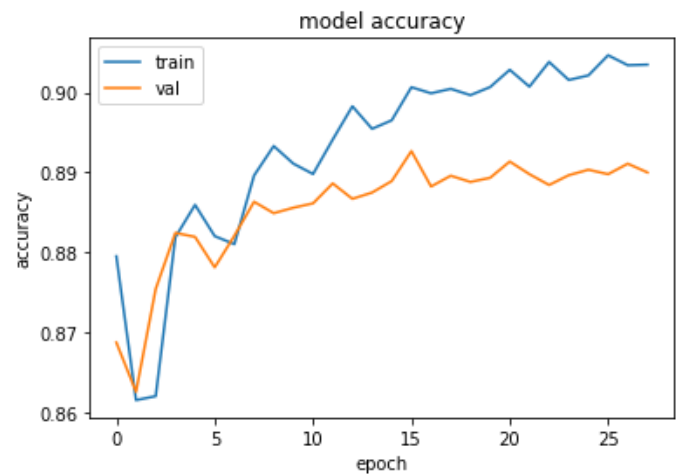


Figure 5.68: Model Accuracy of second 28

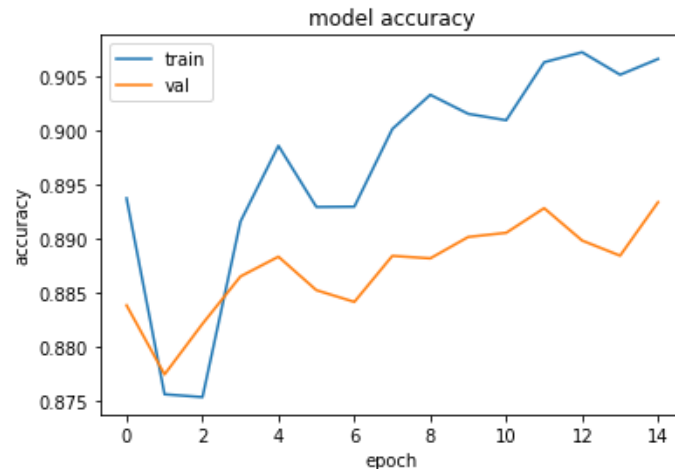


Figure 5.68: Model Accuracy of the last 15 epoch

The model improves slightly so we will try to go with transfer learning and see the effect.

5.5.7 Experiment five

A lot of tuning happened but the model still needs to be improved so next we will use weights of “imagenet” dataset as initialization to the model (VGG16), we back to VGG16 to see first what will be the improve with this “simple” model.

The changes presented next:

- First the images standardized as in the previous experiment, use zoom, vertical flip and horizontal flip as augmentation, back to constant learning rate instead of cycling learning rate and VGG16 model with three dense layers with dropout.

The next figures present what happen after train 26 epoch with this model.

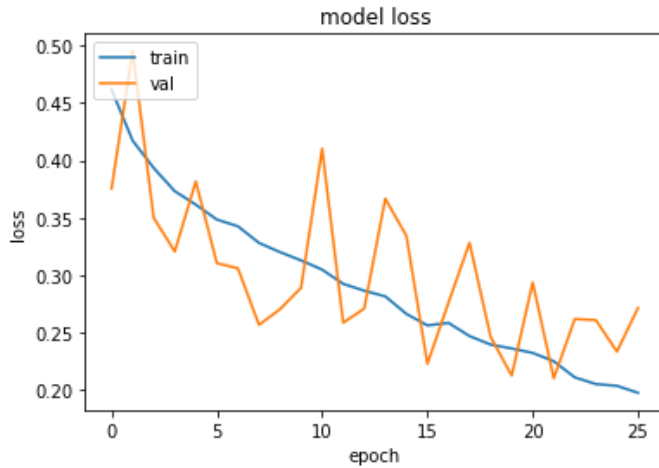


Figure 5.69: Model Loss

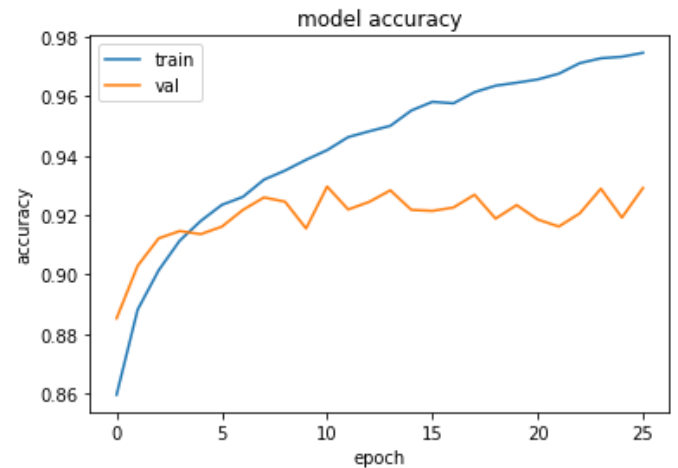


Figure 5.70: Model Accuracy

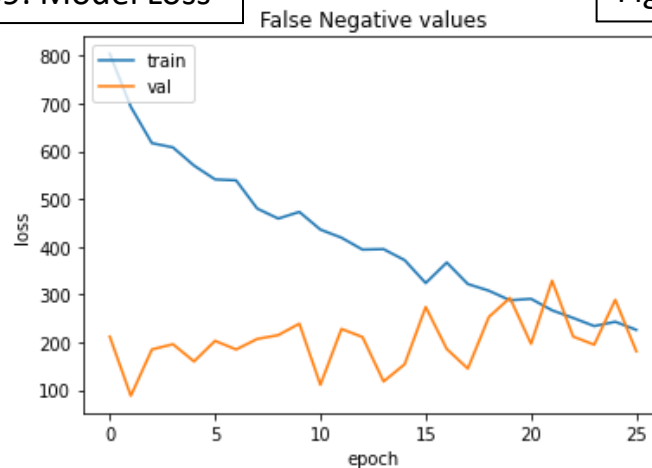


Figure 5.71: Model False negative values

The training and validation accuracy reach to 97.5% and 92.4% respectively, the Loss in training is so good but not in validation.

The test accuracy in the competition is 0.8837 as private score and 0.9149 as public score.

So, the model reaches to a good point, but we have to treat this overfitting.

- Add rotation, shear range, width shift range and height shift range augmentation.

The next figures present what happen after train 26 epoch with this model.

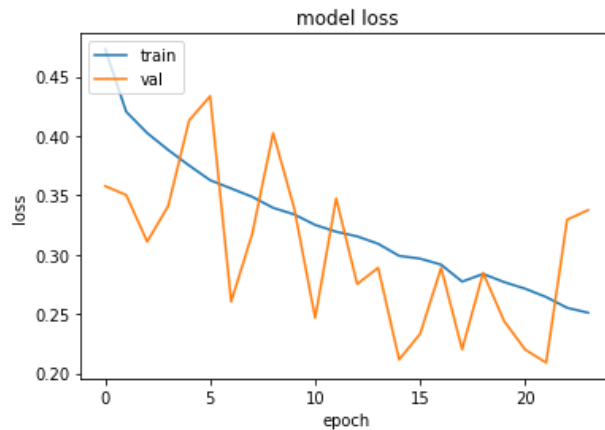


Figure 5.72: Model Loss

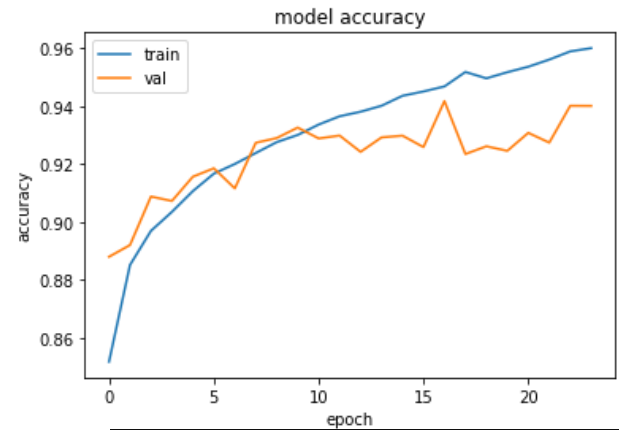


Figure 5.73: Model Accuracy

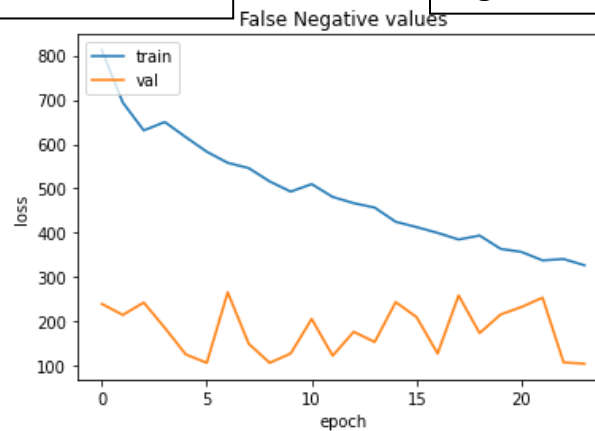


Figure 5.74: Model False negative values

The model improved as shown in the figures the training accuracy decreased and validation accuracy increased, the Loss of validation is better than the previous and in False negative values its so good comparing to the previous.

The test accuracy in the competition is 0.8813 as private score and 0.8890 as public score.

The test accuracy of the previous model is slightly better than this model so next we will try to merge the two models.

- Take the weights of the first model and use it as initialization with the second model augmentation.

The next figures present what happen after train 27 epoch with this model.

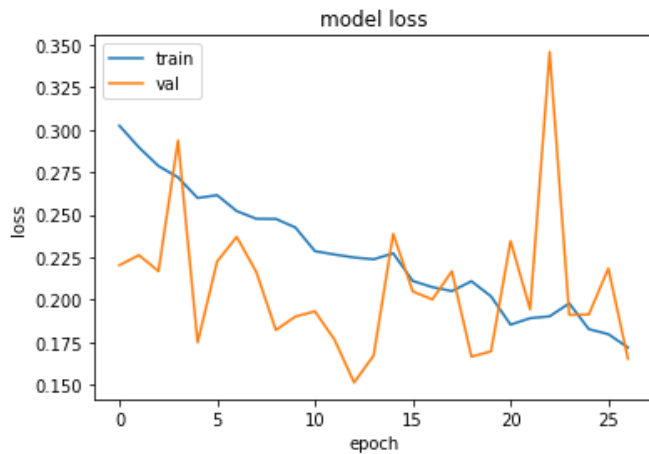


Figure 5.75: Model Loss

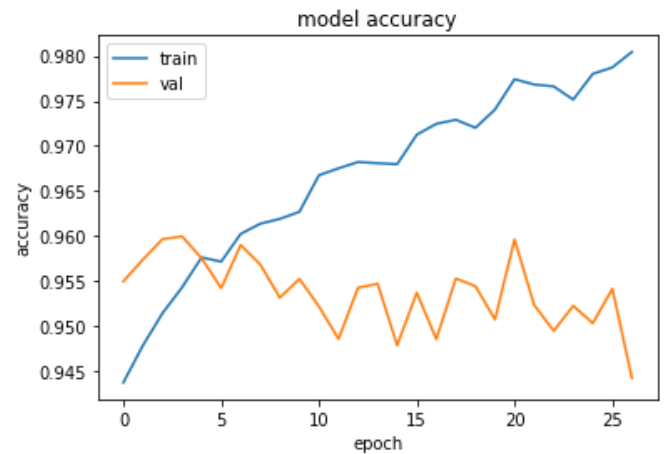


Figure 5.76: Model Accuracy

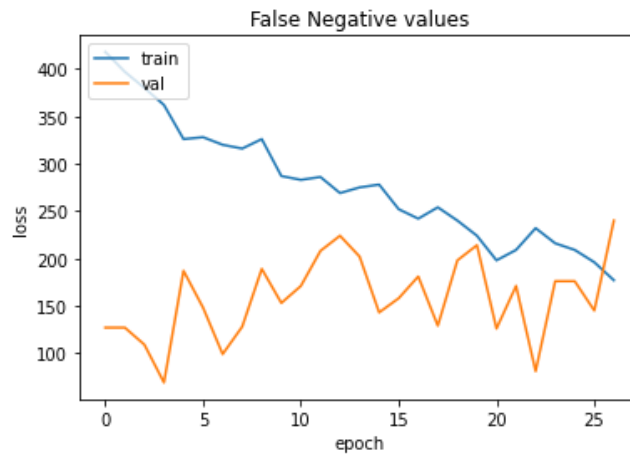


Figure 5.77: Model False negative values

From these figures and the previous experiments when retrain the model with another model weights it always overfit the data may this because the weights still not as good to use as initialization.

5.5.8 Experiments Six

The main idea in this experiment is using the “imagenet” dataset weights as initialization to the model with a lot of tuning to reach to the best model.

Firstly, we will use VGG16 model and do a lot of tuning then go with a more complex model.

Next figures present the results of using VGG16 with imagenet weights with the last model.

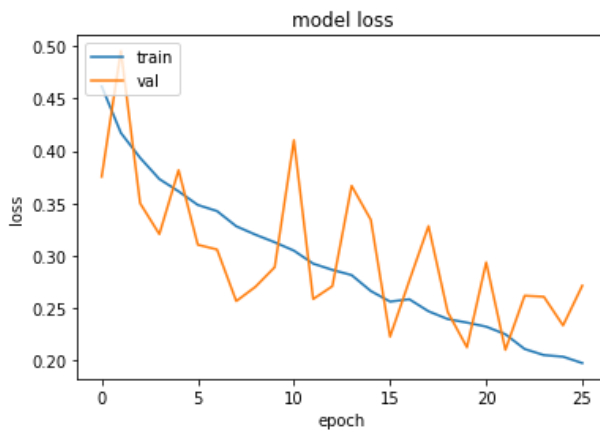


Figure 5.78: Model Loss

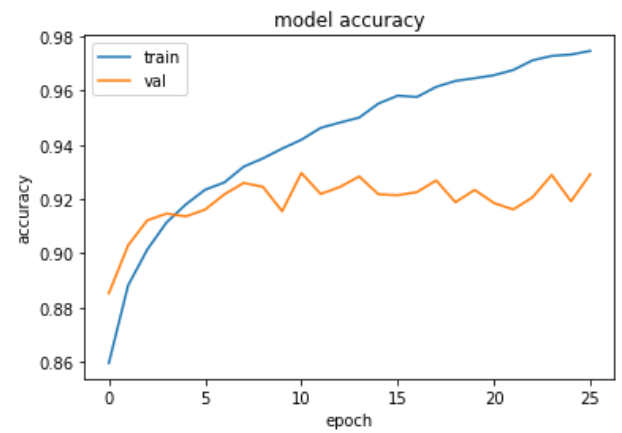


Figure 5.79: Model Accuracy

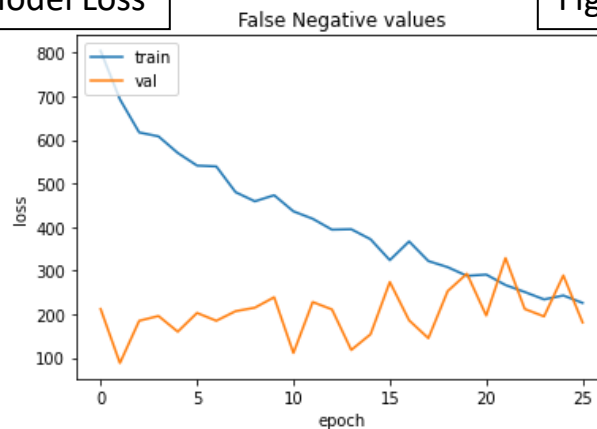


Figure 50: Model False negative values

From the above model is a very satisfied for us in the training accuracy needs to some improve in validation we will try to that next.

The test results of this model n the competition is 0.8837 as private score and 0.9149 as public score.

Like the other experiments we see that augmentation reduces the overfitting so we will add some augmentation like rotation range, width shift range, shear range.

Next figures represent the model results after adding the above augmentation and run with 24 epochs.

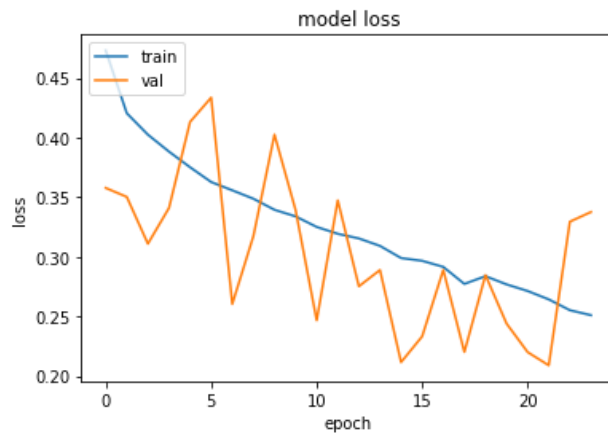


Figure 51: Model Loss

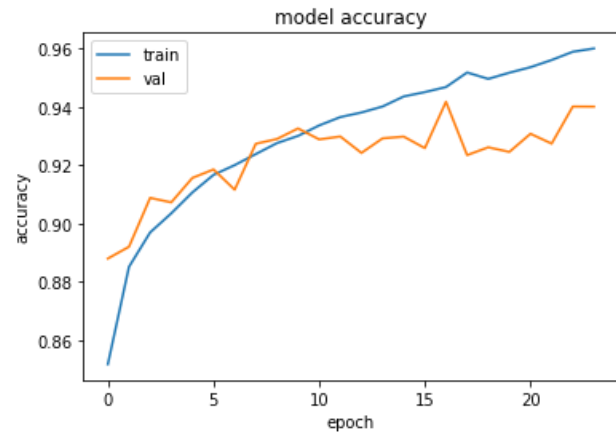


Figure 52: Model Accuracy

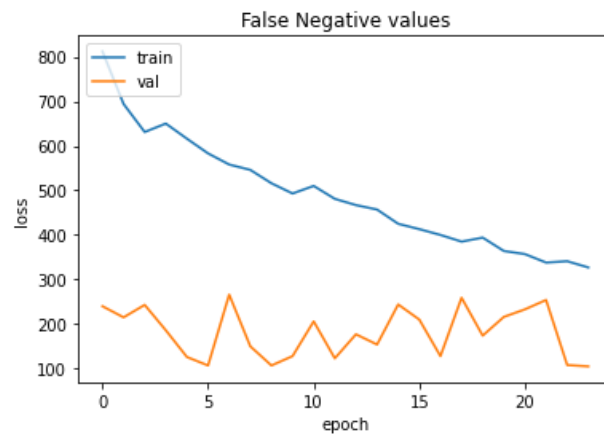


Figure 53: Model False negative values

The model train accuracy decreased, and the validation increased causing a great performance in reducing overfitting.

But we are computing the model performance according to the competition, so this model test results in the competition is 0.8813 as private score and 0.8890 as public score, again it reduces the overfitting.

Next, we will try a bigger network architecture (VGG19), trying to increase the accuracy.

The next figures present this model results after using a bigger network.

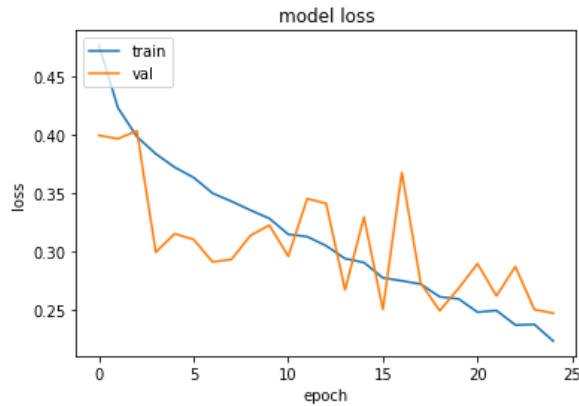


Figure 54: Model Loss

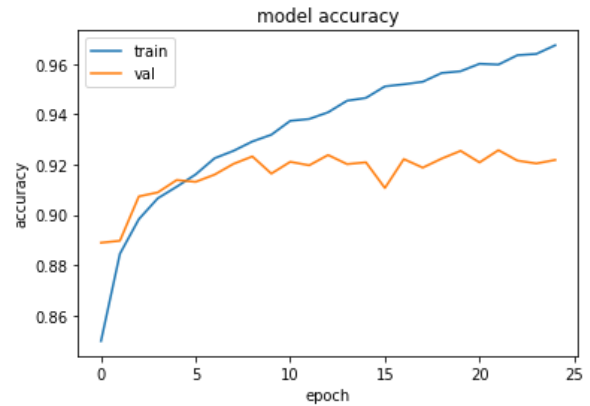


Figure 55: Model Accuracy

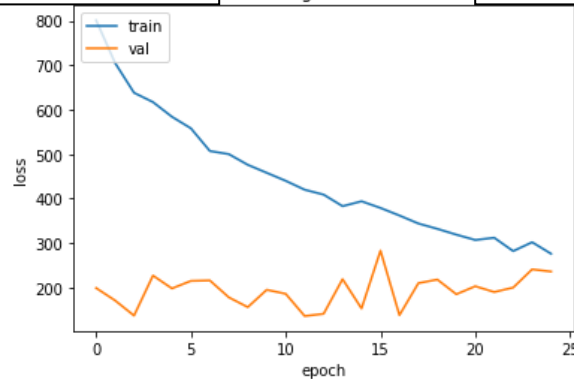


Figure 56: Model False negative values

The model performance goes down in both train and validation and return to the previous overfitting.

We will back the previous model and tune learning rate to reach to one close to the optimal learning rate, we try 12 values up and down the 10^{-5} learning rate and find that $10^{-5.5}$ is the closest of what we try this is according to the competition results.

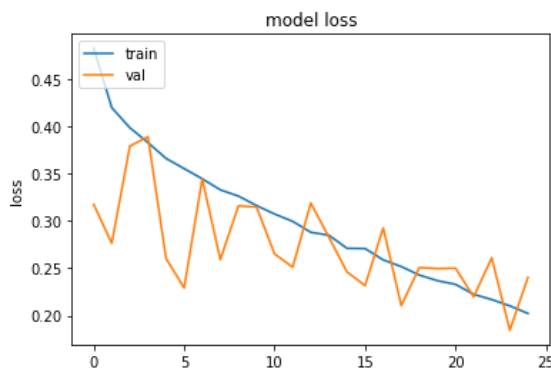


Figure 57: Model Loss

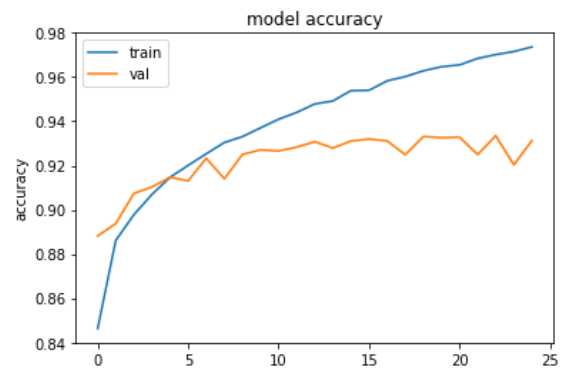


Figure 58: Model Accuracy

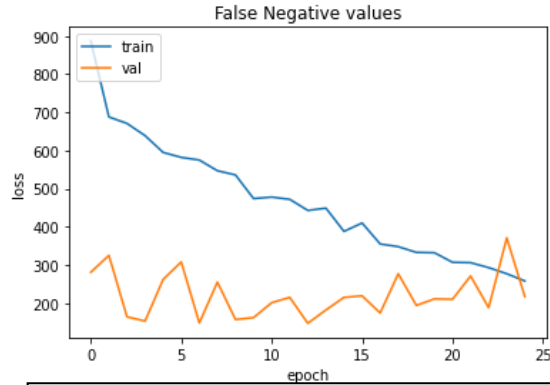


Figure 59: Model False negative values

The competition test results of this model are 0.8912 as private score and 0.8912 as public score.

Let's try removing dropouts and organize the values of dense layers units, we use this experiment with both learning rates 10^{-5} and $10^{-5.5}$ and let's see the results.

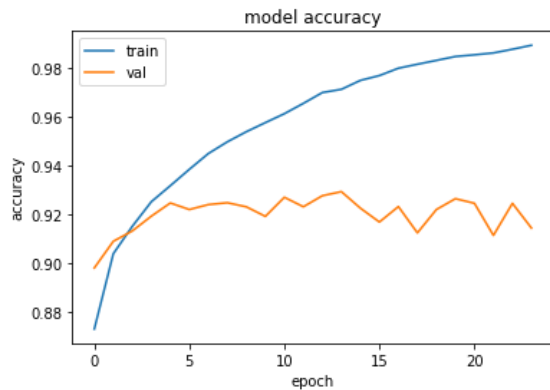


Figure 61: Model Accuracy with $LR=10^{-5.5}$

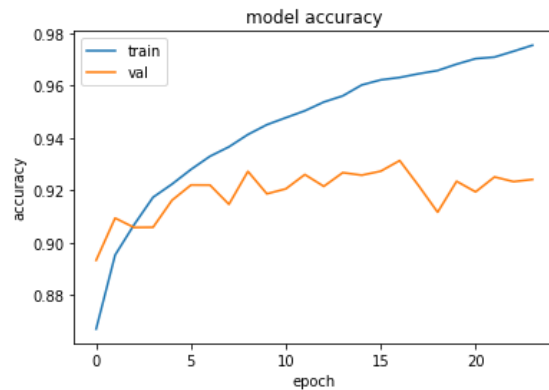


Figure 60: Model Accuracy with $LR=10^{-5}$

The test results of the first is 0.8873 as private score and 0.9092 as public score, for the second model private score is 0.8826 and public score is 0.8977.

From that we found that adding the dense layers and removing dropout changes the optimal learning rate for this model.

Next, we will try bigger networks like ResNet50, ResNet101 and DenseNet169, these network architectures share the complexity of connecting the previous layers with the next.

We will try the three architectures in parallel with different augmentation and with 10^{-5} learning rate.

The different augmentation back to the size of each network layers and the learning rate 10^{-5} because it was tried with architectures like this before and did will, next we will try the learning rate that tuned above.

- Use ResNet50 with imagenet weights and zoom, horizontal flip and vertical flip augmentation.

Next, figures represent the mode performance in 27 epochs.

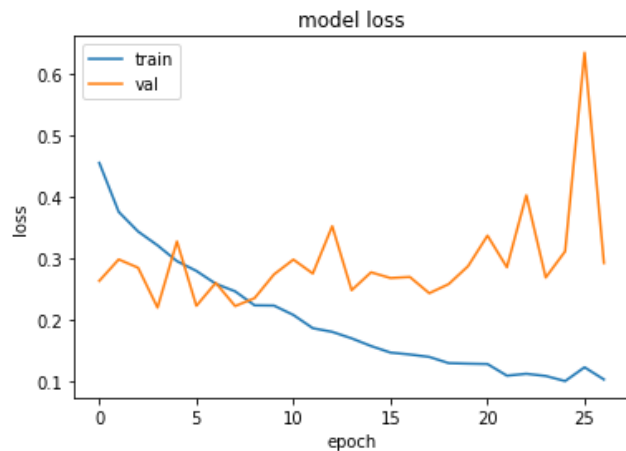


Figure 62: Model Loss

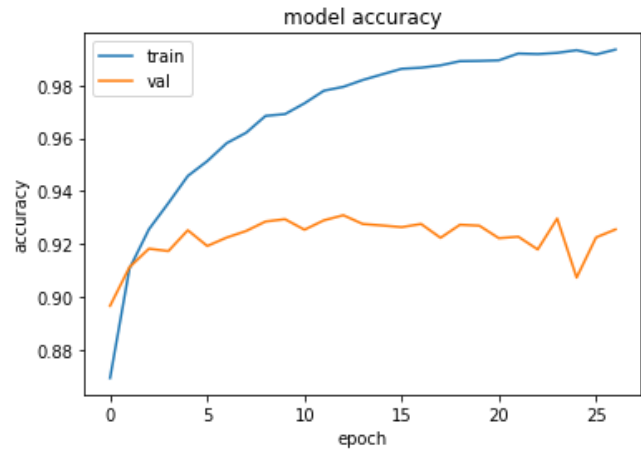


Figure 63: Model Accuracy

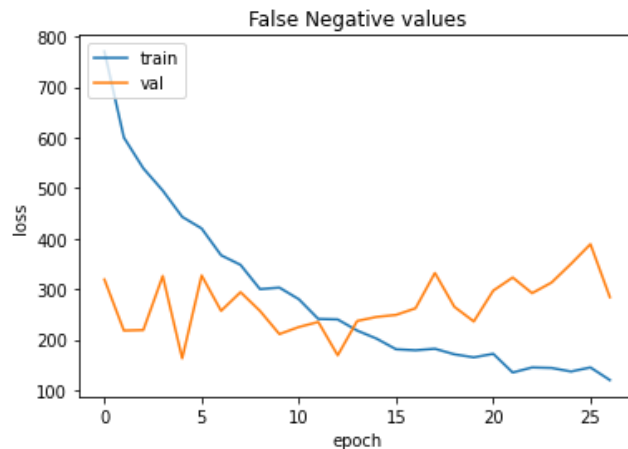


Figure 64: Model False negative values

The model performance is not good training goes in a very good way, but validation does not improve in the two accuracy matrices and loss function.

This model test results in the competition is 0.8695 as private score and 0.8755 as public score, the low performance defines in the competition results.

- Use ReseNet101 with imagenet weights, zoom, horizontal flip, vertical flip augmentation, width shift and height shift augmentation and add more dense layers.

Next, figures represent the mode performance in 27 epochs.

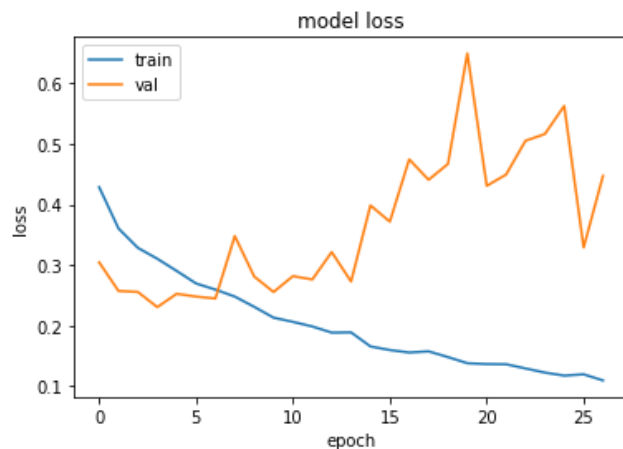


Figure 65: Model

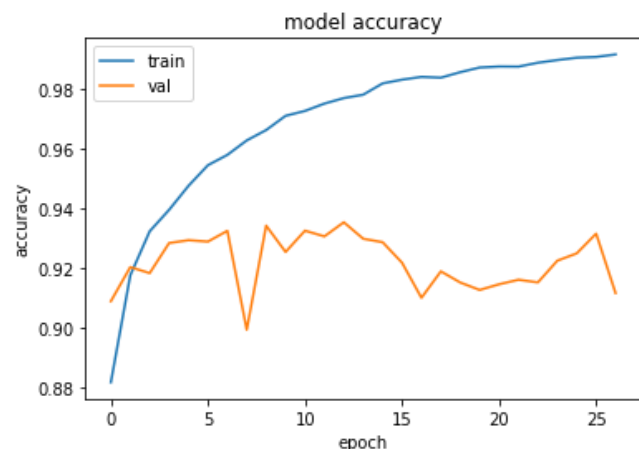


Figure 66: Model

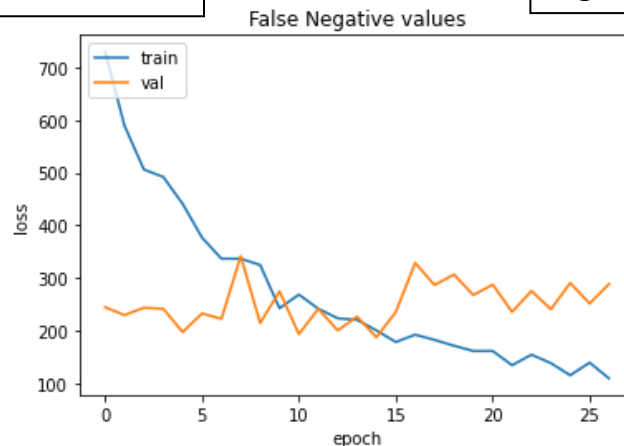


Figure 67: Model False negative values

The model doing very good in training and start good in validation at a certain epoch does not increase may decrease but we save the model with the highest accuracy may this affect the testing performance.

This model test results in the competition is 0.8901 as private score and 0.8645 as public score, the effect of saving the model with the highest accuracy presented in the competition results.

- Use DenseNet169 with imagenet weights, zoom, horizontal flip, vertical flip augmentation, width shift, height shift, shear augmentation and add more dense layers.

Next, figures represent the mode performance in 28 epochs.

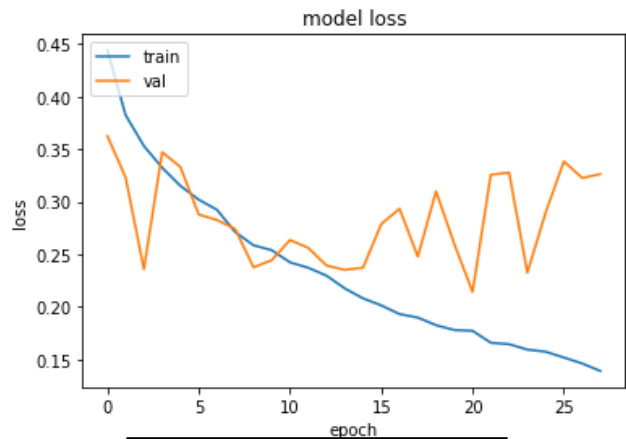


Figure 68: Model Loss

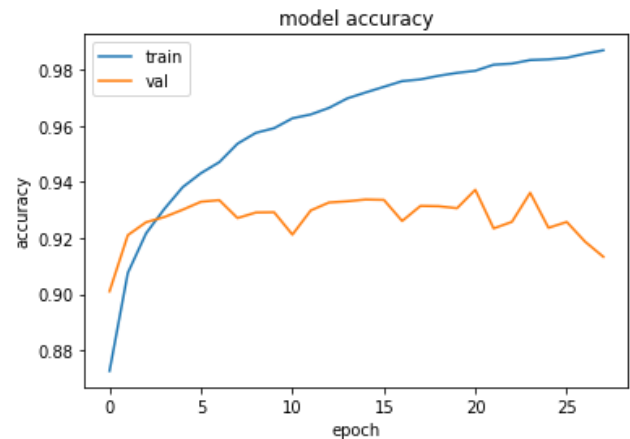


Figure 69: Model Accuracy

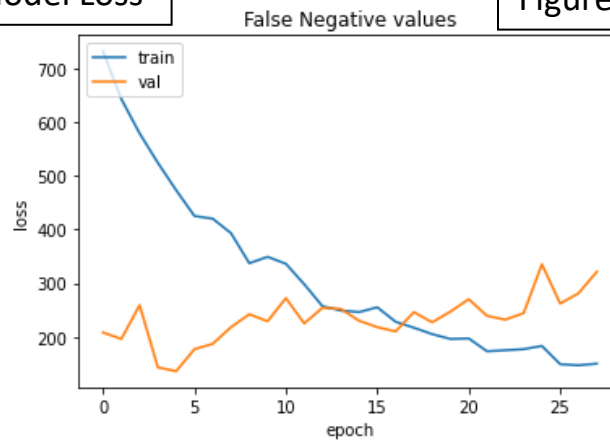


Figure 70: Model False negative values

Like the previous model this model doing very good in training and start good in validation at a certain epoch does not increase may decrease but we save the model with the highest accuracy may this affect the testing performance.

This model test results in the competition is 0.8819 as private score and 0.8984 as public score, the previous model doing good in the test results than this model.

The best model of the three in the test results is ResNet101, so we will us $10^{-5.5}$ learning rate to see the effect on this model.

Next, figures represent the mode performance in 26 epochs.

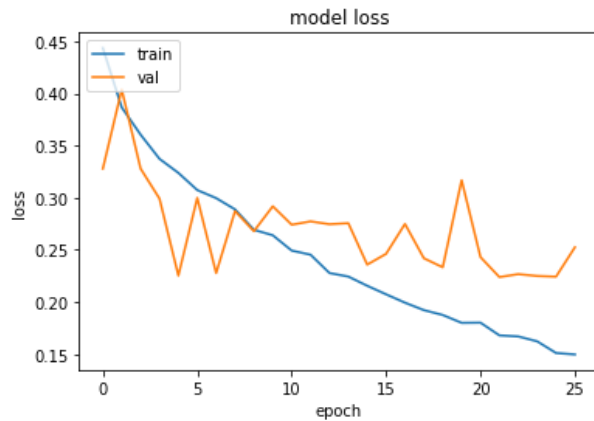


Figure 71: Model Loss

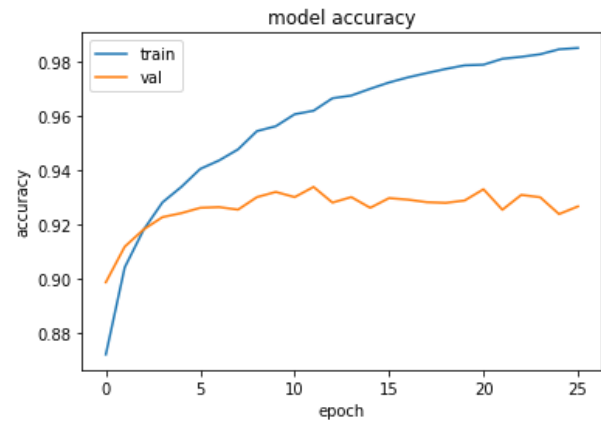


Figure 72: Model Accuracy

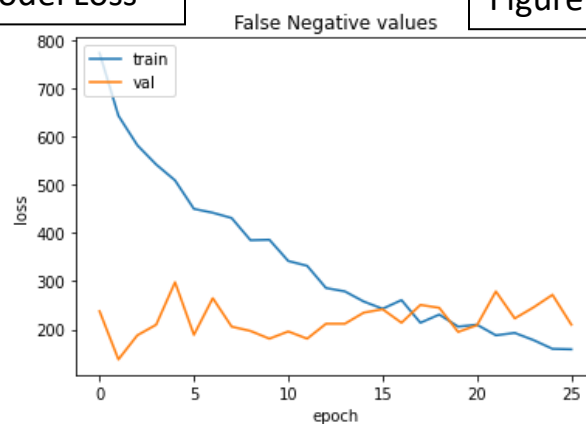


Figure 73: Model False negative values

The improves in the matrices and Loss Function comparing to it before changing the learning rate.

This model test results in the competition is 0.9002 as private score and 0.8784 as public score.

Now we are with the best model of all that we try in the experiments but we facing a problem that the size of the model as very big this is not good so we try a lot of techniques to reduce the make model with this accuracy and small size, some of the tried techniques are L1 regularization, L2 regularization, Drop out regularization...etc, all of these techniques reduces the size but the accuracy become so bad until we reach to using 1x1 convolution instead of dense layers and we did a lot of experiments on it until reaching to the next model.

Next, figures represent the mode performance in 28 epochs.

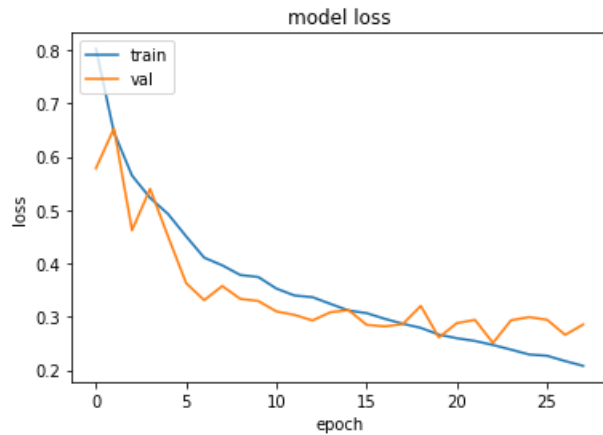


Figure 74: Model Loss

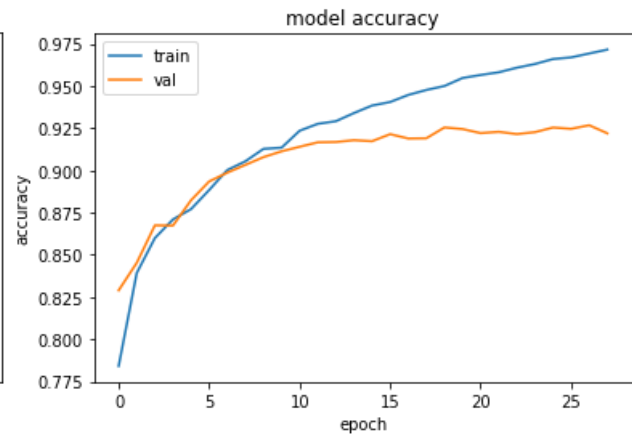


Figure 75: Model Accuracy

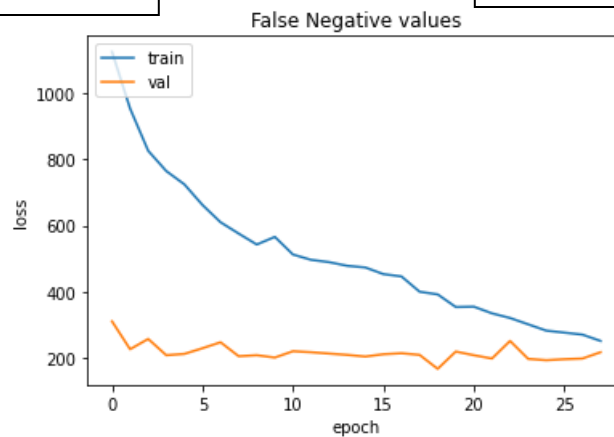


Figure 76: Model False negative values

As shown in the figure this model has a very good performance comparing to the last models.

This model test results in the competition is 0.9017 as private score and 0.8863 as public score, so we tend to go with this model and deploy it to connect with the embedded device.

We can say that the process of improving the model stopped now as its will not stop ever we enter the market to save lives in we will always work on this goal.

5.6 Deployment

5.6.1 API deployment

This section will concern on deploy Flask web server to Azure cloud services.

Firstly, we will create a Azure Virtual Machine (VM) which is one of several types of on-demand, scalable computing resources that Azure offers.

We will use here Ubuntu VM with Standard B2s size to deploy our Python application. Once we created the VM, we connect to the machine using SSH protocol. After that we cloned the API from GitHub.

```
$ git clone https://github.com/GProjet/Melanosizer.git
```

We navigate to API directory and run the following commands

```
$ sudo apt-get update
```

```
$ python3 -m venv env
```

```
$ source env/bin/activate
```

```
$ pip3 install -r requirements.txt
```

```
$ export FLASK_APP=run.py
```

```
$ python run.py
```

Now, our application is running at port **5000**.

Next, we will use **Nginx** reverse proxy server that takes a client request, passes it on to servers, and subsequently delivers the server's response back to the client. We install the nginx using the following command:

```
$ sudo apt install nginx -y
```

One last thing, we will use **Gunicorn** which is built so many different web servers can interact with it. It also does not really care what you used to build our web application.

To install Gunicorn use the following command:

```
$ pip3 install gunicorn
```

And we will navigate to API directory with running the following command

```
$ gunicorn --bind 0.0.0.0:5000 app:app
```

Now, our web server is running on port 80 based on Nginx reverse proxy within the VM.

5.6.2 Raspberry PI deployment

Use TensorFlow lite to convert the model into TensorFlow lite model to run the model on raspberry pi.

5.7 Embedded

5.7.1 Design Constraints:

In this chapter we discuss the basic design constraints, the hardware and software environments.

5.7.2 Hardware Environment

A melanosizer is a device drawn or powered by external electricity or connected to the pc or laptop.

it is controlled from the Raspberry Pi TOUCH SCREEN LCD HDMI that connects to the Raspberry Pi Camera Module to your Raspberry Pi and take pictures, record video, and apply image effects.

5.7.3 Raspberry Pi TOUCH SCREEN LCD HDMI

1. Raspberry pi touch display

The Raspberry Pi Touch Display is an LCD display which connects to the Raspberry Pi through the DSI connector. In some situations, it allows for the use of both the HDMI and LCD displays at the same time (this requires software support)

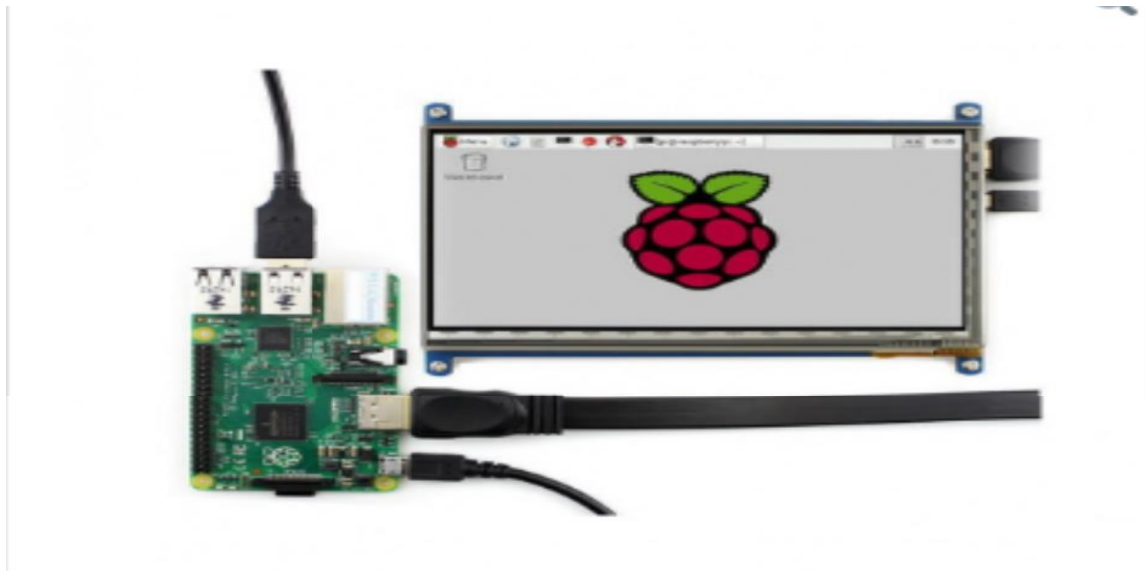


Figure 5. 9 Raspberry Pi TOUCH SCREEN LCD HDMI

2. Board support

The DSI display is designed to work with all models of Raspberry Pi, however early models that do not have mounting holes (the Raspberry Pi 1 Model A and B) will require additional mounting hardware to fit the HAT-dimensioned stand-offs on the display PCB

3. Physical Installation

The following image shows how to attach the Raspberry Pi to the back of the Touch Display (**if required**), and how to connect both the data (**ribbon cable**) and power (**red/black wires**) from the Raspberry Pi to the display. If you are not attaching the Raspberry Pi to the back of the display, take extra care when attaching the ribbon cable to ensure it is the correct way round. The black and red power wires should be attached to the GND and 5v pins respectively.

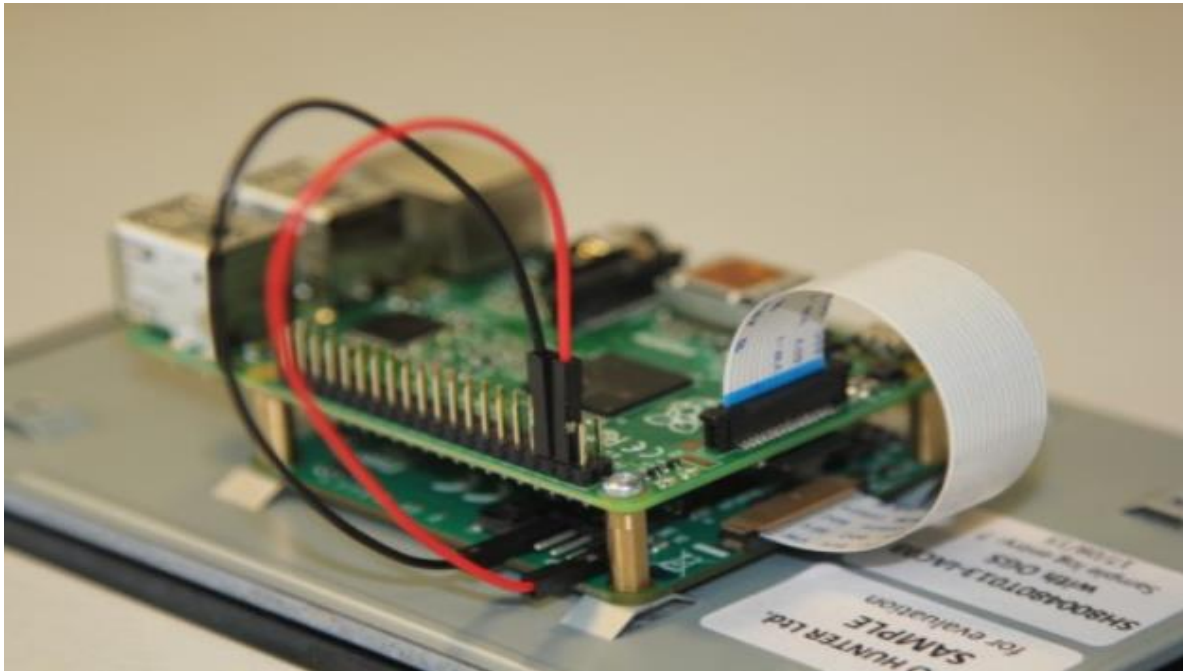


Figure 5. 10 Raspberry Pi to the back of the Touch Display

4. specifications

- 800×480 RGB LCD display
- 24-bit colour
- Industrial quality: 140-degree viewing angle horizontal, 130-degree viewing angle vertical
- 10-point multi-touch touchscreen
- PWM backlight control and power control over I2C interface
- Metal-framed back with mounting points for Raspberry Pi display conversion board and Raspberry Pi
- Backlight lifetime: 20000 hours
- Operating temperature: -20 to +70 degrees centigrade
- Storage temperature: -30 to +80 degrees centigrade
- Contrast ratio: 500
- Average brightness: 250 cd/m2

- Power requirements: 200mA at 5V typical, at maximum brightness.

5.8 Raspberry Pi computer with a Camera Module port

All current models of Raspberry Pi have a port for connecting the Camera Module.

There are two versions of the Camera Module:

- The standard version, which is designed to take pictures in normal light
- The NoIR version, which doesn't have an infrared filter, so you can use it together with an infrared light source to take pictures in the dark

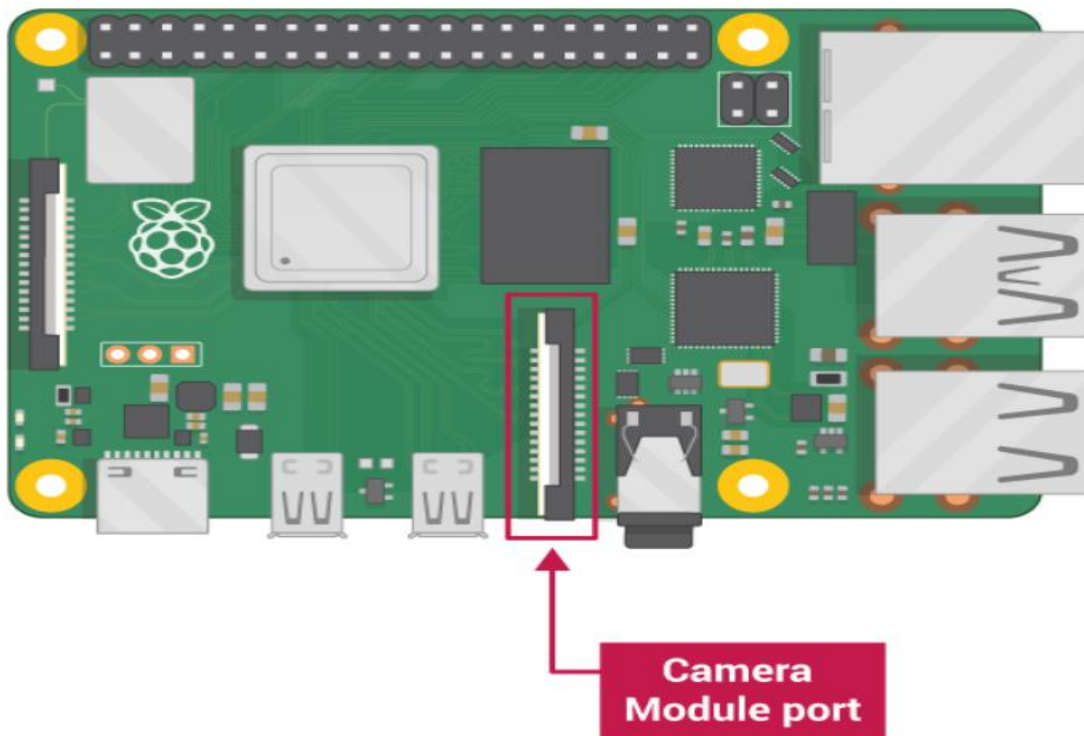


Figure 5. 11 camera module port

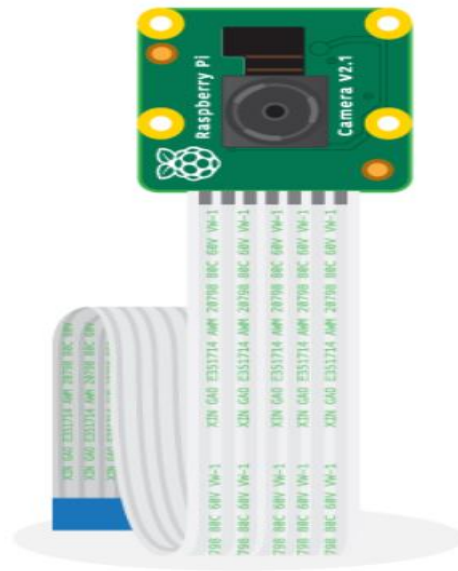


Figure 5. 12 Raspberry Pi Camera Module

1. Connect the Camera Module

Ensure your Raspberry Pi is turned off.

1. Locate the Camera Module port
2. Gently pull up on the edges of the port's plastic clip
3. Insert the Camera Module ribbon cable; make sure the connectors at the bottom of the ribbon cable are facing the contacts in the port.
4. Push the plastic clip back into place

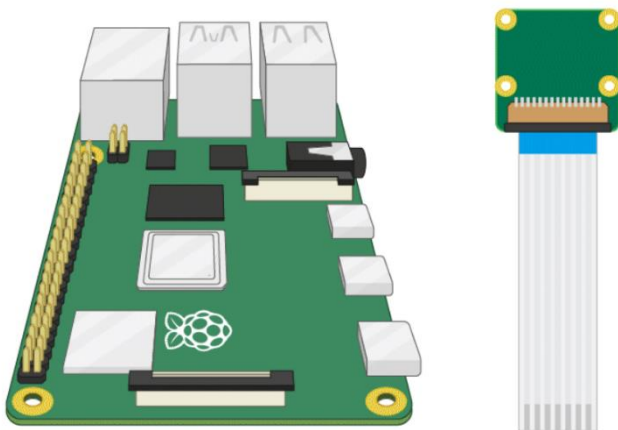


Figure 5. 13 shows step 1,2



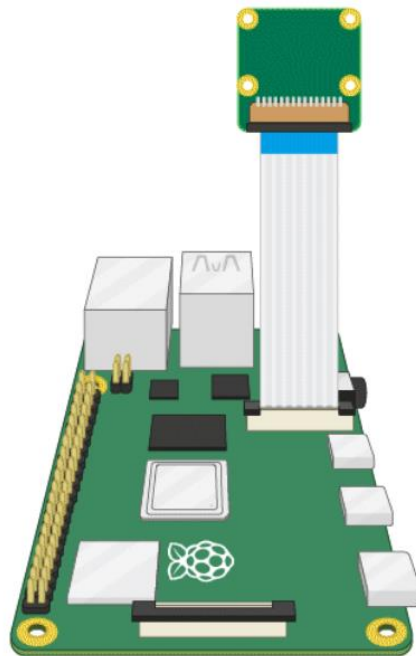


Figure 5. 14 show steps 3,4

- Start up your Raspberry Pi.
- Go to the main menu and open the Raspberry Pi Configuration tool.

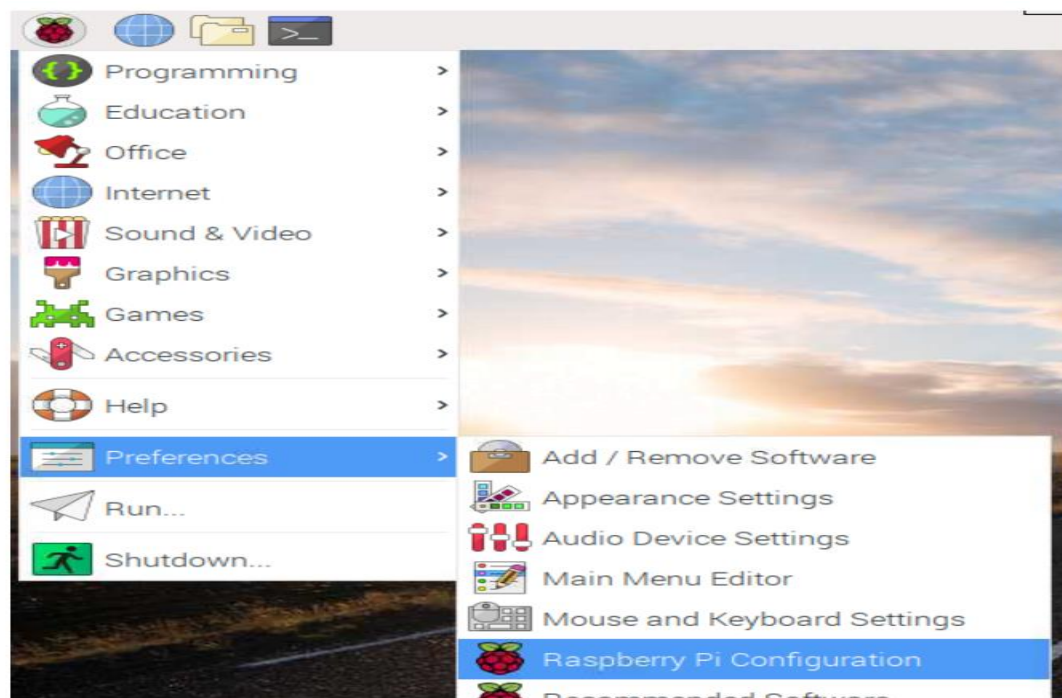

- Select the Interfaces tab and ensure that the camera is enabled:

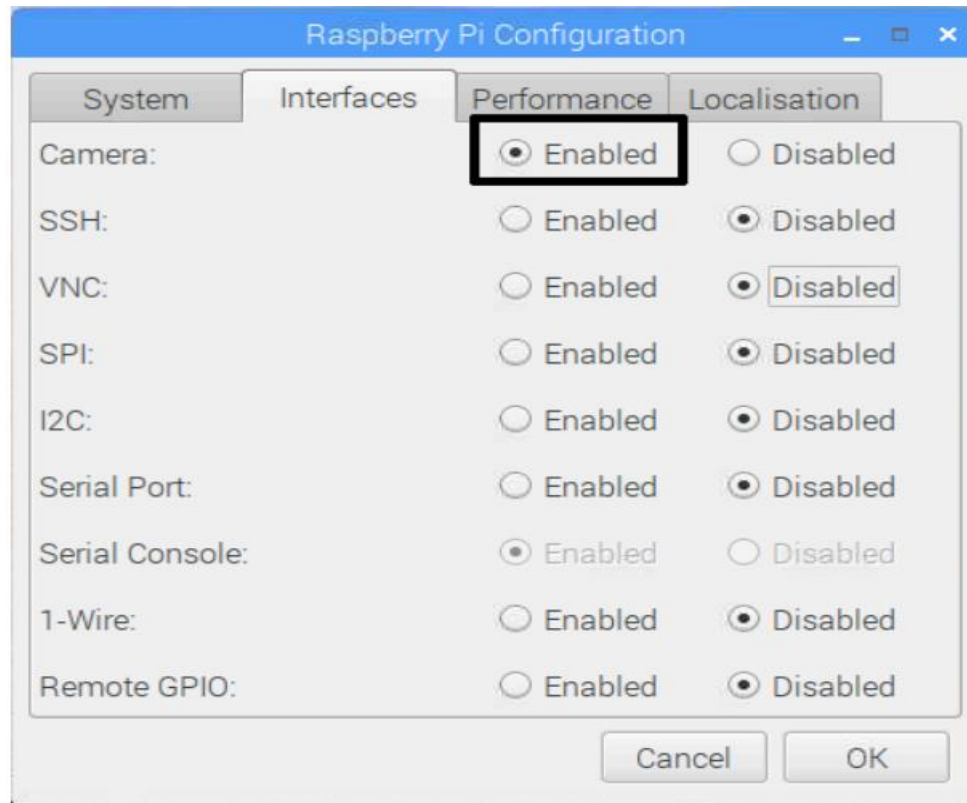

- Reboot your Raspberry Pi.

Raspberry Pi 4 Model B

The Raspberry Pi 4 Model B (Pi4B) is the first of a new generation of Raspberry Pi computers supporting more RAM and with significantly enhanced CPU, GPU and I/O performance; all within a similar form factor, power envelope and cost as the previous generation Raspberry Pi 3B+. The Pi4B is available with either 1, 2 and 4 Gigabytes of LPDDR4 SDRAM.

1. GPIO

A powerful feature of the Raspberry Pi is the row of GPIO (general-purpose input/output) pins along the top edge of the board. A 40-pin GPIO header is found on all current Raspberry Pi boards (unpopulated on Raspberry Pi Zero, Raspberry Pi Zero W and Raspberry Pi Zero 2 W). Prior to the Raspberry Pi 1 Model B+ (2014), boards comprised a shorter 26-pin header.

Voltages


Two 5V pins and two 3.3V pins are present on the board, as well as a number of ground pins (0V), which are unconfigurable. The remaining pins are all general purpose 3.3V pins, meaning outputs are set to 3.3V and inputs are 3.3V-tolerant.

Outputs

A GPIO pin designated as an output pin can be set to high (3.3V) or low (0V).

Inputs

A GPIO pin designated as an input pin can be read as high (3.3V) or low (0V). This is made easier with the use of internal pull-up or pull-down resistors. Pins GPIO2 and GPIO3 have fixed pull-up resistors, but for other pins this can be configured in software.

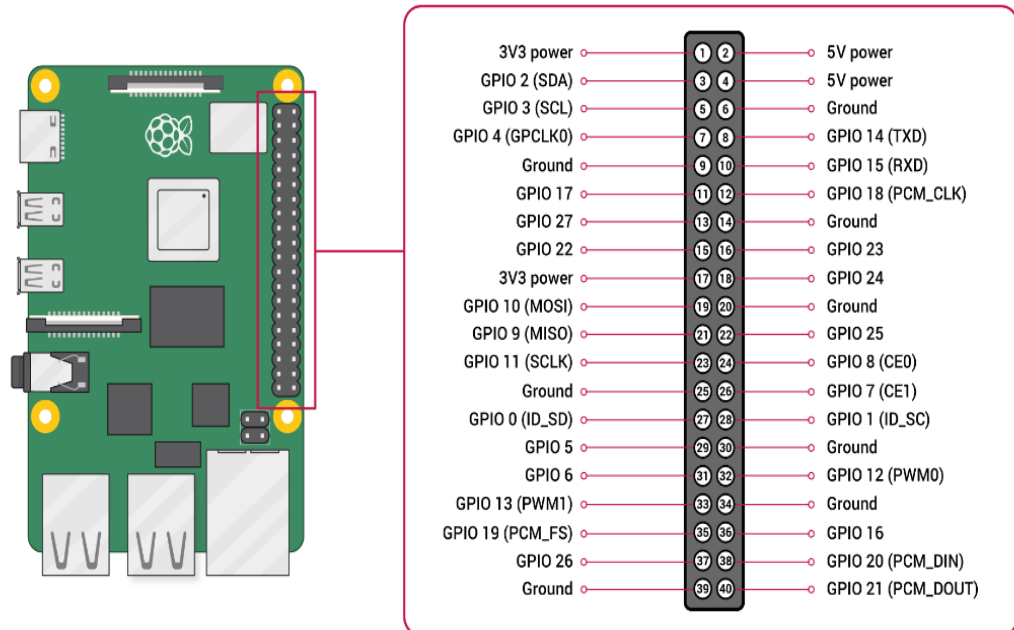

5.9 Software Environment

5.9.1 Operating system

Raspberry Pi OS is a free operating system based on Debian, optimised for the Raspberry Pi hardware, and is the recommended operating system for normal use on a Raspberry Pi. The OS comes with over 35,000 packages: precompiled software bundled in a nice format for easy installation on your Raspberry Pi. Raspberry Pi OS is under active development, with an emphasis on improving the stability and performance of as many Debian packages as possible on Raspberry Pi.

5.9.2 Updating Raspberry Pi OS

It's important to keep your Raspberry Pi up to date. The first and probably the most important reason is security. A device running Raspberry Pi OS contains

millions of lines of code that you rely on. Over time, these millions of lines of code will expose well-known vulnerabilities. The only way to mitigate these exploits as a user of Raspberry Pi OS is to keep your software up to date.

The second reason, related to the first, is that the software you are running on your device most certainly contains bugs. Some bugs are CVEs, but bugs could also be affecting the desired functionality without being related to security. By keeping your software up to date, you are lowering the chances of hitting these bugs.

5.9.3 Using APT

The easiest way to manage installing, upgrading, and removing software is using APT (Advanced Packaging Tool) from Debian. To update software in Raspberry Pi OS, you can use the apt tool from a Terminal window.

5.9.4 Using rpi-update

rpi-update is a command line application that will update your Raspberry Pi OS kernel and VideoCore firmware to the latest pre-release versions.

6 Chapter Six (Testing and Validation)

This chapter introduces you to the testing setup and how we tested each module in isolation and the end-to-end integration that we planned and followed to prove to a high certainty the correctness of the system.

Testing of the project is done to check whether the actual results match what is expected, ensure that each module is doing what it is intended to do, and to ensure that the system as a whole is working correctly. Our project depends on real world skin cancer images, so, most of our testing was done through testing the modules on different images and checking whether the output is satisfactory compared to the results of similar projects.

In the end, we test the workflow of the project to check whether the results are clear enough to interpret or not.

6.1 Module Testing

Model testing is done sequentially into three steps: Model testing which includes testing metrics of deep neural network. Hardware testing is concerned with Raspberry Pi device components testing separately. Web server testing includes testing the API (Application Programming Interface) which built using Flask to accelerate the development of web applications.

6.1.1 Model Testing

- Model metrices

Two matrices are used in this model AUC (Accuracy Under the Curve) and False Negative values.

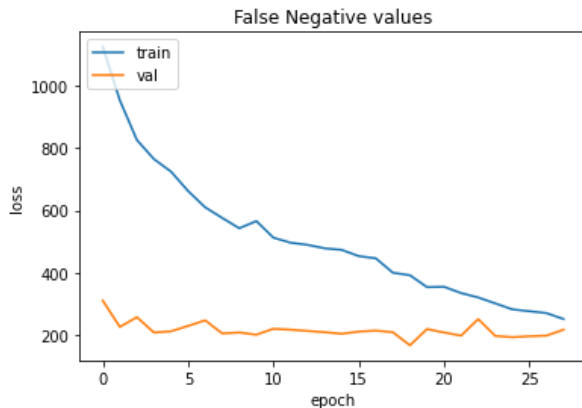


Figure 6.1: False Negative values matrix

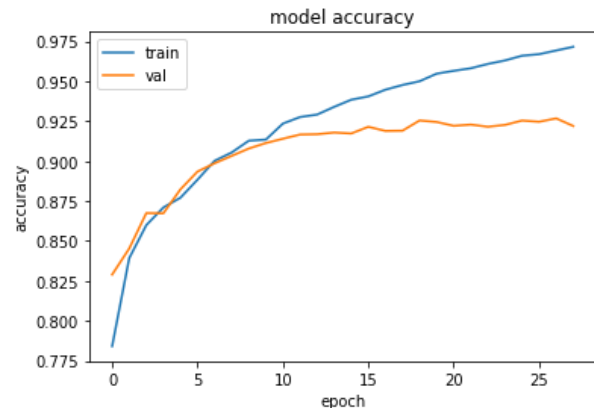


Figure 6.2: AUC matrix

6.1.2 Web Server Testing

Since the web server is built using Flask, we document the API with Swagger UI framework which depends on OpenAPI technology to exchanges the requests between the client and the web server.

To be able to predict an image whether the tumor is benign or malignant. We have created an end point called predict which takes an image as input and return string results as output. The post method is allowed only images as input and return an error if else.

Figure 6.3 shows the Swagger UI, it has a documentation of the API which has an example of predict endpoint.

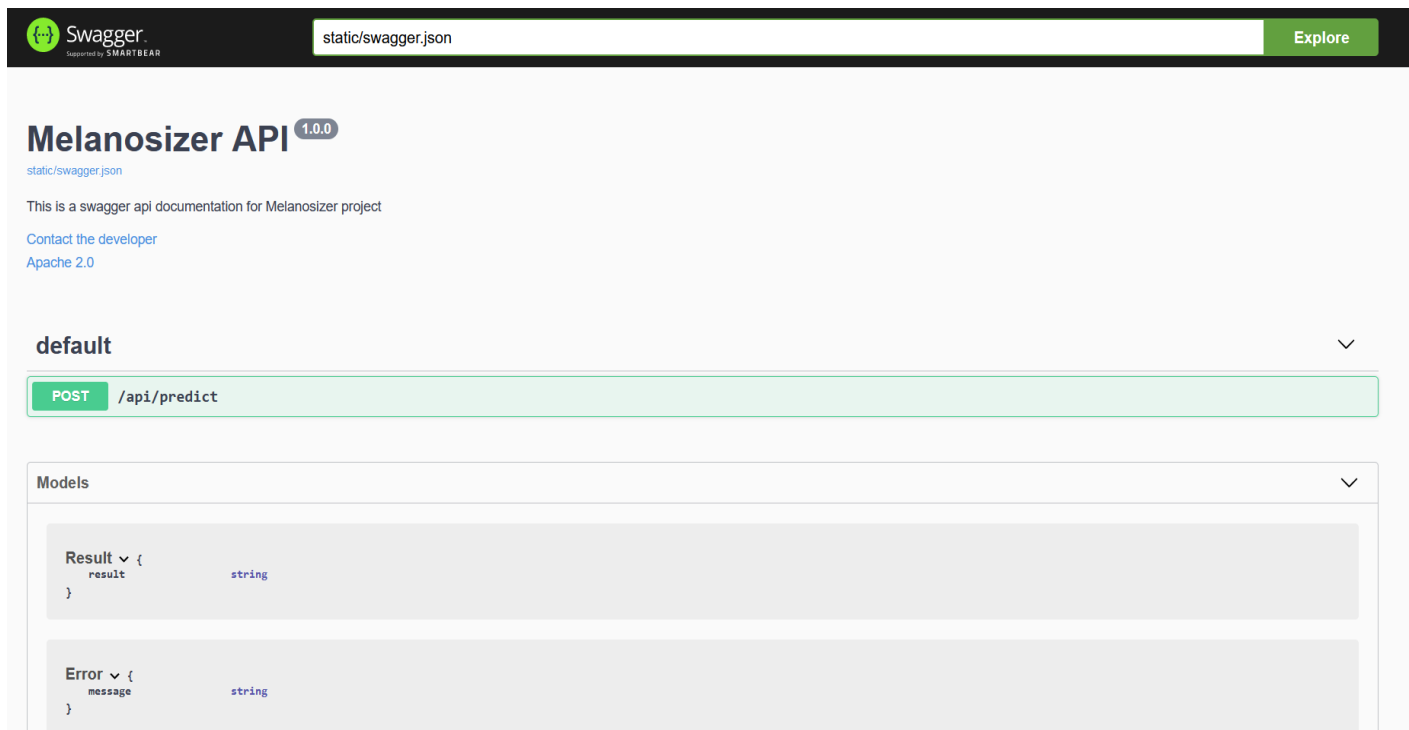


Figure 6. 3 Swagger UI documentation of API

Figure 6.4 shows the prediction result of tumor image to determine whether is benign or malignant.

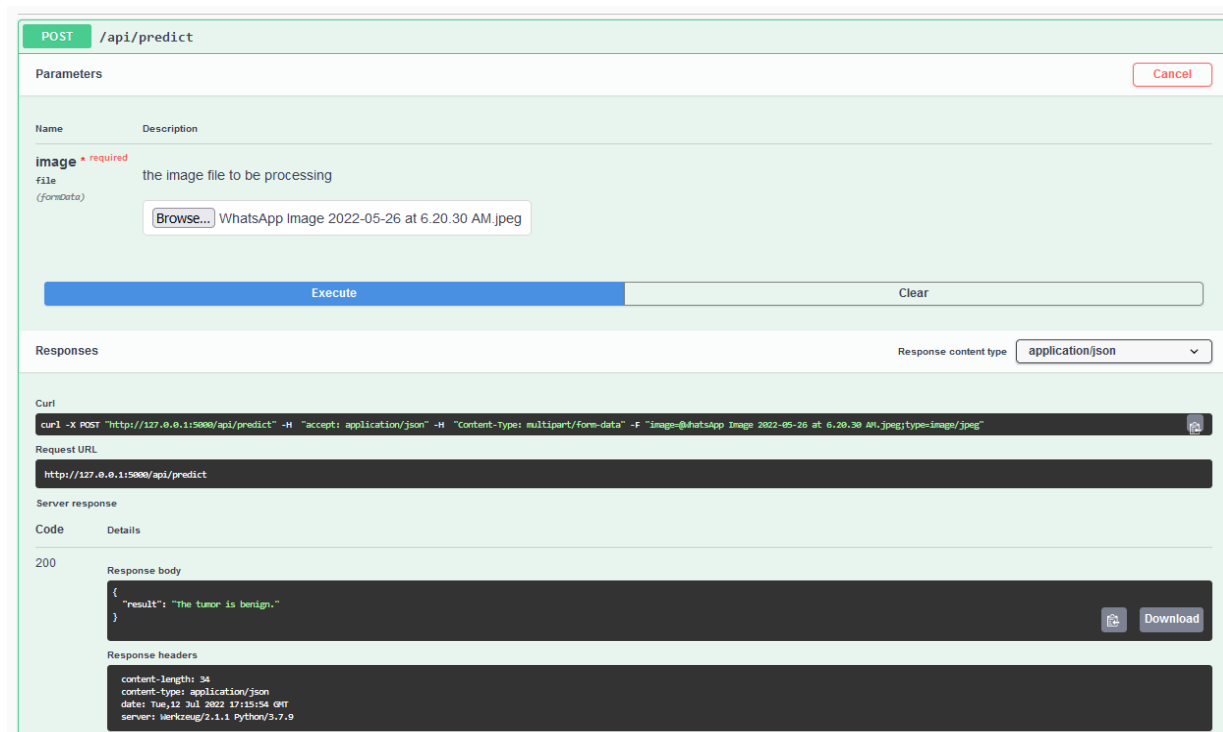


Figure 6. 4 Example of skin image prediction

6.2 Overall Testing

After merging all modules together and testing the whole project integrated with Raspberry PI device, we started with takes an image of tumor by camera and send to the web server to be predicted and return the results finally to the raspberry pi screen.

Figure 6.5 shows the input sample taken by Raspberry PI device



Figure 6. 5 Screen on input image

Figure 6.6 shows the result/output of prediction

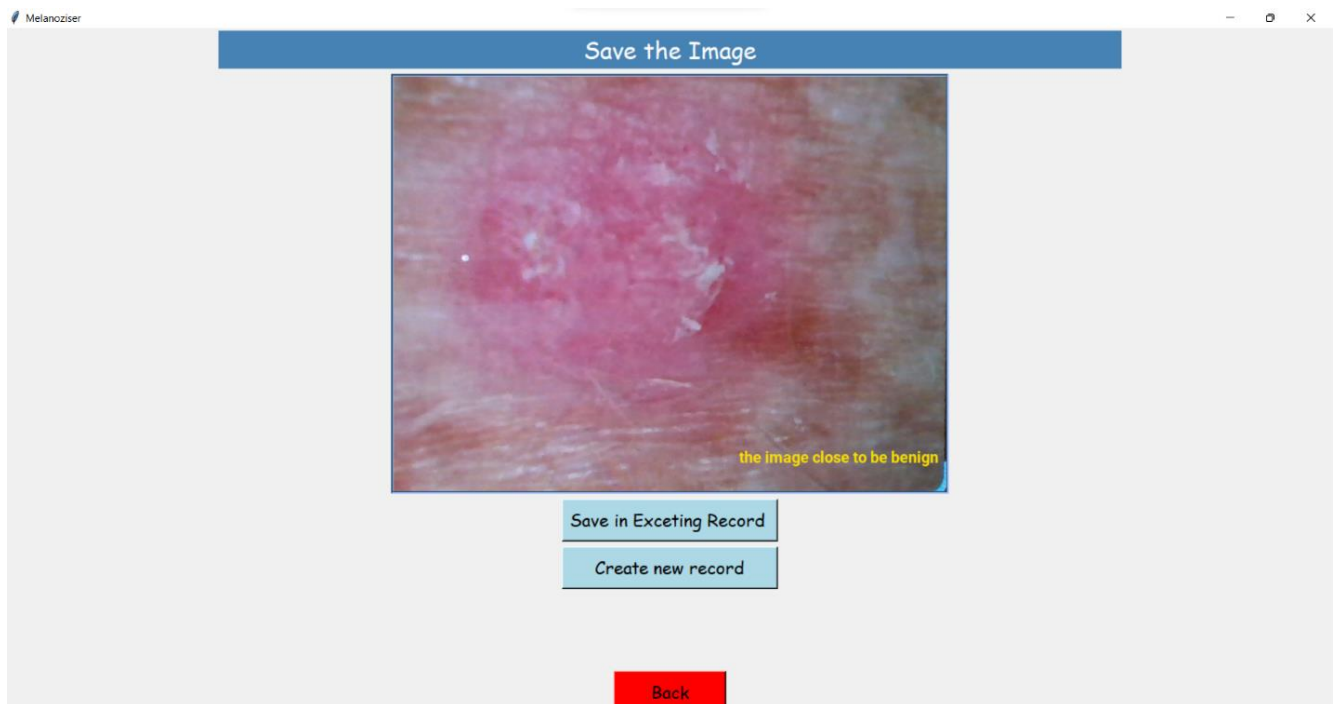


Figure 6. 6 Screen on output prediction

7 Chapter Seven (Tools and Technologies)

in this chapter we Describe any state-of-the-art tools and technologies you used in the project.

7.1 Libraries and environment used in ML model

Normally, a library is a collection of books or is a room or place where many books are stored to be used later. Similarly, in the programming world, a library is a collection of precompiled codes that can be used later on in a program for some specific well-defined operations. Other than pre-compiled codes, a library may contain documentation, configuration data, message templates, classes, and values, etc.

1. libraries

➤ TensorFlow compat .v1

TensorFlow is an open-source library used for high-level computations It is also used in machine learning and deep learning algorithms. It contains a large number of tensor operations.

We use TensorFlow compatible version 1 because it contains functions that we will need but another library TensorFlow doesn't have some required functions.

➤ **Matplotlib**

it is an open-source library, is responsible for plotting numerical data. And that's why it is used in data analysis

➤ **Pandas**

it is an open-source machine learning library that provides flexible high-level data structures and a variety of analysis tools

➤ **Numpy**

It is a popular machine learning library that supports large matrices and multi-dimensional data, Used for numerical operations

works with **SciPy** to handle complex computations

➤ **Scikit-learn**

works in association with Numpy and SciPy.

➤ **Keras**

It's TensorFlow API

➤ **Random**

➤ **TKinter**

open source, portable graphical user interface (GUI) library designed for use in Python scripts.

➤ **Pillow**

Imaging Library (expansion of PIL) is the de facto image processing package for Python language. It incorporates lightweight image processing tools that aids in editing, creating and saving images.

➤ **Datetime**

supplies classes to work with date and time

➤ **OpenCV**

library of Python bindings designed to solve computer vision problems.

2. environments

➤ **Colab**

Python in google colab for online GPU as environment

➤ **Kaggle Notebooks**

an online community platform for data scientists and machine learning enthusiasts as environment.

➤ **AZURE ML studio**

An online platform provides a powerful and trusted machines to run your model, and do many other things.

➤ **THonny IDE**

Python IDE supported by Raspberry Pi Foundation

8 Chapter eight (Future work)

8.1 Conclusion

Finally, we get an embedded device that can be used by a doctor to capture skin picture from skin cancer patient and detect whether it is melanoma or other type of skin cancer, and produce deep learning model to implement this process and diagnose the image and find use 1x1 convolution with the model to reduce power consumption.

8.2 Future Work

we have to work on improving our model more and more to enter the market and make the doctors and the society happy to see this improvement in melanoma diagnosis devices.

We will work on improving some points:

1. Integrate the device more with the medical devices by making the device work with Reflectance confocal microscopy (RCM) instead of a camera.

Reflectance confocal microscopy (RCM) enables imaging of skin lesions at cellular level resolution at the bedside (in vivo) or in freshly excised tissue (ex vivo).

RCM features of common melanocytic and non-melanocytic skin neoplasms such as melanoma, actinic keratosis/squamous cell carcinoma, basal cell carcinoma, and nevi have been well defined and show good correlation with dermoscopic and histopathologic findings. Due to its technical properties, RCM is especially suitable for the examination of flat skin lesions.

This will take the performance to another level as the details that the Reflectance confocal microscopy (RCM) image will enable us to work with it.

2. Collect more images to the data set to improve the model.

As long as the device work in a Radiology and analysis center we will see what we will have to do to take the pictures diagnosed on the device to do two things collect more image to the data set and do error analysis on these images to improve model performance by seeing the four matrices to all image the mode diagnosis (True Positive, False Positive, True Negative, False Negative).

3. Work with many models instead of one model.

Work on the image with many techniques before the recognition, this needs a good processing unit to process these models so we will work on improving the work of the processor by integrating it with a GPU.

Try to integrate with A powerful cloud machine this need to handle the financial problem by collaborating with some medical companies or other ways to reach the target.

4. Try to go with feature engineering deep learning instead of end-to-end deep learning.

End-to-end deep learning requires a huge data set which not available now we work on improve the data set but the collecting of data set in medical field is not an easy way.

Feature engineering deep learning did a very good work in these cases, but it needs a hard work on engineering the features (split the end-to-end model into some features and work with many steps instead of only one step (taking image and produce output)).

9 References

1. "Skin Cancer Treatment (PDQ®)". NCI. 25 October 2013. Archived from the original on 5 July 2014. Retrieved 30 June 2014.
2. Stahl W, Sies H (November 2012). "β-Carotene and other carotenoids in protection from sunlight". *The American Journal of Clinical Nutrition*. 96 (5): 1179S–84S.
3. Kuschal C, Thoms KM, Schubert S, Schäfer A, Boeckmann L, Schön MP, Emmert S (January 2012). "Skin cancer in organ transplant recipients: effects of immunosuppressive medications on DNA repair". *Experimental Dermatology*.
4. Saladi RN, Persaud AN (January 2005). "The causes of skin cancer: a comprehensive review". *Drugs of Today*. 41 (1): 37–53.
5. Lozano R, Naghavi M, Foreman K, Lim S, Shibuya K, Aboyans V, et al. (December 2012). "Global and regional mortality from 235 causes of death for 20 age groups in 1990 and 2010: a systematic analysis for the Global Burden of Disease Study 2010" (PDF). *Lancet*. 380 (9859): 2095–128.
6. Schulz, Hannes; Behnke, Sven (1 November 2012). "Deep Learning". *KI - Künstliche Intelligenz*. 26 (4): 357–363.
7. Bengio, Yoshua; LeCun, Yann; Hinton, Geoffrey (2015). "Deep Learning". *Nature*. 521 (7553): 436–444.
8. Ciresan, D.; Meier, U.; Schmidhuber, J. (2012). "Multi-column deep neural networks for image classification". *2012 IEEE Conference on Computer Vision and Pattern Recognition*. pp. 3642–3649.
9. Cireşan, Dan Claudiu; Meier, Ueli; Gambardella, Luca Maria; Schmidhuber, Jürgen (21 September 2010). "Deep, Big, Simple Neural Nets for Handwritten Digit Recognition". *Neural Computation*. 22 (12): 3207–3220.
10. Wallach, Izhar; Dzamba, Michael; Heifets, Abraham (9 October 2015). "AtomNet: A Deep Convolutional Neural Network for Bioactivity Prediction in Structure-based Drug Discovery".
11. Littleton, Cynthia (March 15, 2019). "CNN Original Series Ride News Tide to Multiplatform Success". *Variety*. Archived from the original on September 29, 2019. Retrieved May 5, 2019.
12. Petski, Denise (April 11, 2018). "CNN Adds Six New Original Series To 2019 Slate; Projects From Sanjay Gupta, Vox Media, More". *Deadline*. Archived from the original on May 5, 2019. Retrieved May 5, 2019.

13. Robbins, Stephanie. "TV Week September 6, 2007 CNN HD Debuts".
Tweek.com. Archived from the original on October 15, 2013. Retrieved
October 12, 2013.
14. Alfonso, Fernando. "CNN Center in Atlanta damaged during protests". CNN.
Archived from the original on May 30, 2020. Retrieved January 16, 2021.
15. "CNN Center". CNN. Archived from the original on January 21, 2021. Retrieved
January 16, 2021.
16. Dickson, Glen (December 15, 2008). "CNN Gets New Graphic Look".
17. Dekel, O. Shamir, and L. Xiao. Learning to classify with missing and corrupted
features. Machine Learning, 81(2):149–178, 2010.
18. Krizhevsky. Learning multiple layers of features from tiny images. Technical
report, University of Toronto, 2009.
19. Shrutik Katchi and Pritish Sachdeva, "A Review Paper on Raspberry Pi",
Vol.4, No.6, Dec 2014
20. Raspberry Pi for Dummies by Sean McManus, Mike Cook · 2013
21. Joseph Muniz, Aamir Lakhani, "Penetration Testing with Raspberry Pi", 2015.
22. Khan, Sikandar; Akram, Adeel; Usman, Nighat (2020). Real Time Automatic
Attendance System for Face Recognition Using Face API and OpenCV. Wireless
Personal Communications,
23. Helmi, R., Yusuf, S., & Jamal, A. (2019). Face recognition automatic class
attendance system (FRACAS). In IEEE international conference on automatic
control and intelligent systems (I2CACIS 2019), Selangor, Malaysia, June 29,
2019.
24. Varadharajan, E., Dharani, R., Jeevitha, S., Kavinmathi, B., & Hemalatha, S.
(2016). Automatic attendance management system using face detection. In
Online international conference on green engineering and technologies (IC-
GET).
25. Rosebrock, A. (2018). YOLO object detection with OpenCV. Py Image Search.
Retrieved December 20, 2018.
26. McCluskey, C. P., Bynum, T. S., & Patchin, J. W. (2004). Reducing chronic
absenteeism: An assessment of an early truancy initiative. NCCD News, 50(2),
214–234.
27. Yadav, V., & Bhole, G. P. (2019). Cloud based smart attendance system for
educational institutions. In 2019 international conference on machine learning,
big data, cloud and parallel computing (Com-ITCon), India, Feb 14th–16th, 2019.